

**NICO<sup>2</sup>**  
**NON-INVASIVE  
CARDIAC OUTPUT**

**User's Manual**

**Non-Invasive Cardiac Output Monitor  
Model 7300**

February 1, 2001

Catalog No. 9226-23-05







## Thank you ...

---

Thank you for purchasing the NICO® Non-Invasive Cardiac Output monitor from Novamatrix.

NICO® measures cardiac output through respiratory gas analysis based on the well accepted Fick Principle, providing continuous and accurate display of cardiac output. The monitor also operates in Respiratory Mechanics-only mode, providing the clinician with a respiratory profile of the patient through a combination of capnography, airway flow and pressure, and pulse oximetry.

We expect that you will find the application and use of NICO® extremely simple, making it easy to adopt this exciting technology into your clinical practice. NICO® can provide accurate cardiac output values without the need for invasive procedures, benefitting the patient, the clinician, and the health care system in general.

We appreciate your patronage and look forward to developing a long-term relationship with you and your institution.

Sincerely,



USA TOLL FREE 1-800-243-3444  
PHONE 205-265-7701  
FAX 205-284-0753

WORLD WIDE WEB:  
<http://www.novamatrix.com>

E-MAIL:  
Customer Service  
[sales@novamatrix.com](mailto:sales@novamatrix.com)  
Technical Service  
[techline@novamatrix.com](mailto:techline@novamatrix.com)





## Introduction

---

---

### *About this manual*

This manual is written for clinical personnel using the Novamatrix NICO<sup>®</sup> Non-Invasive Cardiac Output Monitor, Model 7300, and the sensors and accessories intended for use with the monitor.

This document contains information which is proprietary and the property of Novamatrix Medical Systems Inc., and may not be reproduced, stored in a retrieval system, translated, transcribed, or transmitted, in any form, or by any means, without the prior explicit written permission of Novamatrix Medical Systems Inc. Novamatrix reserves the right to change specifications without notice.

### *NICO<sup>®</sup> Monitor Technical Description*

Per requirements of IEC 601-1, the NICO<sup>®</sup> monitor is classified as class II equipment, internally powered, with type BF applied part, and an enclosure protection rating of IPX0. The NICO<sup>®</sup> monitor is Year 2000 compliant.

Transport/Storage: -10 to +55° C (14-131° F), 10-95% R.H. non-condensing  
Operating Conditions: 10 to +40° C (50 to 104° F), 10-90% R.H. non-condensing

The NICO<sup>®</sup> monitor, Model 7300, contains no user serviceable parts. Refer servicing to qualified service personnel. A technical Service Manual is available for use by technical personnel.

### *Manufacturing Quality & Safety*

The Novamatrix Medical Systems Inc. manufacturing facility is certified to both ISO 9001 and EN46001 (MDD93/42/EEC Annex II). Novamatrix' products bear the "CE 0086" mark. The product is certified by Underwriter's Laboratories (UL) to bear the UL mark; and tested by TÜV Rheinland to IEC 601-1/EN60601-1.

### *Declaration of Conformity with European Union Directive*

The Authorized Representative for Novamatrix equipment is:

D.R.M. Green  
European Compliance Services Limited,  
Oakdene House,  
Oak Road,  
Watchfield  
Swindon, Wilts SN6 8TD  
United Kingdom

### *Trademarks and Patents*

CAPNOSTAT CO<sub>2</sub> Sensor and NICO are registered trademarks (®); NICO<sub>2</sub> and the stylized NICO<sub>2</sub> with CO<sub>2</sub> shadow, NICO Sensor, NICO Loop and CObar (cardiac output confidence bar), SuperBright and Y-Sensor are trademarks (™) of Novamatrix Medical Systems Inc. Other trademarks and registered trademarks are the property of their respective owners.

NICO<sup>®</sup> and its sensors and accessories are covered by the following USA patents: 4,859,858, 4,859,859, 4,914,720, 5,146,092, 5,153,436, 5,190,058, 5,206,511, 5,251,121, 5,347,843, 5,369,277, 5,379,650, 5,398,680, 5,535,633, 5,616,923, 5,693,944, 5,789,660, 5,793,044, 5,820,550, 5,891,026, 5,999,834, 6,098,622, 6,126,610, D424,193, 6,179,784. Other patents pending.

### *Manual Revision History*

11-Mar-99 Release at Rev. 00.  
23-Mar-99 Revision 01.  
05-Oct-99 Revision 02. R-N677  
20-Mar-00 Revision 03. R-N741  
19-Oct-00 Revision 04, R-N807  
01-Feb-01 Revision 05, R-N850





# Table of Contents

---

---

<b>Welcome to NICO®</b> .....	1
General Description .....	1
Front Panel .....	1
Rear Panel .....	2
Symbols .....	3
AC/Battery Operation .....	4
NICO® Parameter List .....	5
Principles of Operation .....	8
<b>Navigating in Cardiac Output mode</b> .....	10
Areas of the Display .....	10
Navigating the Display System .....	10
KNOB selectable Monitoring Screens .....	11
MENU key Screen Displays .....	13
DATA ENTRY key Screen Displays .....	14
<b>Navigating in Respiratory Mechanics mode</b> .....	16
Areas of the Display .....	16
Navigating the Display System .....	16
KNOB selectable Monitoring Screens .....	17
MENU key Screen Displays .....	18
DATA ENTRY key Screen Displays .....	19
<b>Safety</b> .....	21
<b>Monitoring Cardiac Output</b> .....	24
Preparing for Use .....	24
Begin NICO® monitoring .....	25
Entering Patient Data .....	28
Entering ABG Data .....	29
C.O. Averaging .....	31
<b>Respiratory Monitoring</b> .....	32
Preparing for Use .....	32
Begin Respiratory Monitoring .....	34
Entering Patient Data .....	34
Entering ABG Data .....	36
<b>Monitoring and Setup Screens</b> .....	39
Cardiac Output mode .....	39
Respiratory Screens (Optional) .....	41
Respiratory Mechanics mode .....	46

<b>Notes on Patient Monitoring</b> .....	53
Automatic Purging .....	53
Manual Purging .....	53
Intra-Aortic Balloon Pump .....	54
Configuration Menu .....	55
Reference Handbooks .....	55
<b>Alerts</b> .....	56
NICO <sup>®</sup> Alert Priorities .....	56
Responding to Alert Audio .....	56
Parameter Limit Alerts .....	57
No Respiration Alert .....	58
<b>NICO Sensor™</b> .....	59
Disposable NICO Sensors™ .....	59
<b>CO<sub>2</sub>/Flow Sensors</b> .....	62
Choosing a CO <sub>2</sub> /Flow Sensor .....	63
Connecting a CO <sub>2</sub> /Flow Sensor .....	63
<b>CAPNOSTAT<sup>®</sup> CO<sub>2</sub> Sensor</b> .....	65
The CAPNOSTAT <sup>®</sup> CO <sub>2</sub> Sensor .....	65
CAPNOSTAT <sup>®</sup> Adapter Zero .....	66
<b>Pulse Oximetry Sensors</b> .....	67
Oximetry Sensors .....	67
Finger Sensor .....	68
Y-Sensor™ .....	68
Single Patient Use Sensors .....	71
<b>Messages</b> .....	73
Message Areas - Cardiac Output mode .....	73
Message Areas - Respiratory Mechanics mode .....	78
<b>External Devices</b> .....	82
GE Medical Systems Solar <sup>®</sup> Interface .....	82
Agilent Technologies VueLink Interface .....	83
ASCII Output .....	85
Analog Output .....	86
Alerts and Messages .....	87
<b>Maintenance</b> .....	88
Cleaning and Sterilization .....	88
Monitor Maintenance Schedules .....	89
Battery Maintenance .....	89
<b>Specifications</b> .....	90
<b>NICO<sup>®</sup> Accessories</b> .....	92
<b>Warranty</b> .....	94



## Welcome to NICO®

### General Description

NICO®, a Non-Invasive Cardiac Output monitor from Novamatrix Medical Systems Inc., non-invasively measures and displays cardiac output (C.O.). The NICO® monitor, Model 7300, also displays cardiac index and stroke volume, as well as various respiratory monitoring parameters including CO<sub>2</sub> elimination (VCO<sub>2</sub>) and alveolar minute ventilation. In Respiratory Mechanics mode, NICO® can be used as a respiratory profile monitor, without cardiac output displayed. In either mode, NICO® provides the clinician with important information to aid in precise and efficient patient management.

### Indications


The NICO® monitor is indicated for use by technically skilled clinical personnel. In Cardiac Output mode, the monitor is used for the monitoring of cardiac output and various respiratory parameters of adult patients receiving mechanical ventilation. In Respiratory Mechanics mode, the NICO® monitor is used for monitoring the respiratory parameters of adult, pediatric and neonatal patients. NICO® is not intended for any other purpose.

### Contraindications

In Cardiac Output mode, use of the NICO® monitor is contraindicated in patients in whom a small rise (3-5 mmHg, 0.4-0.67 kPa) in their PaCO<sub>2</sub> level cannot be tolerated.

### Front Panel

The NICO® monitor's front panel includes a display screen, sensor input connectors, a control knob, and operational push button keys and indicators that are explained below.

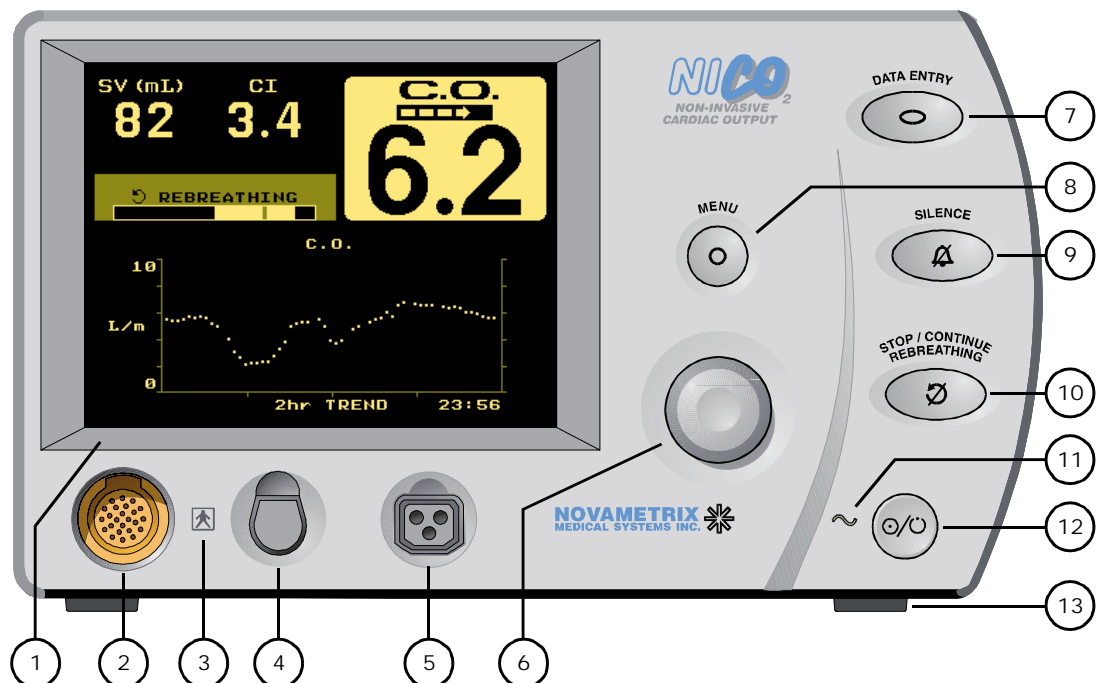


- 1 **Display Screen.** The screen displays NICO® data, respiratory mechanics, trends, waveforms and messages, along with setup and configuration data.
- 2 **CAPNOSTAT® CO<sub>2</sub> Sensor Input Connector.** Connect only a Novametrix CAPNOSTAT® CO<sub>2</sub> Sensor, Catalog Number 9567-00 here.
- 3 **Connector Isolation Icon.** Identifies the connector to either side of this icon as a type BF patient isolation connection.
- 4 **Pulse Oximetry Sensor Input Connector.** Connect only Novametrix pulse oximetry sensors and extension cables approved for use with the NICO® monitor.
- 5 **NICO Sensor™ Input Connector.** Connect only Novametrix NICO Sensors™, Catalog Number 8950-00, 8951-00 and 8952-00 or Novametrix CO<sub>2</sub>/Flow sensors, Catalog Number 9765-00, 9766-00 and 9767-00.
- 6 **KNOB.** The **KNOB** is used to select monitoring screens, scroll through menus and to change or enter values. The **KNOB** is generally turned to access different monitoring screens and to highlight menu options, and pressed to accept or change those selections.
- 7 **DATA ENTRY key.** Press to activate the **DATA ENTRY** screen and illuminate the key's green icon. Press the key again to return to the previously displayed screen. From the **DATA ENTRY** screen, you can enter patient information including height, weight and respiratory gas mixture, and access the **ABG DATA ENTRY** screen.
- 8 **MENU key.** Press to activate the **SELECT A SCREEN** menu and illuminate the key's green icon. Press the key again to return to the previously displayed screen. From the **SELECT A SCREEN** menu you can, by turning the knob, highlight the screen you wish to display. Press the **MENU** key or the **KNOB** to display that selected screen.
- 9 **SILENCE key.** The **SILENCE** key is used to mute/prevent audible alerts. It also visually indicates the presence of a "High Priority Alert". The Silence feature operates in two modes; a temporary "2 Minute Silence" mode and an "Audio Disabled" mode.
  - 2 Minute Silence — Press and release to activate or deactivate the two minute silence. The key's icon illuminates amber when active and audible alerts will be muted for two minutes, after which the icon turns off and any audible alert will sound.
  - Audio Disabled — Press and hold for one second to prevent or allow any audible alerts. The key's icon illuminates and flashes amber to indicate that all audible alerts are being suppressed.
  - High Priority Alerts — (See "NICO® Alert Priorities" on page 56) The **SILENCE** key's icon illuminates and flashes red to indicate High Priority Alert is active. The icon alternately flashes red and amber if the audio is disabled and a High Priority Alert is active.
- 10 **STOP/CONTINUE REBREATHING key.** Press to start NICO® monitoring and the automatic rebreathing process. Subsequent presses will stop (amber indicator illuminated) or continue (amber indicator off) the rebreathing process. Press and hold for two seconds to reset the NICO® algorithm; the C.O. value and averaging filter will be cleared. The **STOP/CONTINUE REBREATHING** key will be amber and inactive in Respiratory Mechanics mode.
- 11 **AC Mains Power Indicator.** This icon illuminates green to indicate AC Mains power is applied to the monitor. To illuminate the icon, the monitor must be plugged into the AC outlet and the monitor's rear panel power switch must be on ("I").
- 12 **OPERATE/STANDBY key.** Press this key to turn the monitor on. If connected to the AC outlet, the monitor uses AC power, otherwise it powers up using its internal battery (provided the battery is charged). Press the **OPERATE/STANDBY** key again to put the monitor into Standby mode (if using AC power) or to turn it off (if using battery power).
- 13 **Kickstand (front and rear).** The NICO® monitor can be positioned for better viewing from above or below by extending the kickstand at the front or rear of the monitor.

## Rear Panel

The NICO® monitor's rear panel includes an AC Mains power input module, three RS232 serial communications ports, an analog input/output port, equipotential connector, fan and ventilation slots, and the monitor's serial number label. These are explained below.



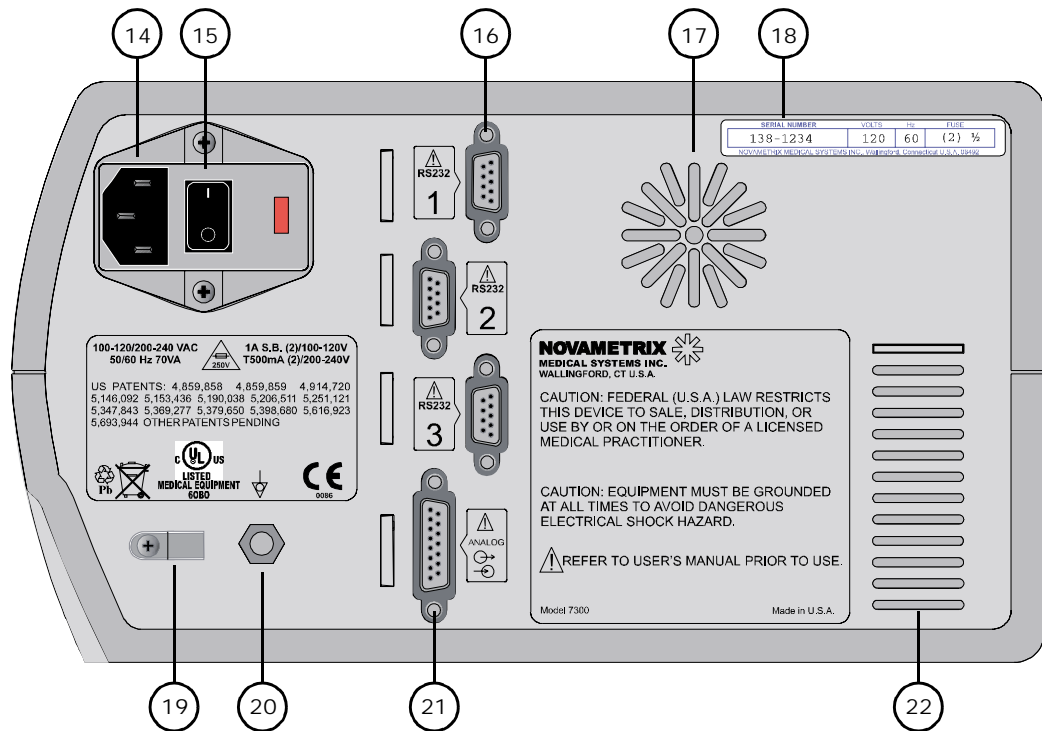

- 14 AC Mains Power Cord Connection.** Connect only approved hospital-grade line cords to this connector.
- 15 AC Mains Power Switch.** This switch controls the flow of AC current into the NICO® monitor. Press the “|” portion of the switch to supply the monitor with AC power, or the “O” portion of the switch to interrupt the flow of AC power. If supplied with AC power, the monitor illuminates the front panel AC Mains Power Indicator, energizes the fan and recharges the internal battery.
- 16 RS232 Communications Ports.** Three 9-pin serial communications ports provide for digital communications with the NICO® monitor. (See “RS232 Communications” ).
- 17 Fan.** The fan draws air in through the monitor. Do not block the fan’s air intake slots.
- 18 Serial Number Label.** The NICO® monitor’s serial number is shown here. Refer servicing to qualified personnel.
- 19 Power cord retaining clip.** If desired, remove the screw, slip the cord through the clip and insert and tighten the screw. Use only the supplied screw to secure the clip.
- 20 Equipotentiality.** Connection to the monitor’s chassis (earth ground system).
- 21 Analog Input/Output Port.** This 15-pin connector provides analog signal output capability for the NICO® monitor (Input reserved for future use).
- 22 Ventilation Slots.** Do not block the air ventilation slots.

## Symbols

These symbols may be found on the NICO® monitor, its sensors, accessories and documentation.



**Attention**  
Consult manual for detailed information.



**Patient Isolation**  
Identifies the patient isolation connection as type BF.



**Single Patient Use**  
Treat in accordance with protocol for “single patient use” items.



**AC Mains Power Switch**  
“|” ON-connection to mains;  
“O” OFF-disconnection from mains



**Mains Fuse Rating**  
Mains rating for replacement fuses



**Separate collection**  
Take appropriate steps to ensure that spent batteries are collected separately when disposed of. This symbol is found on the internal battery and the monitor enclosure.



**Recyclable item**  
This symbol is found on the internal battery and the monitor enclosure.



**Equipotentiality**  
Connection to monitor's chassis.



**Heavy Metal Content**  
Indicates heavy metal content, specifically lead. This symbol is found on the internal battery and the monitor enclosure.


## AC/Battery Operation

The NICO® monitor is designed to be operated from AC Mains power. An internal battery provides uninterrupted monitoring and trending capability for short periods (no more than 45 minutes) if the AC power is removed.

### AC Mains Operation

To operate NICO® from AC Mains power:

- 1 Plug the line cord into the rear panel connector and the AC Mains power outlet.
- 2 Set the rear panel power switch to the On "I" position.
  - The front panel AC Mains Power indicator illuminates.
  - The monitor's fan turns on.
  - The internal battery starts to recharge.

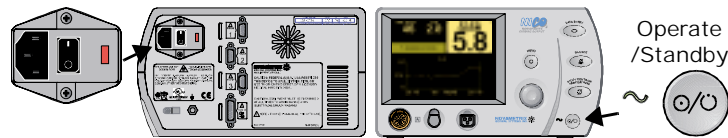

- 3 Press the front panel **Operate/Standby** key to turn the monitor on and off.

### Battery Operation



The NICO® monitor automatically switches to battery power from AC Mains power if the AC source is removed. A fully charged battery will power the monitor for up to 45 minutes. While on battery power, NICO® displays a battery icon that "drains" as power is consumed.

The battery icon starts to flash when approximately 5 minutes of battery power remain. An audible alert tone also sounds.

Reconnect to AC Mains power or the monitor will automatically shut off. A depleted battery may require 12-16 hours to fully recharge.

To operate NICO® on battery power:


- 1 Unplug the line cord or set the rear panel power switch to the Off "O" position.
  - The front panel AC Mains Power indicator turns off.
- 2 Press the front panel **Operate/Standby** key to turn the monitor on and off.

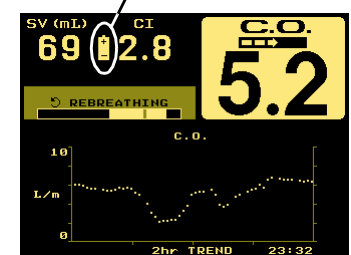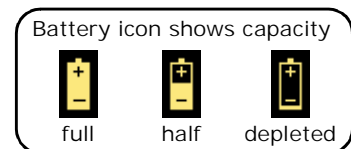

## NICO® Parameter List


**Cardiac Output mode** The NICO® monitor displays the parameters described in this table.

Label	Parameter	Range/Units	Description	Screen Display
C.O.	Cardiac Output	0.5-19.9 L/min	Volume of blood pumped by the heart each minute	All
CO-a	Average Cardiac Output	0.5-19.9 L/min	C.O. averaged value, displayed when the fast-mode cardiac output mode is chosen for large display.	Last Completed Cycle rebreathing curve & Tabular Data
CO-f	Fast-mode Cardiac Output	0.5-19.9 L/min	C.O. unaveraged value, displayed when the average cardiac output mode is chosen for large display.	Last Completed Cycle rebreathing curve & Tabular Data
Cdyn	Dynamic Compliance	0-500 ml/cmH <sub>2</sub> O	Volume the lungs expand for a given pressure  Note that if the ventilator is set for an inspiratory pause that is detected by NICO, Cdyn becomes Cstat.	Numerics & Tabular Data
CI	Cardiac Index	0-9.9 L/min/m <sup>2</sup>	C.O. divided by body surface area	All
ETCO <sub>2</sub>	End Tidal Carbon Dioxide	0-150 mmHg 0-20.0 % 0-20.0 kPa	Maximum CO <sub>2</sub> plateau value at the end of the breath (reflects alveolar CO <sub>2</sub> )	Numerics, Respiratory Numerics, CO <sub>2</sub> /SpO <sub>2</sub> , SBCO <sub>2</sub> , Tabular Data & rebreathing curves
Insp CO <sub>2</sub>	Inspired Carbon Dioxide	3-50 mmHg 0.4-6.7 % 0.4-6.7 kPa	Maximum CO <sub>2</sub> value observed during the baseline portion of the Inspiratory phase of the breath (baseline shift above zero point)	General Message Area (if above 3 mmHg for 10 sec (0.4 % or kPa))
MAP	Mean Airway Pressure	0-100 cmH <sub>2</sub> O	Mean (average) pressure in the airway throughout the breath	Numerics & Tabular Data
MV	Minute Volume	2-40 L/min Adult	Volume (in liters) of gas delivered to the patient per minute	Numerics, Respiratory Numerics & Tabular Data
MV alv	Alveolar Minute Volume	0.05-16 L/min	MV less deadspace (wasted) ventilation	Numerics, Respiratory Numerics & Tabular Data
PCBF	Pulmonary Capillary Blood Flow	0.5-19.9 L/min	Portion of the cardiac output that is effective in gas exchange	3 min Cycle in Progress rebreathing curve & Tabular Data
PeCO <sub>2</sub> /FeCO <sub>2</sub>	Mixed expired CO <sub>2</sub>	0-100 mmHg, 0-13.2 kPa or %	Volume weighted average CO <sub>2</sub> in the breath	Respiratory Numerics & Tabular Data
PEEP	Positive End Expiratory Pressure	0-99 cmH <sub>2</sub> O	Pressure in the lungs at the end of expiration	Numerics, Flow/Pressure & Tabular Data
PEF	Peak Expiratory Flow	2-180 L/min	Highest absolute flow rate during expiration	Tabular Data
PIF	Peak Inspiratory Flow	2-180 L/min	Highest absolute flow rate during inspiration	Tabular Data
PIP	Peak Inspiratory Pressure	0-120 cmH <sub>2</sub> O	Peak (highest) pressure in the airway during inspiration	Numerics, Flow/Pressure & Tabular Data

Label	Parameter	Range/Units	Description	Screen Display
	Pulse Rate	30-250 bpm	Number of pulse beats per minute	Numerics, CO <sub>2</sub> /SpO <sub>2</sub> & Tabular Data
Raw	Airway Resistance	0-100 cmH <sub>2</sub> O/L/sec	Pressure required to cause gas flow at a given rate	Numerics & Tabular Data
RR	Respiration Rate	2-120 br/min	Number of breaths per minute	Numerics, CO <sub>2</sub> /SpO <sub>2</sub> & Tabular Data
RSBI	Rapid Shallow Breathing Index	0-1000 br/min/L	Respiratory rate divided by average spontaneous tidal volume (only calculated when RR < 57)	Respiratory Numerics & Tabular Data
SpO <sub>2</sub>	Oxygen Saturation	0-100 %	Oxyhemoglobin as a percentage of total hemoglobin less dysfunctional hemoglobin	Numerics, CO <sub>2</sub> /SpO <sub>2</sub> , Tabular Data & rebreathing curves
SV	Stroke Volume	0-250 ml	Volume of blood pumped by the heart each beat	All
SVI	Stroke Volume Index	0-125 ml	Stroke volume divided by body surface area	Tabular Data
SVR	Systemic Vascular Resistance	0-5000 dynes sec/cm <sup>5</sup>	Resistance exerted by the blood vessels on blood flow and is an indicator of left ventricular afterload.	SVR Calculation & Tabular Data
SVRI	Systemic Vascular Resistance Index	0-5000 dynes sec/cm <sup>5</sup>	SVR normalized to body surface area	Tabular Data
VCO <sub>2</sub>	Carbon Dioxide Elimination	0-3000 ml/min	Volume of CO <sub>2</sub> eliminated through the breath each minute	Numerics, Respiratory Numerics, SBCO <sub>2</sub> , Tabular Data & rebreathing curves
Vd Aw	Airway deadspace	0-500 ml	Includes added mechanical deadspace proximal to the flow sensor	Respiratory Numerics, SBCO <sub>2</sub> , Tabular Data
Vd/Vt	Deadspace to tidal volume ratio	0-1.00 ml	(PaCO <sub>2</sub> -PeCO <sub>2</sub> )/PaCO <sub>2</sub>	Respiratory Numerics & ABG Data Entry
Vd alv	Alveolar deadspace	0-500 ml	Difference between physiologic and airway deadspace	Respiratory Numerics & ABG Data Entry
Vt alv	Alveolar tidal volume	0-2400 ml	Tidal volume less airway deadspace	Respiratory Numerics, SBCO <sub>2</sub> & Tabular Data
Vte	Expired Tidal Volume	200-3000 ml	Volume of gas exhaled per breath	Respiratory Numerics, Flow/Pressure & SBCO <sub>2</sub>
Vte-m	Expired Tidal Volume - mechanical	200-3000 ml	Volume of mechanically exhaled gas, per breath	Tabular Data
Vte-s	Expired Tidal Volume - spontaneous	200-3000 ml	Volume of spontaneously exhaled gas, per breath	Tabular Data
Vti	Inspired Tidal Volume	200-3000 ml	Volume of gas inhaled per breath	Respiratory Numerics, Flow/Pressure & SBCO <sub>2</sub>
Vti-m	Inspired Tidal Volume - mechanical	200-3000 ml	Volume of mechanically inhaled gas, per breath	Tabular Data
Vti-s	Inspired Tidal Volume - spontaneous	200-3000 ml	Volume of spontaneously inhaled gas, per breath	Tabular Data

**Respiratory  
Mechanics mode**

The NICO® monitor displays the parameters described in this table.

Label	Parameter	Range/Units	Description	Screen Display
Cdyn	Dynamic Compliance	0-500 ml/cmH <sub>2</sub> O	Volume the lungs expand for a given pressure	Numerics & Tabular Data
ETCO <sub>2</sub>	End Tidal Carbon Dioxide	0-150 mmHg 0-20.0 % 0-20.0 kPa	Maximum CO <sub>2</sub> plateau value at the end of the breath (reflects alveolar CO <sub>2</sub> )	All
I:E	I:E Ratio	1:9.9 or 4:1	Ratio of inspiratory time (ti) to expiratory time (te)	Flow/Pressure
MAP	Mean Airway Pressure	0-100 cmH <sub>2</sub> O	Mean (average) pressure in the airway throughout the breath	Numerics & Tabular Data
MV	Minute Volume	0.4-40 L/min adult 0.06-30 L/min pedi. 0.01-5 L/min neonatal	Volume (in liters) of gas delivered to the patient per minute	All except Flow/Pressure & Loops
MV alv	Alveolar Minute Volume	0-16 L/min adult 0-8 L/min pediatric 0-4 L/min neonatal	MV less deadspace (wasted) ventilation	Numerics & Tabular Data
NIP	Negative Inspiratory Pressure	0 to -120 cmH <sub>2</sub> O	Maximum negative pressure during inspiratory cycle	Loops
PeCO <sub>2</sub> / FeCO <sub>2</sub>	Mixed expired CO <sub>2</sub>	0-100 mmHg, 0-13.2 kPa or %	Volume weighted average CO <sub>2</sub> in the breath	Numerics, SBCO <sub>2</sub> & Tabular Data
PEEP	Positive End Expiratory Pressure	0-99 cmH <sub>2</sub> O	Pressure in the lungs at the end of expiration	Numerics, Flow/Pressure & Tabular Data
PEF	Peak Expiratory Flow	2-180 L/min adult 0.5-100 L/min pedi. 0.25-25 L/min neo.	Highest absolute flow rate during expiration	Loops & Tabular Data
PIF	Peak Inspiratory Flow	2-180 L/min adult 0.5-100 L/min pedi. 0.25-25 L/min neo.	Highest absolute flow rate during inspiration	Loops & Tabular Data
PIP	Peak Inspiratory Pressure	0-120 cmH <sub>2</sub> O	Peak (highest) pressure in the airway during inspiration	Numerics, Flow/Pressure & Tabular Data
	Pulse Rate	30-250 bpm	Number of pulse beats per minute	All
Raw	Airway Resistance	0-100 cmH <sub>2</sub> O/L/sec adult/pediatric  0-500 cmH <sub>2</sub> O/L/sec neonatal	Pressure required to cause gas flow at a given rate	Numerics & Tabular Data
RR	Respiration Rate	2-120 br/min adult 2-150 br/min pedi. 10-150 br/min neo.	Number of breaths per minute	All
RSBI	Rapid Shallow Breathing Index	0-1000 br/min/L (adult only)	Respiratory rate divided by average spontaneous tidal volume (only calculated when RR < 57)	Loops
SpO <sub>2</sub>	Oxygen Saturation	0-100 %	Oxyhemoglobin as a percentage of total hemoglobin less dysfunctional hemoglobin	All

Label	Parameter	Range/Units	Description	Screen Display
VCO <sub>2</sub>	Carbon Dioxide Elimination	1-3000 ml/min adult/pediatric  0-300 ml/min neonatal	Volume of CO <sub>2</sub> eliminated through the breath each minute	All except Flow/Pressure & Loops
Vd Aw	Airway deadspace	0-500 ml	Includes added mechanical deadspace proximal to the flow sensor	Numerics, SBCO <sub>2</sub> , Tabular Data
Vd/Vt	Deadspace to tidal volume ratio	0-1.00 ml	(PaCO <sub>2</sub> -PeCO <sub>2</sub> ) / PaCO <sub>2</sub>	Numerics, SBCO <sub>2</sub> & ABG Data Entry
Vd alv	Alveolar deadspace	0-500 ml	Difference between physiologic and airway deadspace	Numerics, SBCO <sub>2</sub> & ABG Data Entry
Vt alv	Alveolar tidal volume	0-2400 ml adult 0-1200 ml pediatric 0-160 ml neonatal	Tidal volume less airway deadspace	Numerics, SBCO <sub>2</sub> & Tabular Data
Vte	Expired Tidal Volume	200-3000 ml adult 30-400 ml pediatric 1-100 ml neonatal	Volume of gas exhaled per breath	All except Flow/Pressure & Loops
Vte-m	Expired Tidal Volume - mechanical	200-3000 ml adult 30-400 ml pediatric 1-100 ml neonatal	Volume of mechanically exhaled gas, per breath	Tabular Data
Vte-s	Expired Tidal Volume - spontaneous	200-3000 ml adult 30-400 ml pediatric 1-100 ml neonatal	Volume of spontaneously exhaled gas, per breath	Tabular Data
Vti	Inspired Tidal Volume	200-3000 ml adult 30-400 ml pediatric 1-100 ml neonatal	Volume of gas inhaled per breath	All except Flow/Pressure & Loops
Vti-m	Inspired Tidal Volume - mechanical	200-3000 ml adult 30-400 ml pediatric 1-100 ml neonatal	Volume of mechanically inhaled gas, per breath	Tabular Data
Vti-s	Inspired Tidal Volume - spontaneous	200-3000 ml adult 30-400 ml pediatric 1-100 ml neonatal	Volume of spontaneously inhaled gas, per breath	Tabular Data

## Principles of Operation

### Non-Invasive Cardiac Output (NICO)

NICO<sup>®</sup> calculates cardiac output (C.O.) non-invasively based on respiratory gas analysis, using a technique known as “differential Fick partial rebreathing.” The key to this technique is a NICO Sensor<sup>™</sup>, consisting of a rebreathing valve and a combined CO<sub>2</sub>/Flow sensor placed in the breathing circuit. The NICO Sensor<sup>™</sup> is placed into the ventilator circuit between the patient elbow and ventilator wye. The rebreathing valve is automatically controlled by the monitor. When the valve is activated, the flow of the inspired and expired gas is diverted through a rebreathing NICO Loop.<sup>™</sup> When the valve is deactivated, this additional rebreathing volume is bypassed and normal ventilation resumes. Every three minutes, a baseline, rebreathing and stabilization phase occurs. (See “The NICO<sup>®</sup> Cycle” on page 27.) A non-invasive cardiac output calculation is made following the end of each three minute cycle. The calculation is based on the changes induced in CO<sub>2</sub> elimination and end tidal CO<sub>2</sub> in response to the rebreathing volume. The increase in end tidal CO<sub>2</sub>, which reflects the increase in PaCO<sub>2</sub>, is usually 3-5 mmHg (0.4-0.67 kPa) and returns to baseline in less than 30 seconds.

The Fick equation using CO<sub>2</sub> as an indicator states that cardiac output is equal to CO<sub>2</sub> elimination divided by the venous-arterial difference in the CO<sub>2</sub> content:  $VCO_2 / (CvCO_2 - CaCO_2)$ . The partial rebreathing method yields a differential form of the Fick equation, eliminating the need to measure mixed venous CO<sub>2</sub> (assumed constant during the rebreathing period and therefore cancels out of the equation). This indirect Fick method is then corrected for shunt,

based on the Nunn's iso-shunt curves using SpO<sub>2</sub> (or entered PaO<sub>2</sub>) and a user-entered value for FiO<sub>2</sub> (INSP O<sub>2</sub>).

### Carbon Dioxide (CO<sub>2</sub>)

NICO® uses the CAPNOSTAT® CO<sub>2</sub> Sensor to measure CO<sub>2</sub> by using the infrared absorption technique. The principle is based on the fact that CO<sub>2</sub> molecules absorb infrared (IR) light energy of specific wavelengths, with the amount of energy absorbed being directly related to the CO<sub>2</sub> concentration. When an IR beam is passed through a gas sample containing CO<sub>2</sub>, the electronic signal from the photodetector (which measures the remaining light energy) can be obtained. This signal is then compared to the energy of the IR source and calibrated to accurately reflect CO<sub>2</sub> concentration in the sample. To calibrate, the photodetector's response to a known concentration of CO<sub>2</sub> is stored at the factory in the monitor's memory. A reference channel accounts for optical changes in the sensor, allowing the system to remain in calibration without user intervention.

### Flow and Pressure

Flow and pressure measurements in the NICO® monitor are made by a fixed orifice differential pressure pneumotachometer. Respirated gas flowing through the flow sensor causes a small pressure drop across the two tubes connected to the sensor. This pressure drop is transmitted through the tubing to a differential pressure transducer located inside the monitor, and is correlated to flow according to the factory stored calibration. User calibration is not required due to the ability of the plastic injection mold to repeatedly produce precision flow sensors. The pressure transducer is automatically "zeroed" to correct for changes in ambient temperature and electronics. The NICO® monitor system software compensations allow accurate flow and volume measurements in the presence of high oxygen concentrations, anesthetic gases and helium-oxygen mixtures. When compensated, gas density and viscosity effects do not cause significant errors in flow measurement.

### Carbon Dioxide Elimination (VCO<sub>2</sub>)

Carbon dioxide elimination (VCO<sub>2</sub>) is a key measurement for NICO® calculations. It is calculated based on a mathematical integration of the measured flow and CO<sub>2</sub> signals. These signals are obtained from practically the same point at the patient's airway, thereby insuring optimal accuracy. Both the flow and CO<sub>2</sub> sensors are integral components of the NICO Sensor™.

### Oxygen Saturation (SpO<sub>2</sub>) & Pulse Rate

Oxygen saturation (SpO<sub>2</sub>) is used by the NICO® monitor to calculate the shunt correction of the NICO® calculation, and the pulse rate is used to calculate stroke volume.

SpO<sub>2</sub> is determined using sensors containing red and infrared light emitting diodes (LEDs). The light from each LED is beamed through a pulsating vascular bed such as the patient's finger or toe. The remaining light not absorbed by the tissue reaches a photodiode light receptor in the sensor. Oxygen saturated blood absorbs different amounts of light at each wavelength as compared to unsaturated blood. Therefore, the amount of light absorbed by the blood in each pulse can be used to calculate oxygen saturation.

NICO® is calibrated to display "functional" saturation. This differs from the "fractional" saturation value displayed by most co-oximeters. Functional saturation represents the amount of oxyhemoglobin as a percentage of the hemoglobin that can be oxygenated. Dysfunctional hemoglobins, (COHb and METHb) are not included in the measurement of functional saturation.

- Functional Saturation =  $HbO_2 / 100 - (COHb + METHb)$ ; HbO<sub>2</sub> is fractional hemoglobin, COHb is carboxyhemoglobin, and METHb is methemoglobin.

Pulse Rate, derived from the pulse oximetry sensor, is calculated by measuring the time interval between the peaks of the infrared light waveform. The inverse of this measurement is displayed as pulse rate.



## Navigating in Cardiac Output mode

### Areas of the Display


The major sections of the Cardiac Output mode screen are identified below.

Stroke Volume (SV) and Cardiac Index (CI) are displayed in all views. A battery icon also appears if on battery power.

Cardiac Output (C.O.) value is displayed in all views. The CObar™ confidence indicator is replaced by **FAST MODE** when in fast mode.

A General Message area for status, alert and error information. Displayed in all views (here shown blank).


A Cardiac Output Message area for C.O. information is displayed in all views. The Rebreathing Bar is displayed during the rebreathing portion of the cycle.

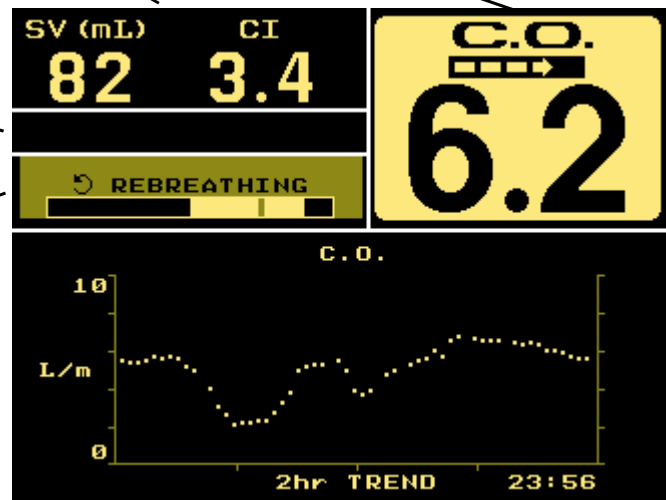

The lower half of the display presents trend, waveform, respiratory and numeric data to the user. Various data entry, setup and alert menus are also presented here. Use the **KNOB** and the **MENU** and **DATA ENTRY** keys to select the various displays.

### Navigating the Display System


Use the **KNOB**, **MENU**, and **DATA ENTRY** keys to navigate the NICO® display system (as outlined in the following sections).

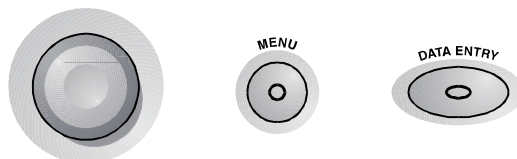



## KNOB selectable Monitoring Screens


The **KNOB** is used to page through monitoring screens, scroll through menus and make selections, and to change or enter values. The **KNOB** is generally turned to access different monitoring screens and to highlight menu options, and pressed to accept or change selections.

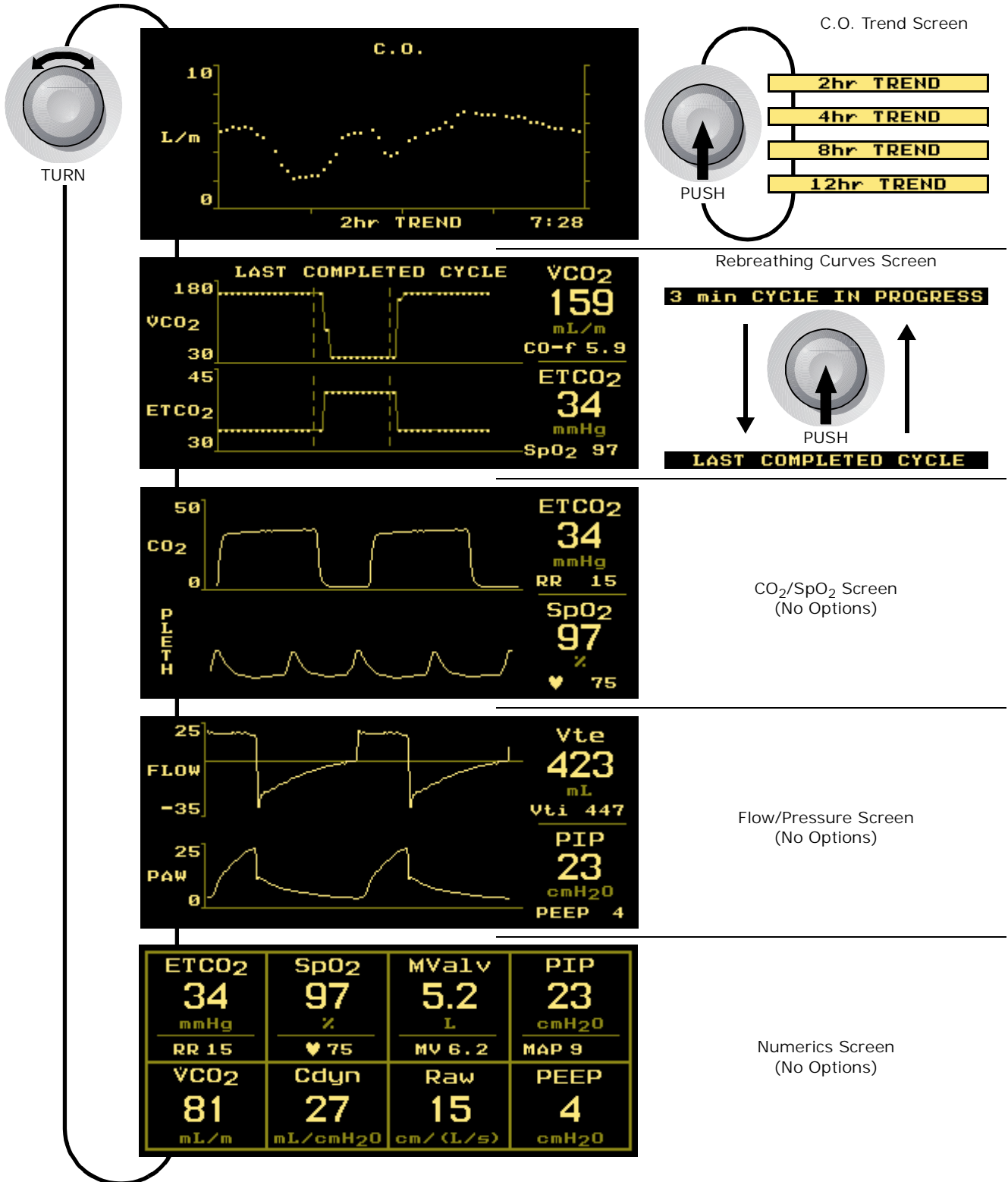


**KNOB Selectable  
Respiratory Screens**

The following Respiratory Screens are available in the monitoring mode **only** when activated. Press the **MENU** key and select **RESP SCREENS** from the **SELECT A SCREEN** menu by turning and then pressing the **KNOB**.

From the **ADD RESPIRATORY SCREENS** menu, highlight and then select which screens will appear in the base monitoring mode by turning and then pressing the **KNOB**. When enabled, selected screens can be displayed by turning the **KNOB** while viewing any monitoring screen.

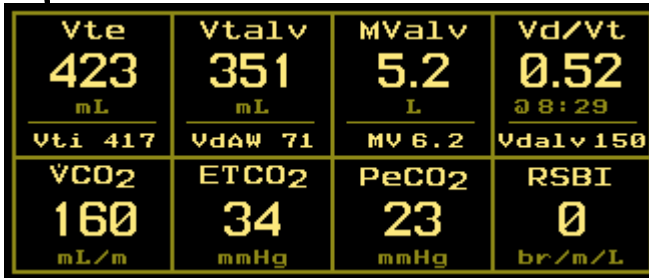


Respiratory Numerics

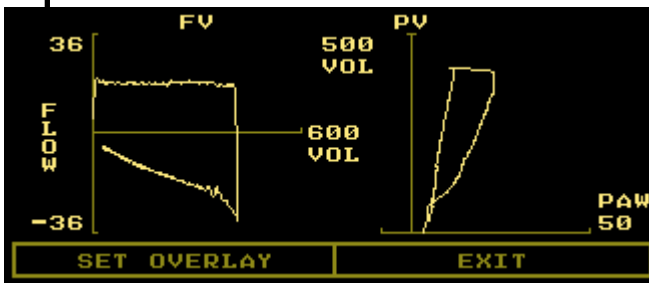

Flow Volume Loop and  
Pressure Volume Loop




Single Breath CO<sub>2</sub> Waveform

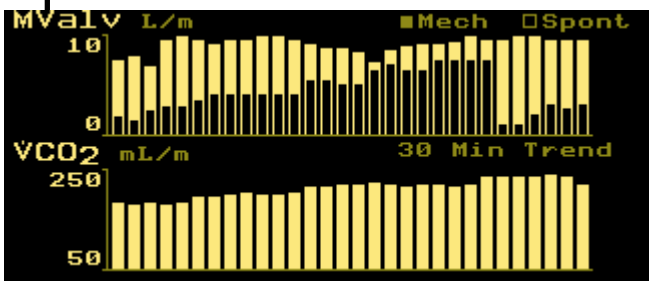

VCO<sub>2</sub>/MValv Trend


## MENU key Screen Displays

Press the **MENU** key to activate the **SELECT A SCREEN** menu and illuminate the key's green icon. Press the key again to return to the previously displayed screen. From the **SELECT A SCREEN** menu turn the **KNOB** to highlight the screen you wish to display. Press the **MENU** key or the **KNOB** to display that selected screen.

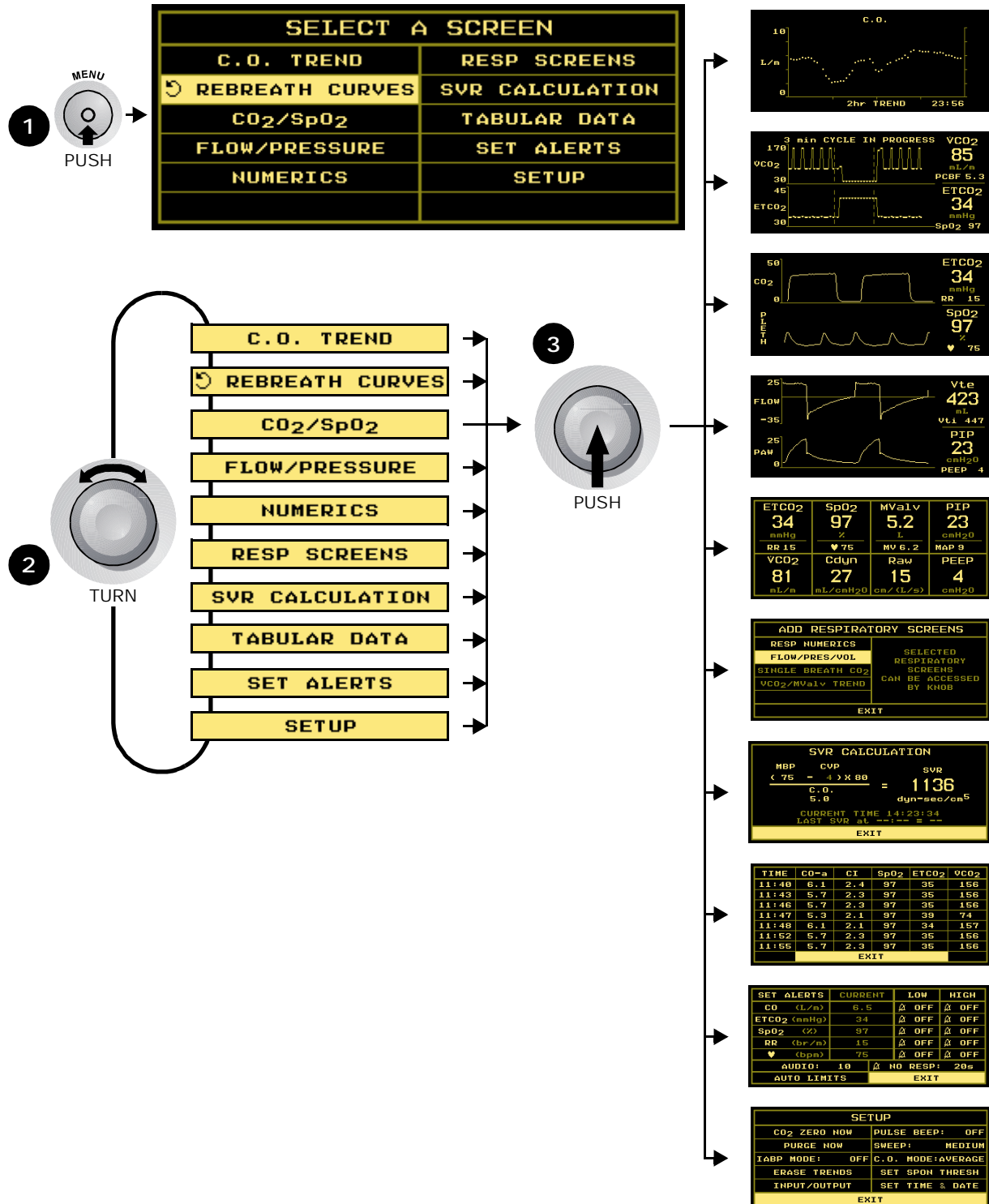


## DATA ENTRY key Screen Displays

Press the **DATA ENTRY** key to activate the **DATA ENTRY** screen and illuminate the key's green icon. Press the key again to return to the previously displayed screen. From the **DATA ENTRY** screen, you can enter patient information including height, weight and respiratory gas mixture, and access the **ABG DATA ENTRY** screen. (See "Entering Patient Data" on page 28 for details.)

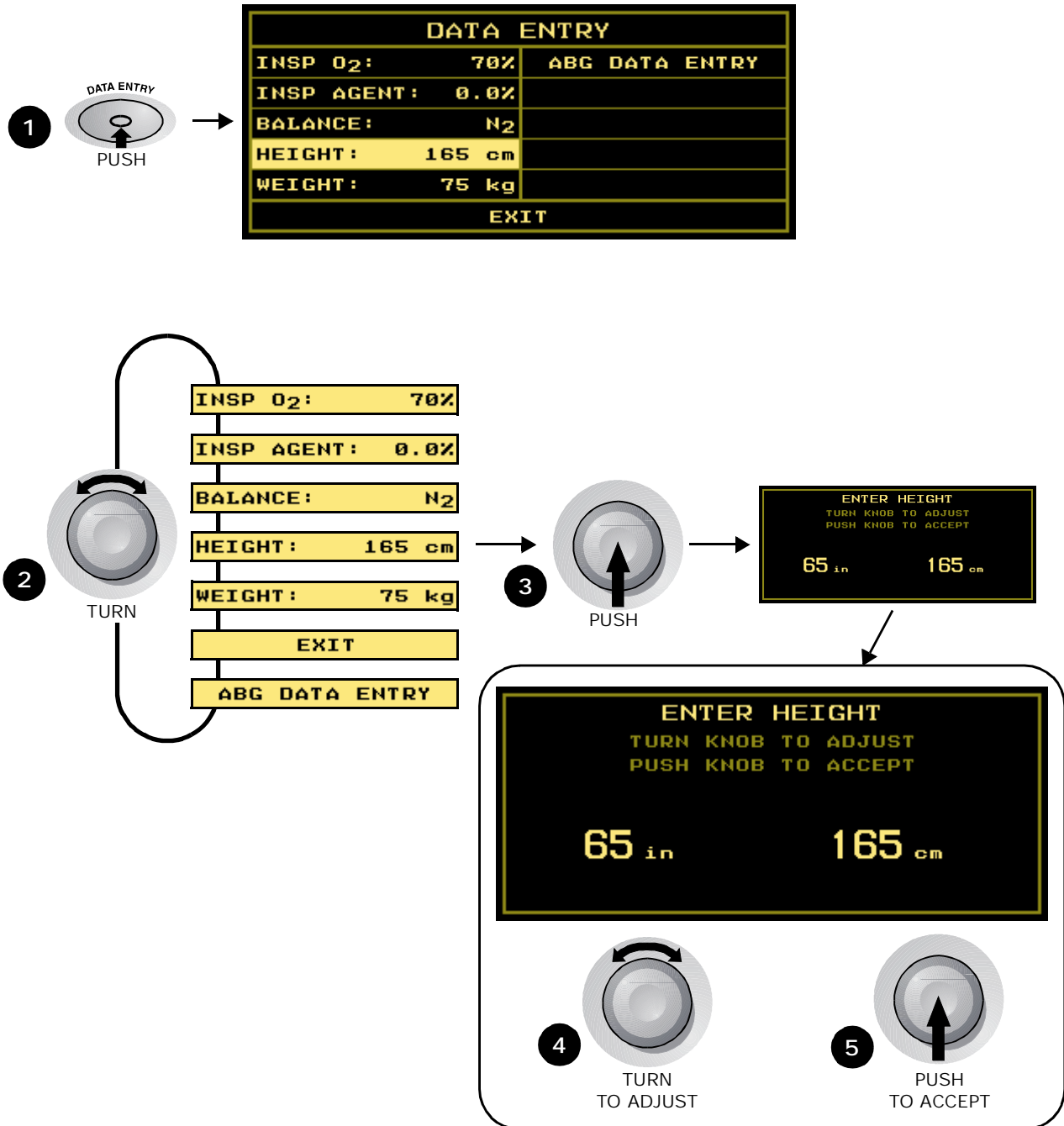

**ABG Data Entry Screens**

From the **DATA ENTRY** screen, select **ABG DATA ENTRY**. Turn and press the **KNOB** to enter PaCO<sub>2</sub>, PaO<sub>2</sub>, and Hemoglobin entry screens. (See "Entering Patient Data" on page 28 for details.)

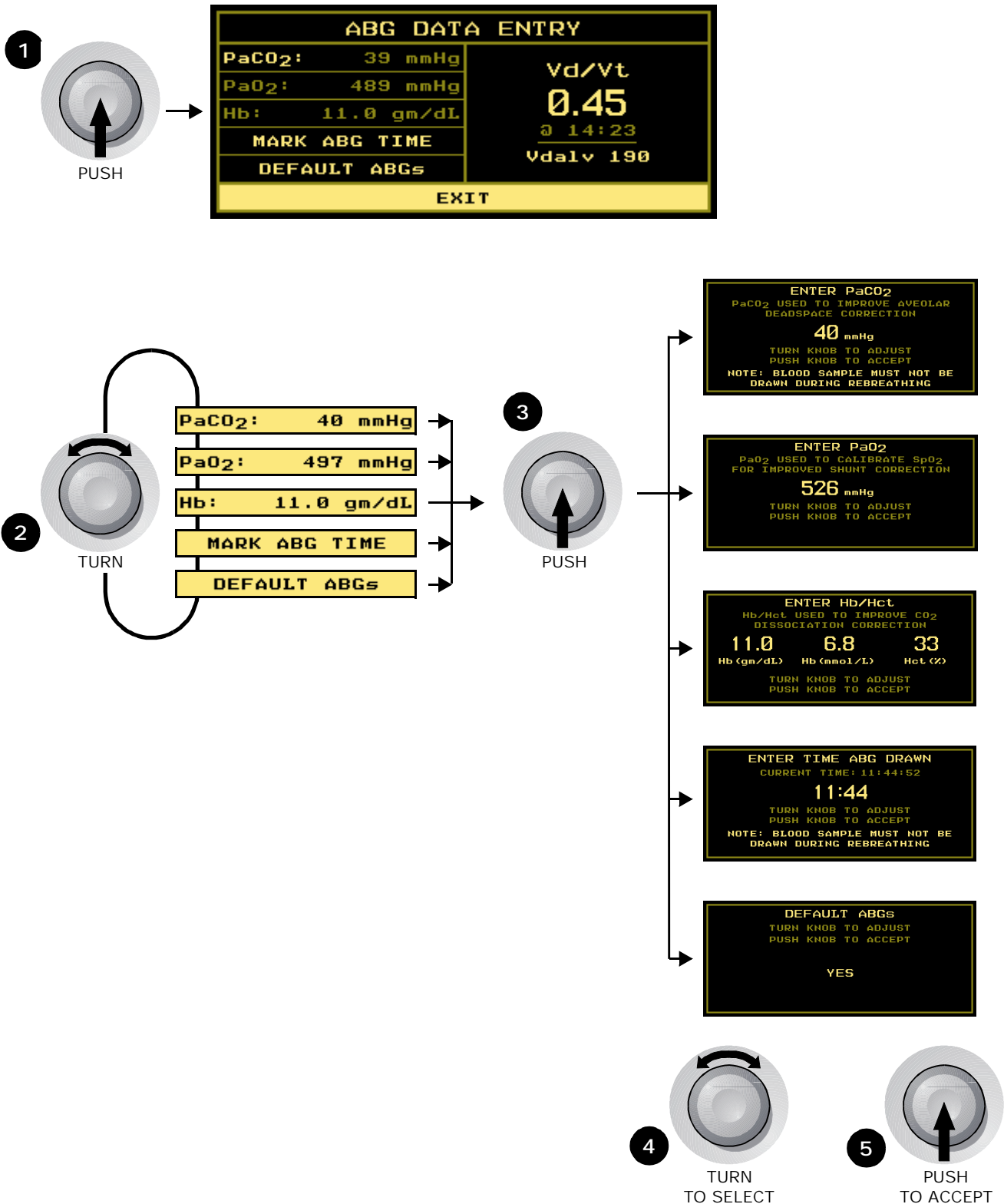




## Navigating in Respiratory Mechanics mode


### Areas of the Display

The major sections of the Respiratory Mechanics mode screen are identified below.

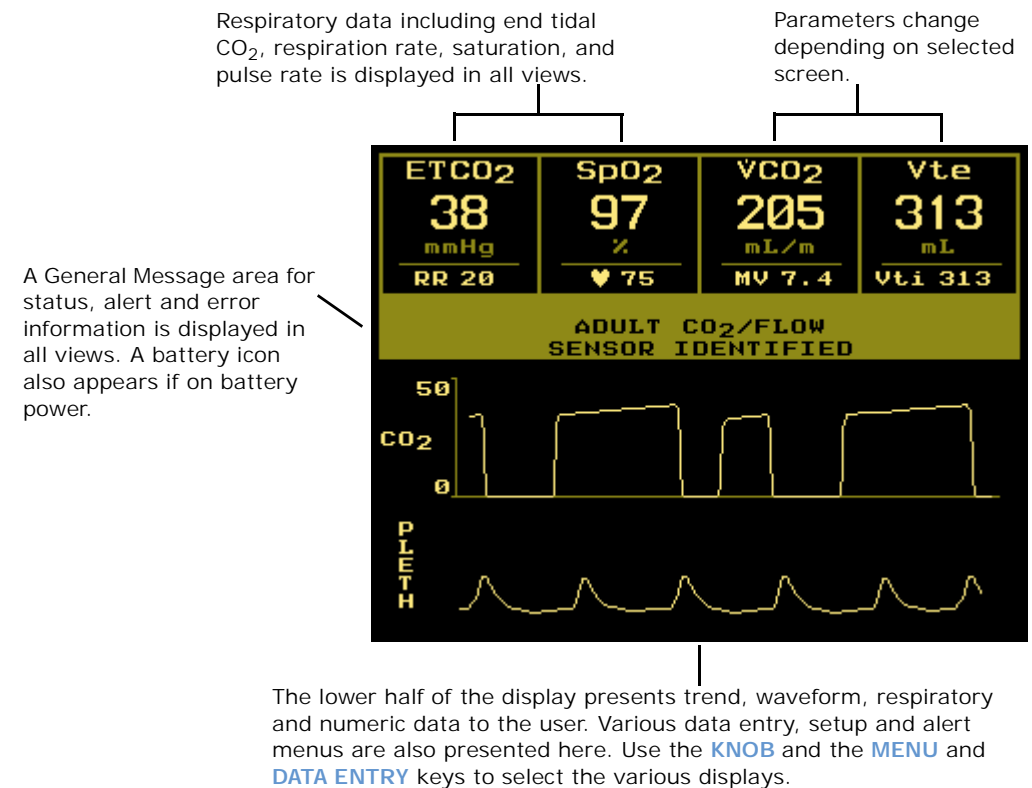


### Navigating the Display System

Use the **KNOB**, **MENU**, and **DATA ENTRY** keys to navigate the NICO® display system (as outlined in the following sections).

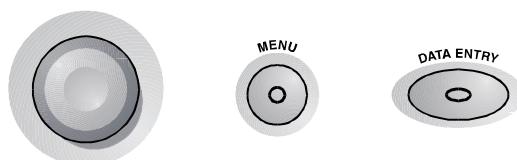


## KNOB selectable Monitoring Screens

The **KNOB** is used to page through monitoring screens, scroll through menus and make selections, and to change or enter values. The **KNOB** is generally turned to access different monitoring screens and to highlight menu options, and pressed to accept or change selections.

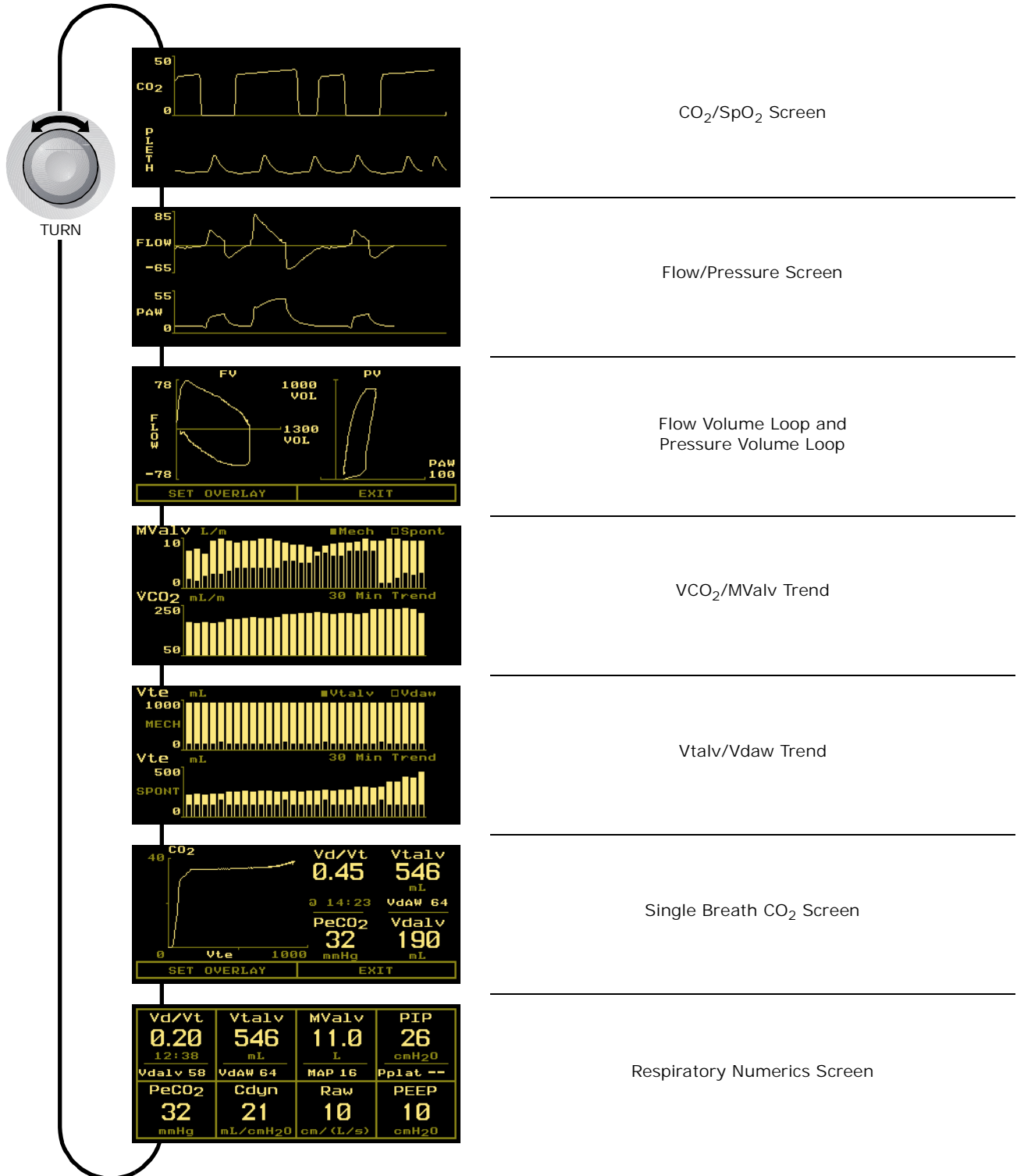

## MENU key Screen Displays


Press the **MENU** key to activate the **SELECT A SCREEN** menu and illuminate the key's green icon. Press the key again to return to the previously displayed screen. From the **SELECT A SCREEN** menu turn the **KNOB** to highlight the screen you wish to display. Press the **MENU** key or the **KNOB** to display that selected screen.

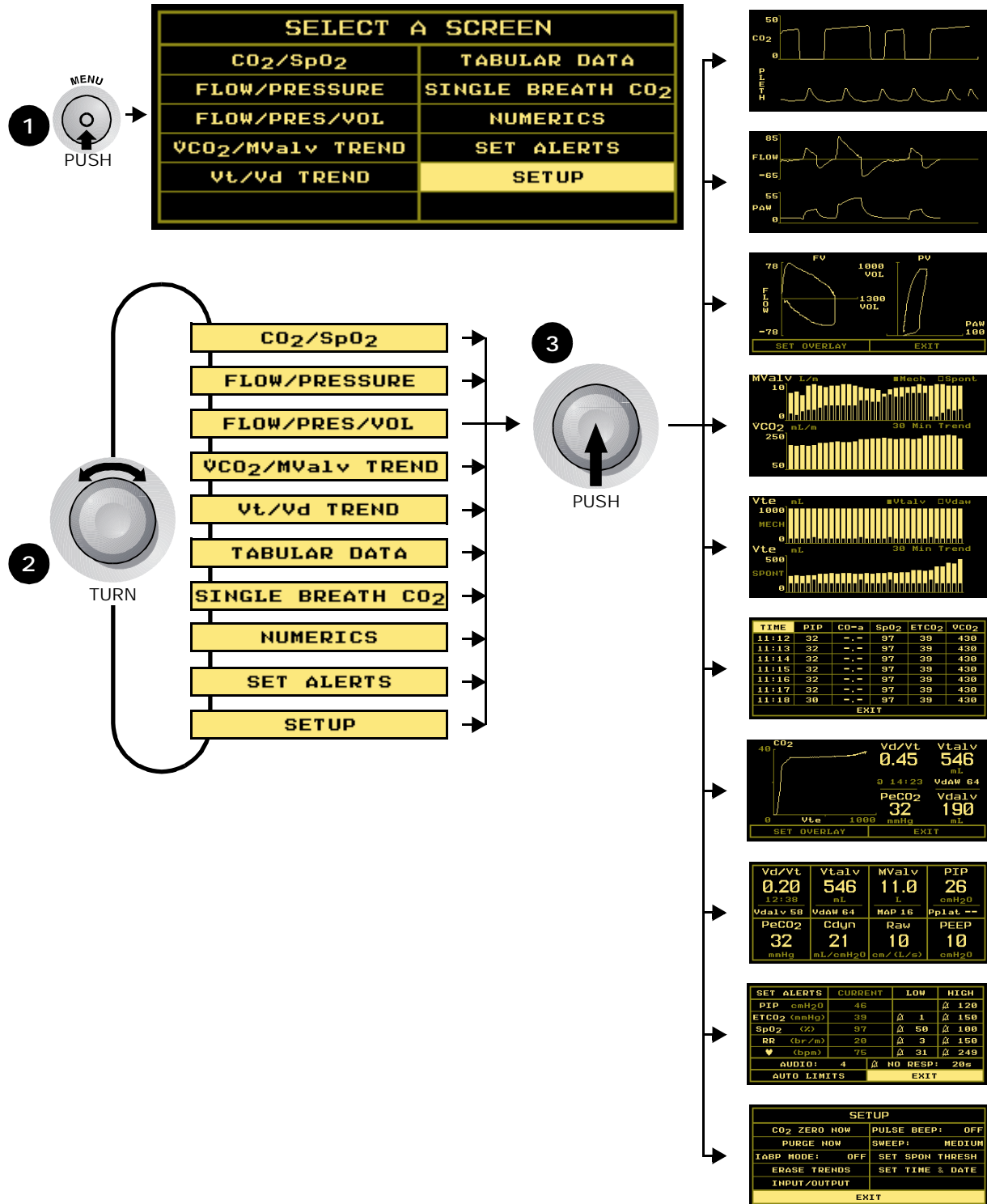




## DATA ENTRY key Screen Displays

Press the **DATA ENTRY** key to activate the **DATA ENTRY** screen and illuminate the key's green icon. Press the key again to return to the previously displayed screen. From the **DATA ENTRY** screen, you can enter patient information including height, weight and respiratory gas mixture, and access the **ABG DATA ENTRY** screen. (See "Entering Patient Data" on page 28 for details.)

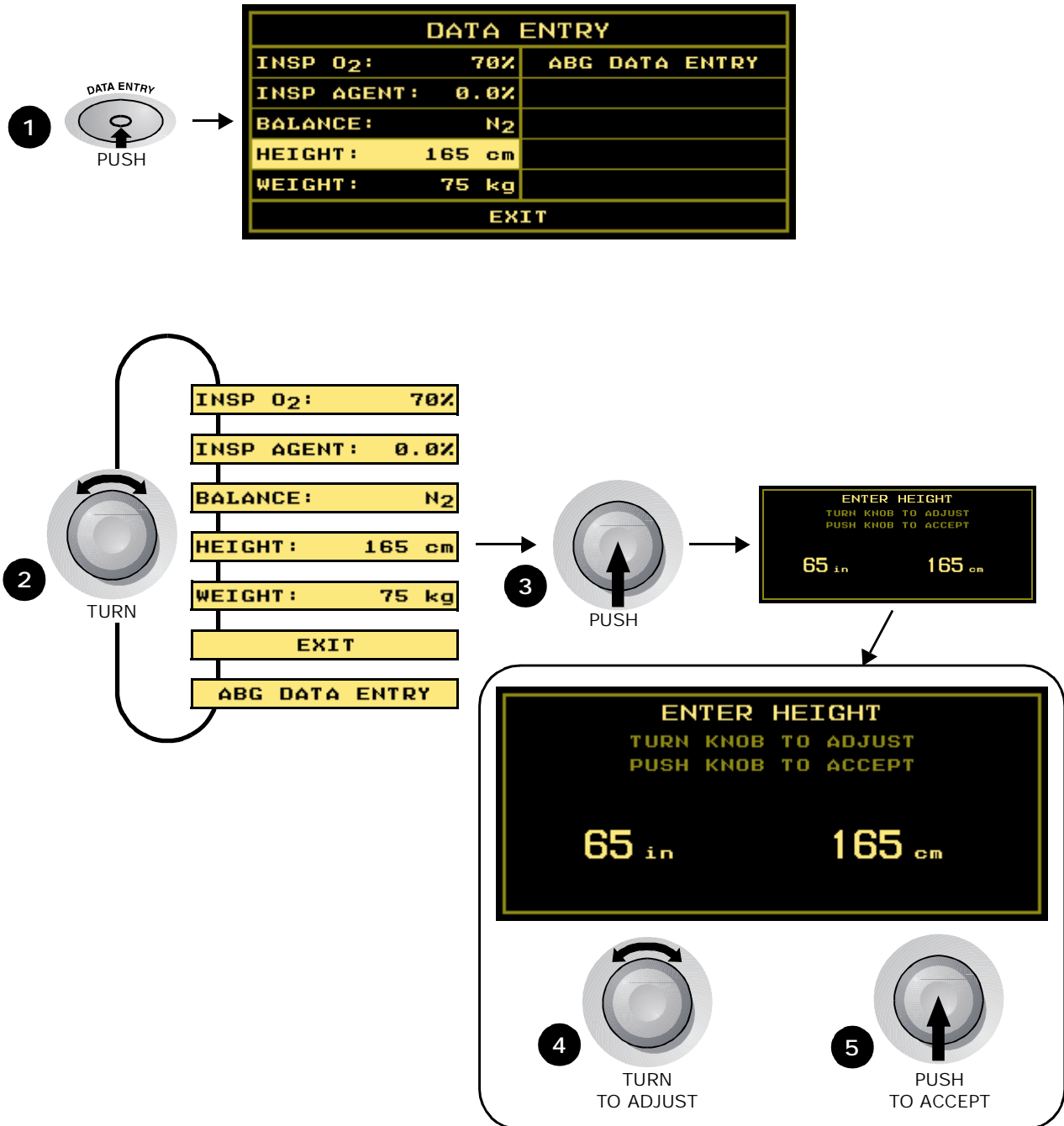

**ABG Data Entry Screens**

From the **DATA ENTRY** screen, select **ABG DATA ENTRY**. Turn and press the **KNOB** to enter PaCO<sub>2</sub>, PaO<sub>2</sub>, and Hemoglobin entry screens. (See "Entering Patient Data" on page 28 for details.)

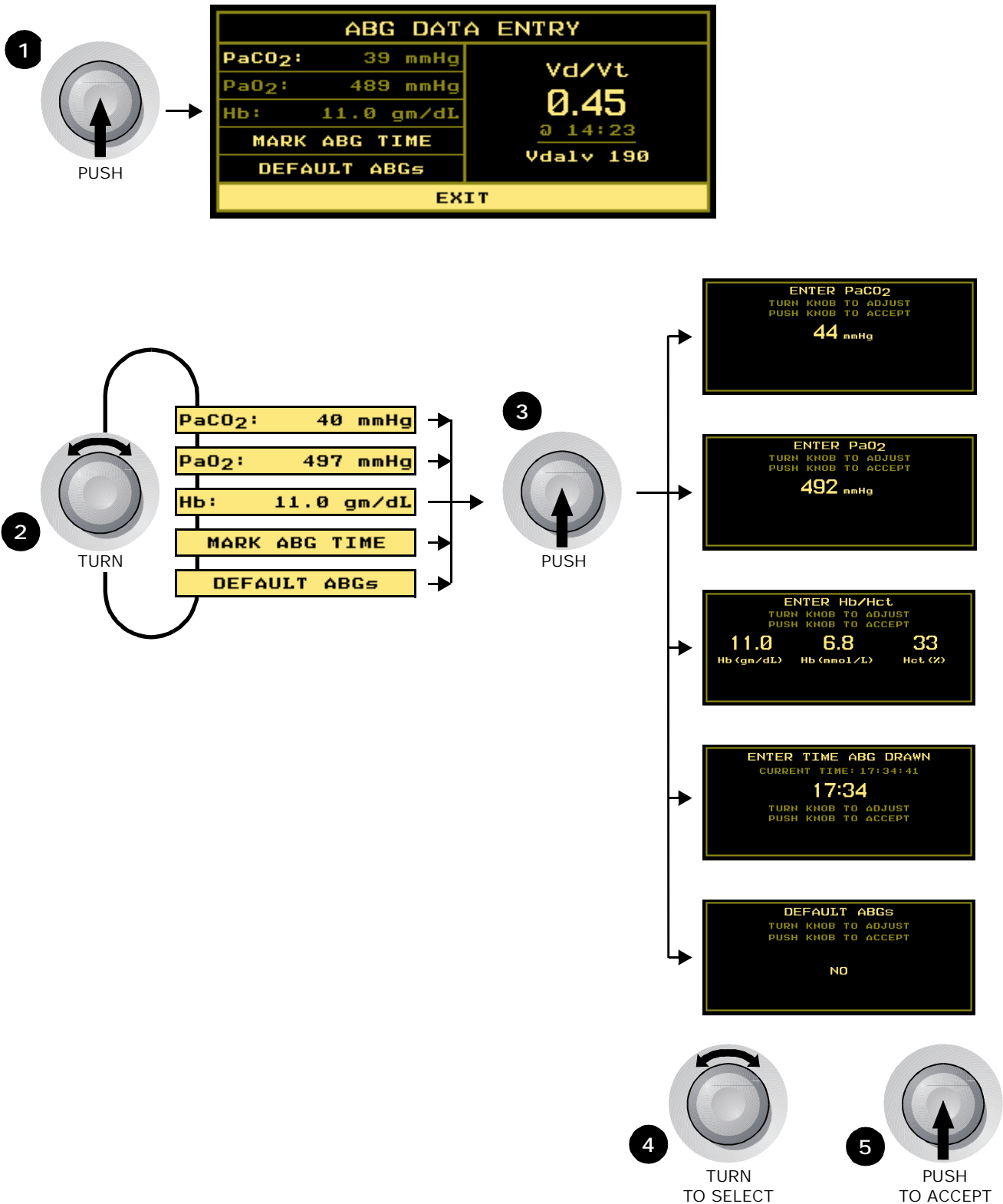



## Safety

For maximum patient and operator safety, observe the following warnings, cautions and notes.

### Warnings



#### WARNING:

Indicates a potentially harmful condition that can lead to personal injury.

- Explosion Hazard: Do not use the NICO<sup>®</sup> monitor in the presence of flammable anesthetics. Use of this instrument in such an environment may present an explosion hazard.
- Electrical Shock Hazard: Always turn the NICO<sup>®</sup> monitor off before cleaning it. Do not use with a damaged external power source. Refer servicing to qualified service personnel.
- Connect the AC Mains power cord to a properly grounded hospital-grade outlet. The NICO<sup>®</sup> monitor should be connected to the same electrical circuit as other equipment in use on the patient. Outlets of the same circuit can be identified by members of the hospital's engineering department.
- Failure of Operation: If the monitor fails to respond as described, do not use it until the situation has been corrected by qualified personnel.
- Reuse (disassembly, cleaning, disinfecting, resterilizing, etc.) of the CO<sub>2</sub>, CO<sub>2</sub>/Flow and NICO Sensors<sup>™</sup> may compromise device functionality and system performance and cause a potential patient hazard. Performance is not guaranteed if a sensor is reused.
- Inspect the CO<sub>2</sub>, CO<sub>2</sub>/Flow, SpO<sub>2</sub> and NICO Sensors<sup>™</sup> prior to use.
  - Do not use if they appear to be damaged or broken.
  - Do not attempt to rotate the NICO Sensor<sup>™</sup> in the breathing circuit by grasping the pneumatic tubes exiting the flow sensor.
  - Do not apply excessive tension to any cable or pneumatic tubing.
  - Periodically inspect sensor tubing lines for kinks.
  - Replace the CO<sub>2</sub>, CO<sub>2</sub>/Flow, or NICO Sensor<sup>™</sup> if excessive moisture or secretions are observed in the tubing.
- NICO<sup>®</sup> automatically identifies the type of sensor (small, standard or large NICO Sensor<sup>™</sup>, or neonatal, pediatric or adult CO<sub>2</sub>/Flow sensor ) when it is connected. If a sensor identification message is not displayed when a sensor is first connected, DO NOT use the sensor. If the condition persists, refer the monitor to qualified service personnel.
- Do not use the NICO<sup>®</sup> monitor if it is unable to properly identify a CO<sub>2</sub>/Flow sensor or a NICO Sensor<sup>™</sup>. If the condition persists, refer the monitor to qualified service personnel.
- In the event the message **NICO SENSOR FAILURE** is displayed, remove the NICO Sensor<sup>™</sup> from the patient circuit.
- The CO<sub>2</sub>/Flow or NICO Sensor<sup>™</sup> connector should be properly inserted into the front panel receptacle prior to connecting a sensor to the breathing circuit, in order to avoid a circuit leak, or occlusion of sensor tubing.
- NICO Sensors<sup>™</sup> increases airway deadspace by 35 cc (minimum). At low tidal volumes, compensatory changes to ventilation protocol should be considered.
- NICO Sensors<sup>™</sup> are not for pediatric use.

- Patient Safety: Care should be exercised to assure continued peripheral perfusion distal to the SpO<sub>2</sub> sensor site after application.
  - Inspect the SpO<sub>2</sub> sensor site for adequate circulation at least once every four hours.
  - When applying sensors take note of patient's physiological condition. For example, burn patients may exhibit more sensitivity to heat and pressure and therefore additional consideration such as more frequent site checks may be appropriate.
- Periodically check sensors and tubing for excessive moisture or secretion build up. Although NICO<sup>®</sup> automatically purges the lines, excessive moisture or secretions may still remain.
- While using the sensors, a system leak, such as that caused by uncuffed endotracheal tubes or a damaged sensor may significantly affect flow related readings. These include flow, volume, pressure, deadspace, CO<sub>2</sub> production and other respiratory mechanics parameters.
- Do not position sensor cables or tubing in any manner that may cause entanglement or strangulation.
- The NICO<sup>®</sup> monitor is not intended to be used as an apnea monitor.
- The NICO<sup>®</sup> monitor has no protection against the ingress of water.

## Cautions



### CAUTION:

Indicates a condition that may lead to equipment damage or malfunction.

- Use only Novamatrix approved sensors and accessories with the NICO<sup>®</sup> monitor.
- Do not operate the NICO<sup>®</sup> monitor when it is wet due to spills or condensation.
- Do not operate the product if it appears to have been dropped or damaged.
- Never sterilize or immerse the monitor in liquids.
- Do not sterilize or immerse sensors except as directed in this manual.
- No tension should be applied to any sensor cable or tubing.
- To avoid the effects of excessive moisture in the NICO Sensor™, insert it in the ventilator circuit with the pneumatic tubes upright. Excessive moisture in the NICO Sensor™ may affect the accuracy of the measurements.
- To avoid the effects of excessive moisture in the measurement circuit, insert the CO<sub>2</sub>/Flow sensor in the ventilator circuit with the tubes upright. Improper placement may result in erroneous data.
- Excessive moisture in the CO<sub>2</sub>/Flow sensor tubing may affect the accuracy of the measurements.
- It is recommended that CO<sub>2</sub>/Flow or NICO Sensors™ be removed from the circuit whenever an aerosolized medication is delivered. This is due to the increased viscosity of the medications which may contaminate the sensor windows, causing the sensor to fail prematurely.
- Operate the monitor at temperatures between 10 to +40° C (50 to 104° F), 10-95% R.H. non-condensing.
- Avoid storing the monitor at temperatures less than -10° C or greater than +55° C (<14° F or >131° F) 10-95% R.H. non-condensing
- Observe precautions for electrostatic discharge (ESD) and electromagnetic interference (EMI) to and from other equipment.
- Where electromagnetic devices (i.e., electrocautery) are used, patient monitoring may be interrupted due to electromagnetic interference. Electromagnetic fields up to 3 V/m will not adversely affect system performance.
- Caution: Federal (U.S.A.) law restricts this device to sale, distribution, or use by or on the order of a licensed medical practitioner.

**NOTE**

A point of particular interest or emphasis intended to provide more efficient or convenient operation.

- In order to ensure proper monitoring of oxygenation and ventilation:
  - The use of pulse oximetry is recommended during NICO® monitoring.
  - Setting of ETCO<sub>2</sub> and SpO<sub>2</sub> alert limits is recommended.
- A “NO RESPIRATION” alert is not generated when both the CAPNOSTAT® CO<sub>2</sub> sensor and the NICO Sensor™ or CO<sub>2</sub>/Flow sensor are disconnected from the NICO® monitor.
- Be certain that the monitor is not in Demo mode while monitoring. Demo mode can be identified by the flashing **DEMO MODE** label in the General Message area of the display. To exit Demo mode and return to normal monitoring mode, turn the power off and back on.
- The NICO® monitor contains no user serviceable parts. Refer servicing to qualified service personnel. A technical Service Manual is available for use by technical personnel.
- Do not attach an SpO<sub>2</sub> sensor distal to a blood pressure cuff. Valid data cannot be processed when the cuff is inflated. Attach the sensor to the limb opposite to the site used for the blood pressure cuff.
- This product and its accessories which have patient contact are free of latex.
- The NICO® monitor is Year 2000 compliant.
- Data Validity: Inaccurate SpO<sub>2</sub> and Pulse Rate values may be caused by:
  - Incorrect application or use of a sensor
  - Significant levels of dysfunctional hemoglobin; carboxyhemoglobin or methemoglobin
  - Significant levels of indocyanine green, methylene blue, or other intravascular dyes
  - Exposure to excessive illumination such as surgical lamps—especially ones with a xenon light source, or direct sunlight
  - Excessive patient movement
  - Venous pulsations
  - Electrosurgical interference
  - Use of an IABP.
- NICO® measurements will occur provided the following conditions are met:
  - The NICO Sensor™ assembly is properly installed in the patient’s breathing circuit.
  - Valid flow and CO<sub>2</sub> signals are detected with no significant signal artifact.
  - VCO<sub>2</sub> is greater than 20 mL/min.
  - ETCO<sub>2</sub> is between 15 and 85 mmHg (2.0 - 11.5 kPa or %) during baseline
  - ETCO<sub>2</sub> is between 15 and 100 mmHg (2.0 - 13.5 kPa or %) during rebreathing
  - The tidal volume is greater than 200ml (small and standard sizes)
  - The tidal volume is greater than 400 ml (large size).
  - The respiratory rate is between 3 and 60 br/min.
  - The **STOP/CONTINUE REBREATHING** key is not illuminated.
  - NICO® is not paused by the monitor for any other reason (displayed in the C.O. message area)
- When a new CAPNOSTAT® CO<sub>2</sub> sensor is attached to the monitor, or is moved from one monitor to another, it must be zeroed before use.
- After the life cycle of the equipment and accessories has been met, disposal should be accomplished following national and local requirements.



## Monitoring Cardiac Output

This section details the steps needed to begin patient monitoring with the NICO® monitor.

### Preparing for Use

#### Inspect


Before monitoring, take a few moments to inspect the NICO® monitor and its sensors. Check that all items are clean, dry, and physically intact with no broken or damaged components.

#### Turn on the monitor

Turn the NICO® monitor on.

- 1 Press the **Operate/Standby** key to turn the monitor on and off.

- NICO® can operate from its internal battery or from the AC Mains. (See “AC/Battery Operation” on page 4 for details.)

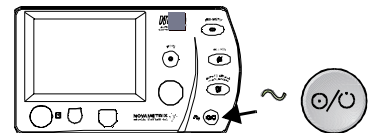

- 2 The monitor performs a quick self-test.

- An audible tone sounds, the key indicators illuminate, and a **SELF-TEST IN PROGRESS** message is briefly displayed.

- 3 The message **PRESS KNOB TO ERASE STORED TRENDS** is displayed for 5 seconds.

- To erase the contents of the monitor’s trend memory, press the knob. **TRENDS ERASED** is briefly displayed.
- To retain the contents of the monitor’s trend memory, do not press the knob. Wait the 5 seconds and **TRENDS RETAINED** is displayed.
- Note: If the internal battery becomes fully discharged, the message **CHECK DATE/TIME (MENU -> SETUP)** will appear before the **PRESS KNOB TO ERASE STORED TRENDS** message. (See “Setup Screen” on page 45 for details.)

- 4 The power up sequence is completed and a monitoring screen is displayed.

- NICO® displays the screen that was displayed when the monitor was last turned off.
- If the monitor is in Respiratory Mechanics mode, connecting a NICO Sensor™ will cause the monitor to automatically switch to Cardiac Output mode.
- The monitor is in a READY state and the parameters will be dashed and alerts will not be active until parameters are calculated and displayed.
- Parameters will display and their alerts will become active as they are calculated.

#### Connect and apply the sensors

Connect the sensors to the monitor, ventilator circuit, and patient.


- 1 Connect the SpO<sub>2</sub> sensor to the monitor and apply it to the patient. (See “Pulse Oximetry Sensors” on page 67.)


- 2 Connect the CAPNOSTAT® CO<sub>2</sub> Sensor to the monitor. (See “CAPNOSTAT® CO<sub>2</sub> Sensor” on page 65.)

- 3 Select a NICO Sensor™, see “Choosing a NICO Sensor™ size” on page 60.


- 4 Connect the NICO Sensor™ to the monitor and attach a CAPNOSTAT®. (See “NICO Sensor™” on page 59.)

- 5 Use the Initial Adjustment Template as a guide and adjust the NICO Loop™ to match the ventilator’s tidal volume setting, then discard the template. (See Instructions on the NICO® template.

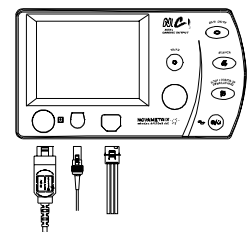

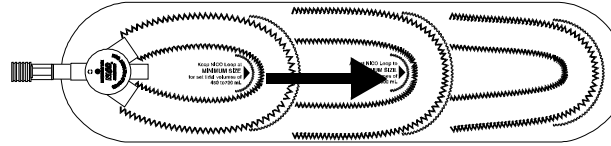


- 6 For optimal results, place the NICO Sensor™ into the ventilator circuit between the endotracheal tube and the ventilator circuit wye.
- Place other devices (HME, filters, etc.) between NICO Sensor™ and the patient connection.
  - The NICO Sensor™ increases airway deadspace by 35 cc (minimum). At low tidal volumes, compensatory changes to ventilation protocol should be considered.
  - Placement of a sidestream gas analyzer sampling port between the NICO Sensor™ and the patient connection may reduce NICO® accuracy at low tidal volumes.
  - Sidestream or mainstream gas analyzers placed between the NICO Sensor™ and the patient circuit “Y” may be inaccurate during the rebreathing phase of the NICO® cycle.
  - Place the sensor so that the triple lumen tubing lines exit from the top of the sensor (to help keep them clear and dry).
  - Keep the sensor clear of accumulations by proper circuit maintenance.

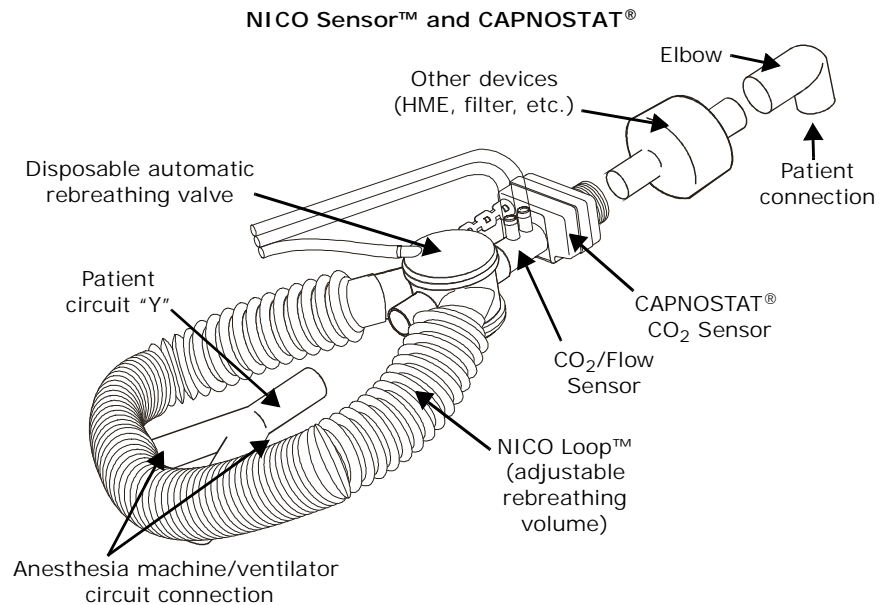

## Begin NICO® monitoring

After the NICO® monitor is turned on and the sensors are properly connected and applied, cardiac output monitoring can begin.

- 1 Press the **STOP/CONTINUE REBREATHING** key to initiate monitoring. Subsequent presses will Stop/Continue the rebreathing process.
  - An icon “⊘” identifies rebreathing status.
  - ⊘ Illuminated: Rebreathing is **DISABLED**.
  - ⊘ Not illuminated: Rebreathing is **ENABLED**.

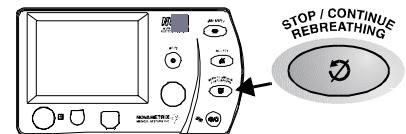

- 2 Enter the patient’s height and weight (for cardiac index calculations) and the delivered oxygen, anesthetic agent and balance gas by pressing the **DATA ENTRY** key. (See “Entering Patient Data” on page 28.)
- 3 If available, enter the patient’s PaCO<sub>2</sub>, PaO<sub>2</sub>, hemoglobin and hematocrit in the **ABG DATA ENTRY** screens, then mark the ABG time before exiting. (See “Entering Patient Data” on page 28.)

NOTE: Vd/Vt Phy, Vd Phy, and Vd alv are calculated based on the PaCO<sub>2</sub> value last entered, and are not updated until the next time the PaCO<sub>2</sub> value is entered or changed.

4 As you begin monitoring with NICO®, please note the following:

- Reliance on Cardiac Output parameters (NICO®, SV & CI) should be taken in context with other monitoring parameters and the physiologic condition of the patient.
- Entry of patient height and weight is required to calculate and display Cardiac Index.
- Pulse oximetry is required to calculate and display Stroke Volume (SV).
- Accuracy of cardiac output and related parameters will be affected by the following:
  - Significant fluctuations in mixed venous CO<sub>2</sub> content or metabolic CO<sub>2</sub> production during any three minute measurement period.
  - Sudden release of CO<sub>2</sub> into the bloodstream, such as when releasing a cross clamp.
  - The presence of excessive moisture or secretions in the NICO® sensor.
  - Entry of blood gas information.
- Due to the periodic rebreathing for NICO® measurements, the patient's effective ventilation will be reduced by typically 10-15% (depending on the rebreathing volume required). This can be offset by increasing the minute ventilation before NICO® monitoring begins.

### Rebreathing Bar

The Rebreathing Bar visually indicates the level of patient rebreathing. Under normal monitoring conditions, the **REBREATHING** icon and Rebreathing Bar appear in the message center each time the automatic rebreathing cycle begins.



- The Rebreathing Bar represents the total possible range of rebreathing from 0-100%.
- The highlighted area within the bar represents the target rebreathing range (50%-90%) for optimal NICO® performance.
- A vertical pointer within the Rebreathing Bar indicates the current rebreathing percentage.
  - The pointer will appear within the highlighted area of the Rebreathing Bar when the NICO Loop™ is properly sized and providing an acceptable percentage of rebreathing.
  - *Note:* The word "OK" will disappear when the pointer falls within that area of the highlighted bar.
- The Rebreathing Bar is not displayed when the ratio of spontaneous to mechanical breaths is greater than or equal to 2:1.

### Expand/Retract NICO Loop™

- If the NICO Loop™ is not sufficiently expanded and the Rebreathing Bar pointer falls below 50%, the message **EXPAND LOOP** is displayed during the rebreathing period.
- If the NICO Loop™ is over-expanded and the Rebreathing Bar pointer is above 90%, the message **RETRACT LOOP** is displayed during the rebreathing period.



To expand or retract the NICO Loop™:

- 1 Grasp the NICO Loop™ with one hand and the automatic rebreathing valve with the other hand so as not to disturb/disconnect the breathing circuit while adjusting the loop.
- 2 Expand or retract the NICO Loop™ 3-6 inches.
  - It may take 2-3 additional breaths before the icon changes.
  - Note that if the loop is still not appropriately sized by the end of the rebreathing period, the message will be removed and may be displayed again during the next rebreathing period.

### Sensor Size Message

If the **EXPAND** or **RETRACT LOOP** message appears for more than three rebreathing cycles, and resizing the NICO Loop™ was not effective, the NICO® monitor will suggest a different sized sensor to correct the condition. See "Status Messages" on page 61.





## The NICO® Cycle

Once rebreathing is enabled, the NICO® monitor automatically repeats a three minute cardiac output measurement cycle. This NICO® cycle has three phases:

- **Baseline:** During the 60-seconds baseline period the rebreathing valve inside the NICO Sensor™ is turned off and the rebreathing volume of the NICO Loop™ is bypassed. During this time, VCO<sub>2</sub>, PaCO<sub>2</sub> and ETCO<sub>2</sub> will be at baseline values.
- **Rebreathing:** The 50-second rebreathing period starts when the monitor turns on the rebreathing valve inside the NICO Sensor™ causing the rebreathing volume of the NICO Loop™ to be added into the circuit. During rebreathing, VCO<sub>2</sub> is reduced, PaCO<sub>2</sub> and ETCO<sub>2</sub> become elevated (3-5 mmHg, typical) and mixed venous CO<sub>2</sub> remains unchanged.

**NOTE:** The rebreathing period will typically induce an increase in PaCO<sub>2</sub> by 3-5 mmHg. An ABG blood sample drawn during this period (☞ **REBREATHING** displayed) or during the first twenty seconds of the stabilization period (where **NEXT ☞** is displayed), may cause PaCO<sub>2</sub> values to reflect higher than normal levels.

- **Stabilization:** After completion of rebreathing, a 70-second stabilization period begins, during which time VCO<sub>2</sub>, PaCO<sub>2</sub> and ETCO<sub>2</sub> return to their baseline values.

NICO® updates the displayed C.O. value following the completion of each three minute NICO® cycle. The CObar™ (cardiac output confidence bar) provides an indication of the system's confidence in the displayed value. (See "CObar™ Confidence Indicator" on page 31.)

## Rebreathing On/Off or Paused

The user can interrupt or resume the rebreathing cycle at any time by pressing the **STOP/CONTINUE REBREATHING** key. The NICO® monitor will not automatically restart the rebreathing cycle—that must be initiated by the user pressing the **STOP/CONTINUE REBREATHING** key. Note that once rebreathing is initiated by the user, the monitor will, under certain conditions, pause rebreathing until a specified condition (see below) is corrected—at which time the monitor will restart the rebreathing.

Rebreathing can be Off, On, or Paused as indicated below:

### Rebreathing OFF (disabled)

- The NICO® monitor starts in this state upon power-up.
- The **STOP/CONTINUE REBREATHING** key is illuminated while in this state.
- The NICO® rebreathing cycle can be placed into this state at any time by pressing the **STOP/CONTINUE REBREATHING** key (rebreathing is immediately disabled).
- The NICO® rebreathing cycle is automatically turned off for certain monitor/sensor conditions. (See "(Alert Class: H-High Priority, M-Medium Priority, L-Low Priority, S-Status Message. See "NICO® Alert Priorities" on page 56 for details." on page 79.)
- Noted in the cardiac output message area as **REBREATHING OFF**.

REBREATHING OFF  
PRESS KEY TO  
CONTINUE

### Rebreathing ON (enabled)

- When the monitor is initially turned on, this state is entered only after pressing the **STOP/CONTINUE REBREATHING** key.
- The **STOP/CONTINUE REBREATHING** key is not illuminated while in this state.
- The cardiac output is calculated and updated while in this state.
- Noted in the cardiac output message area as **☞ REBREATHING** or **NEXT ☞**.

REBREATHING  
OK

NEXT : 1:57

### Rebreathing Paused

- The NICO® monitor automatically pauses the rebreathing cycle and generates a display message under any of these conditions:
  - ETCO<sub>2</sub> is less than 15 mmHg (2.0 kPa or %) or greater than 85 mmHg (11.5 kPa or %)
  - Respiration rate is less than 3 or greater than 60 br/min.
  - VCO<sub>2</sub> is less than 20 mL/min.
- The rebreathing cycle automatically restarts when the condition is corrected.

REBREATHING PAUSED  
WAITING FOR  
RR > 3 br/m

## Entering Patient Data

NICO® monitoring can be enhanced by the entry of key patient specific data including respired gas composition (anesthetic agent, balance gas, and inspired O<sub>2</sub>), patient height and weight, and arterial blood gas data (PaCO<sub>2</sub>, PaO<sub>2</sub>, Hb or Hct). Inclusion of ABG data is especially important when gas exchange impairment is expected (i.e., high shunt or deadspace). **ABG samples should not be obtained during the rebreathing phase of the 3-minute NICO® cycle.**

Patient data should be updated in the **DATA ENTRY** screen whenever possible. The screen may be accessed at any time by pressing the **DATA ENTRY** key.


**DATA ENTRY settings** The following table lists the parameters and ranges accessible in the **DATA ENTRY** screens.

Label	Parameter	Default	Range/Units	Description
INSP O <sub>2</sub>	Inspired Oxygen	70%	21-100 %	Percent of oxygen in the inspired gas. Must be entered in order for NICO® to accurately calculate parameters.
INSP AGENT	Inspired Anesthetic Agent	0%	0-20 %	Percent of anesthetic agent in the inspired gas. Must enter percent delivered in order for NICO® to accurately calculate parameters.
BALANCE	Gas Balance	N <sub>2</sub>	N <sub>2</sub> , He, or N <sub>2</sub> O	N <sub>2</sub> , He or N <sub>2</sub> O. Must select the correct balance in the inspired gas in order for NICO® to accurately calculate parameters.
HEIGHT	Patient Height	--	35-91 in 90-230 cm	Enter patient height for CI calculations.
WEIGHT	Patient Weight	--	55-551 lb 25-250 kg	Enter patient weight for CI calculations.
<b>ABG DATA ENTRY Screen</b>				
PaCO <sub>2</sub>	Arterial Carbon Dioxide	40 mmHg (5.4 kPa or %)  ("--" displayed until an initial value is entered)	0-250 mmHg 0.0-20.0 kPa 0.0-20.0 %	Partial pressure of carbon dioxide in arterial blood. Entering this value can enhance the accuracy of NICO® parameters.
PaO <sub>2</sub>	Arterial Oxygen	FiO <sub>2</sub> (Pb-47 mmHg) ("--" displayed until an initial value is entered)	0-750 mmHg 0.0-99.5 kPa 0.0-99.5 %	Partial pressure of oxygen in arterial blood. Entering this value can enhance the accuracy of NICO® parameters.
Hb	Hemoglobin Concentration or Hematocrit	11.0 gm/dL 6.8 mmol/L 33 % ("--" displayed until an initial value is entered)	Hb: 5.0-20.0 gm/dL Hb: 3.1-12.4 mmol/L Hct: 0-60 %	Concentration of hemoglobin or hematocrit in the blood. Entering this value can enhance the accuracy of NICO® parameters.
MARK ABG TIME	Time when ABG blood sample is drawn	Current Time	hh:mm (hours: minutes)	Enter time ABG is drawn. (Only accepts time since ETCO <sub>2</sub> was first detected.)
DEFAULT ABGs	Blood gas values	PaCO <sub>2</sub> : 5 mmHg (0.7 kPa or %) above the measured ETCO <sub>2</sub> value <sup>a</sup>  PaO <sub>2</sub> : Based on barometric pressure and INSP O <sub>2</sub> value.  Hb: 11.0 gm/dL		Resets blood gas settings to default values.

a. 40mmHg (5.4 kPa or %) if ETCO<sub>2</sub> is not available.

**Entering Patient Data** To enter (or view) patient data:

- 1 Press the **DATA ENTRY** key to activate **DATA ENTRY**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 2 Highlight the desired data by turning the **KNOB**.
- 3 Select the highlighted data by pressing the **KNOB**.
- 4 Turn the **KNOB** to adjust the value as desired.
- 5 Press the **KNOB** to accept the value.
- 6 Repeat these steps for the other settings.

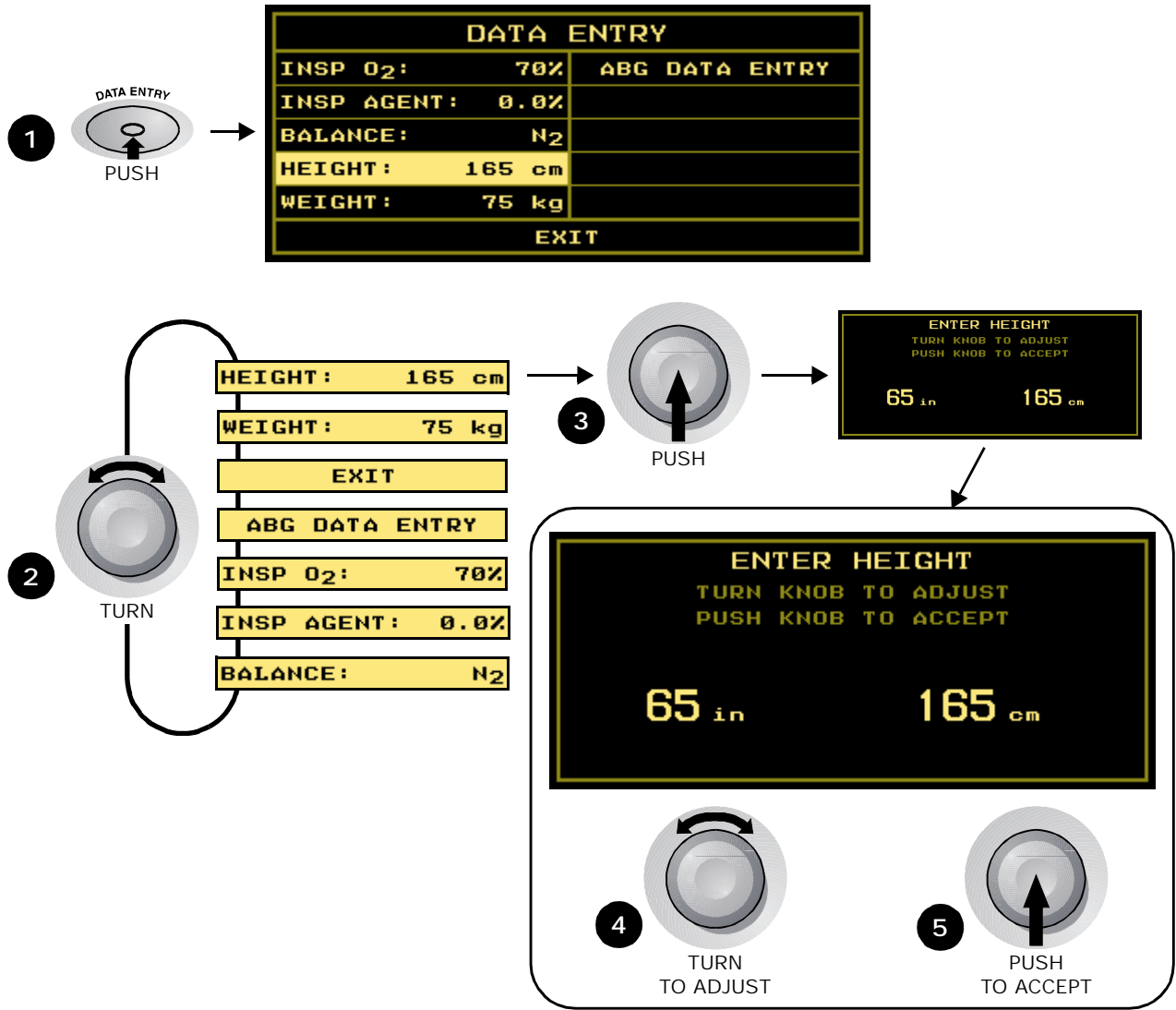




**Entering ABG Data**

To enter ABG data:

- 1 Draw blood sample during the "normal" phase of the NICO® cycle, not during the rebreathing phase.
  - The cardiac output message will display how much time before the next rebreathing phase occurs.

NEXT ⌚: 1:57

- 2 Press the **DATA ENTRY** key to activate **DATA ENTRY**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 3 Highlight **ABG DATA ENTRY** by turning the **KNOB**.
- 4 Select **ABG DATA ENTRY** by pressing the **KNOB**.

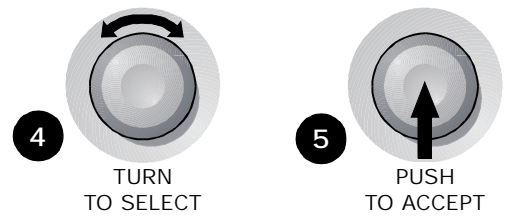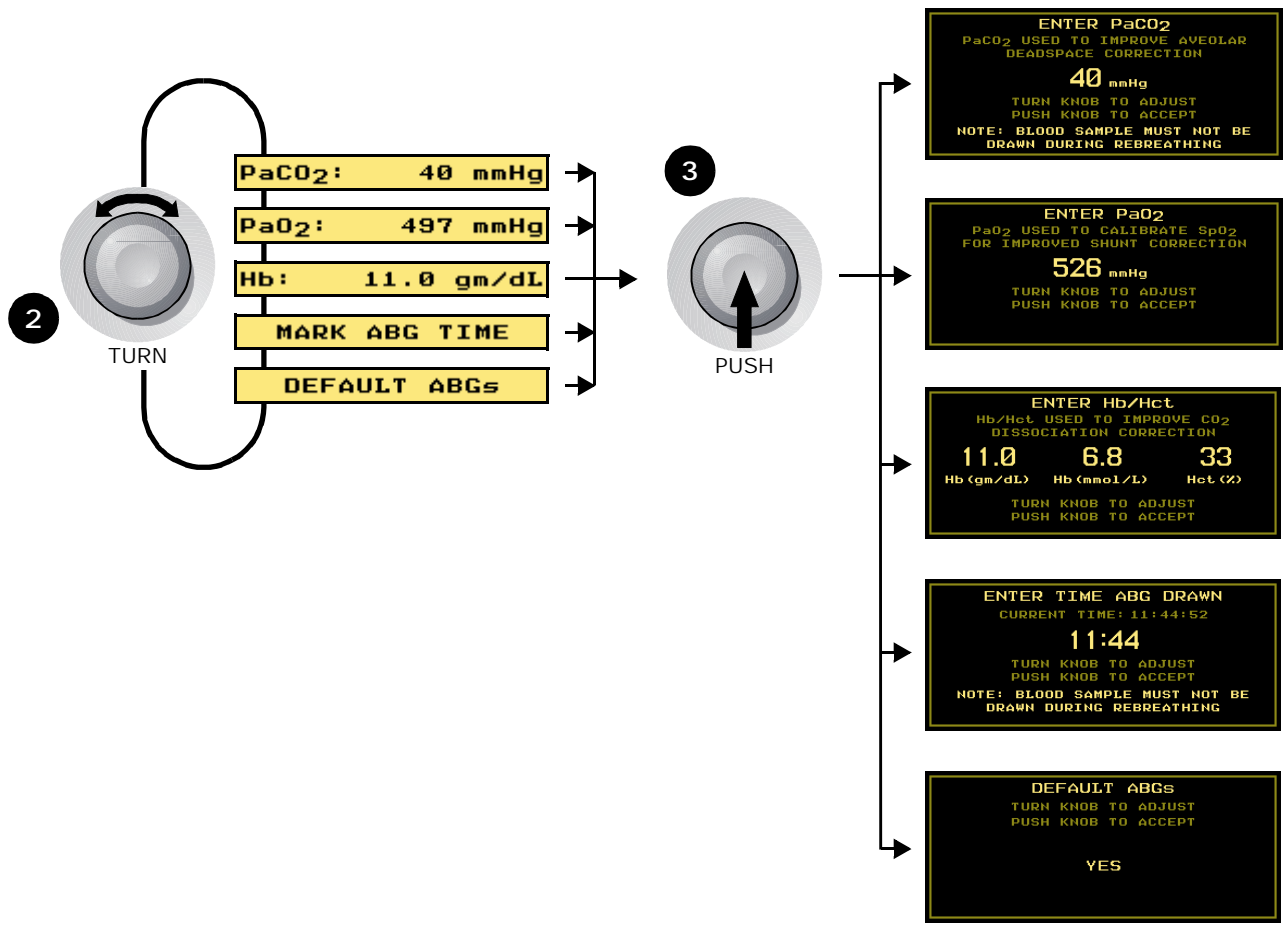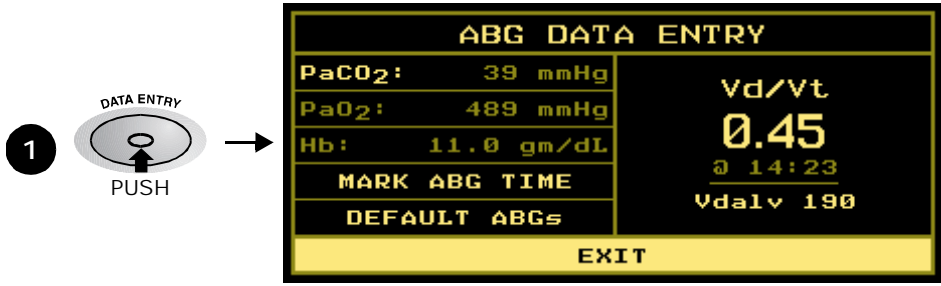

- 5 Turn and press the **KNOB** to select the PaCO<sub>2</sub>, PaO<sub>2</sub>, Hb/Hct, Mark ABG Time or Default ABG screen.
- 6 Turn **KNOB** to adjust value.
  - ABG Time is required for PaCO<sub>2</sub> and PaO<sub>2</sub>.
  - Ensure that the NICO® clock is synchronized to the clock used to determine the draw time.
- 7 Press the **KNOB** to accept the value and return to the **ABG DATA ENTRY** screen.
  - Updated Vd/Vt and Vd<sub>alv</sub> values are displayed.

## C.O. Averaging

### CObar™ Confidence Indicator

The Cardiac Output Confidence Bar, or CObar™, is an indicator of the system's confidence in the currently displayed cardiac output value. The CObar indicator is located above the cardiac output value and can contain up to five segments. The degree of confidence (more segments for higher confidence, fewer segments for lower confidence) is based on multiple factors including ventilatory pattern, sizing of the NICO Loop™, tidal volume, entry of patient data and breathing circuit integrity.



At higher confidence levels the system can more quickly resolve changes in C.O. and is therefore able to display those changes more quickly. At lower confidence levels, the system must work harder at resolving C.O. changes, thereby delaying how quickly that information can be displayed. The delay from detection to display of changes in C.O. can range from:

- 5 segments: 1.5 - 2 minutes
- 4 segments: 3 - 7 minutes
- 3 segments: 5 - 9.5 minutes
- 2 segments: 7 - 9.5 minutes
- 1 segment: 8.5 - 9.5 minutes

A lower number of segments indicates that the displayed cardiac output reading is being averaged more with readings from previous NICO® rebreathing cycles.

If there is no confidence in the signals, the C.O. value is not displayed and no segments are displayed in the CObar.

The CObar is displayed when the **C.O. MODE** is set to **AVERAGE** in the **SETUP** screen.

### C.O. Fast Mode

When the **C.O. MODE** is set to **FAST**, the monitor will display the unfiltered C.O. value (CO-f) rather than the averaged value. The text **FAST MODE** will replace the CObar graphic. The stroke volume and cardiac index will be calculated from this value rather than the averaged value.



The averaged C.O. value can still be viewed as **CO-a** under the VCO<sub>2</sub> value in the **LAST COMPLETED CYCLE** screen.

The **C.O. MODE** option is set in the **SETUP** screen.



## Respiratory Monitoring

This section details the steps needed to begin patient respiratory mechanics monitoring with the NICO® monitor.

### Preparing for Use

#### Inspect


Before monitoring, take a few moments to inspect the NICO® monitor and its sensors. Check that all items are clean, dry, and physically intact with no broken or damaged components.

#### Turn on the monitor

Turn the NICO® monitor on.

- 1 Press the **Operate/Standby** key to turn the monitor on and off.

- NICO® can operate from its internal battery or from the AC Mains. (See “AC/Battery Operation” on page 4 for details.)

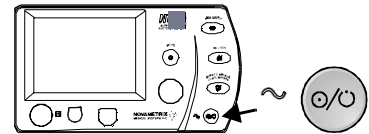

- 2 The monitor performs a quick self-test.

- An audible tone sounds, the key indicators illuminate, and a **SELF-TEST IN PROGRESS** message is briefly displayed.

- 3 The message **PRESS KNOB TO ERASE STORED TRENDS** is displayed for 5 seconds.

- To erase the contents of the monitor’s trend memory, press the knob. **TRENDS ERASED** is briefly displayed.
- To retain the contents of the monitor’s trend memory, do not press the knob. Wait the 5 seconds and **TRENDS RETAINED** is displayed.
- Note: If the internal battery becomes fully discharged, the message **CHECK DATE/TIME (MENU -> SETUP)** will appear before the **PRESS KNOB TO ERASE STORED TRENDS** message. (See “Setup Screen” on page 45 for details.)


- 4 The power up sequence is completed and a monitoring screen is displayed.

- NICO® displays the screen that was displayed when the monitor was last turned off.
- If the monitor is in Cardiac Output mode, connecting a CO<sub>2</sub>/Flow sensor will cause the monitor to automatically switch to Respiratory Mechanics mode.
- The monitor is in a **READY** state and the parameters will be dashed and alerts will not be active until parameters are calculated and displayed.
- Parameters will display and their alerts will become active as they are calculated.


#### Connect and apply the sensors

Connect the sensors to the monitor, ventilator circuit, and patient.


- 1 Connect the SpO<sub>2</sub> sensor to the monitor and apply it to the patient. (See “Pulse Oximetry Sensors” on page 67.)
- 2 Connect the CAPNOSTAT® CO<sub>2</sub> Sensor to the monitor. (See “CAPNOSTAT® CO<sub>2</sub> Sensor” on page 65.)
- 3 Select an appropriate size CO<sub>2</sub>/Flow sensor.
- 4 Connect a CO<sub>2</sub>/Flow Sensor to the monitor and attach a CAPNOSTAT®. (See “CO<sub>2</sub>/Flow Sensors” on page 62.)

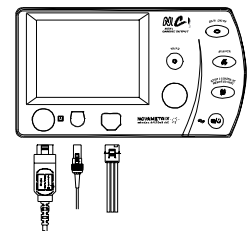

5 Connect the CAPNOSTAT CO<sub>2</sub> sensor to the CO<sub>2</sub>/Flow sensor.

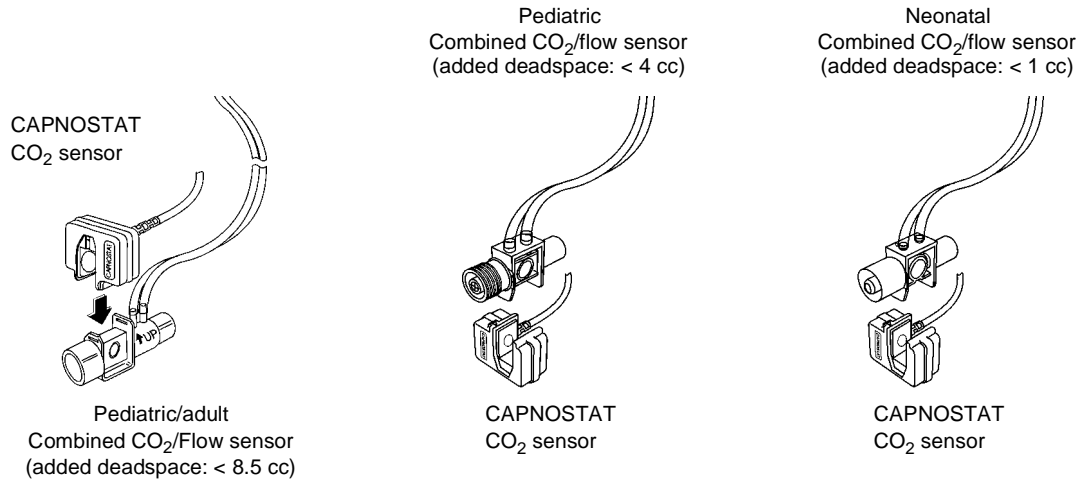

- 6 For optimal results, place the CO<sub>2</sub>/Flow sensor in the ventilator circuit between the endotracheal tube and the ventilator circuit wye.
- Place the sensor proximal to the patient if using other devices in the circuit.
  - Place the sensor so that the tubing lines exit from the top of the sensor (to help keep them clear and dry).
  - Keep the sensor clear of accumulations by proper circuit maintenance.
- 7 Connect the combined CO<sub>2</sub>/Flow sensor to the patient breathing circuit.

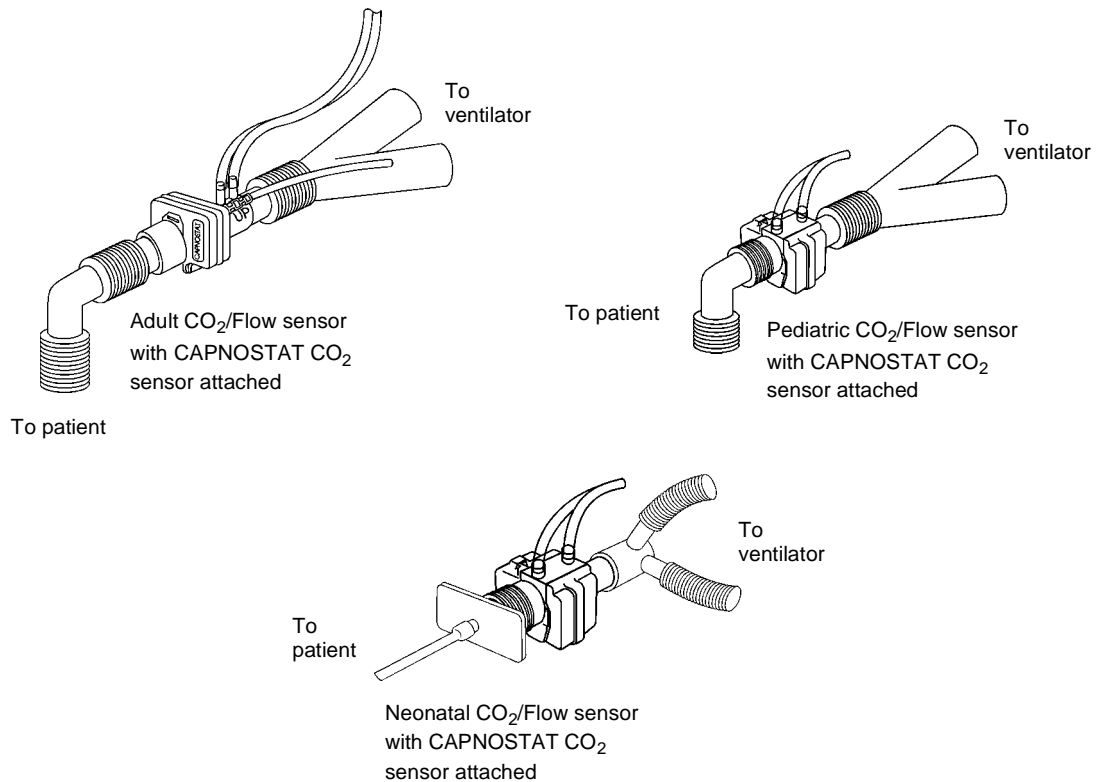


- Do NOT place the CO<sub>2</sub>/Flow sensor between the ET tube and the elbow (pediatric/adult circuit), as this may allow patient secretions to block the adapter windows.
- Position the CO<sub>2</sub>/Flow sensor with its windows in a vertical and NOT a horizontal position: this helps keep patient secretions from “pooling” on the windows.
- To prevent “rain-out” and moisture from draining into the CO<sub>2</sub>/Flow sensor, do NOT place the CO<sub>2</sub>/Flow sensor in a gravity dependent position.
- Periodically check the CO<sub>2</sub>/Flow sensor and tubing for excessive moisture or secretion build up.
- For routine performance of airway care, separate the system between the ET tube and the airway adapter (neonatal circuit), or between the ET tube and elbow (pediatric/adult circuit). Lavage and suctioning of the airway can then be performed without fluids and mucous accumulating on the CO<sub>2</sub>/Flow sensor windows.

## Begin Respiratory Monitoring

After the NICO® monitor is turned on and the sensors are properly connected and applied, Respiratory Mechanics monitoring can begin.

- 1 Enter the delivered oxygen and balance gas or anesthetic agent if present, by pressing the **DATA ENTRY** key. (See “Entering Patient Data” on page 34.)
- 2 As you begin monitoring with NICO®, please note the following:

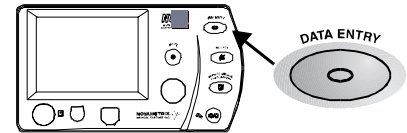

- Water that accumulates in the CO<sub>2</sub>/Flow sensor or the sensor tubing may cause the reported Tidal Volumes to be higher than set volumes. If reported values are higher (or lower) than expected and water is seen in the line or sensor body, purge the lines. If purging does not clear the water, remove the sensor from the circuit and remove the water by shaking the sensor, or by flowing oxygen or compressed air through the tubing or sensor until the water is removed. Do not use high pressure for water removal.
- To minimize the effects of aerosolized medications on the CO<sub>2</sub>/Flow sensor, it is recommended that the CO<sub>2</sub>/Flow sensor be removed from the ventilator circuit prior to the delivery of the medication. ***The decision to remove or not remove the CO<sub>2</sub>/Flow sensor is the responsibility of the clinician.***
- During the purge cycle the pump will be heard.
- Water will condense in the pressure sense lines at a faster rate when used in cooler ambient temperatures.
- Always keep the CO<sub>2</sub>/Flow sensor tubing pointed in an upward position to minimize pooling of water and secretions at the pressure sense line openings.
- The automatic purge mode may not be disabled.
- During a very low battery condition, automatic and manual purging is not allowed.

## Entering Patient Data

Respiratory monitoring can be enhanced by the entry of key patient specific data including respired gas composition (anesthetic agent, balance gas, and inspired O<sub>2</sub>).

Patient data should be updated in the **DATA ENTRY** screen whenever possible. The screen may be accessed at any time by pressing the **DATA ENTRY** key.



**DATA ENTRY settings** The following table lists the parameters and ranges accessible in the **DATA ENTRY** screens.

Label	Parameter	Default	Range/Units	Description
INSP O <sub>2</sub>	Inspired Oxygen	70%	21-100 %	Percent of oxygen in the inspired gas. Must be entered in order for NICO® to accurately calculate parameters.
INSP AGENT	Inspired Anesthetic Agent	0%	0-20 %	Percent of anesthetic agent in the inspired gas. Must enter percent delivered in order for NICO® to accurately calculate parameters.
BALANCE	Gas Balance	N <sub>2</sub>	N <sub>2</sub> , He, or N <sub>2</sub> O	N <sub>2</sub> , He, or N <sub>2</sub> O. Must select the correct balance in the inspired gas in order for NICO® to accurately calculate parameters.
HEIGHT	Patient Height	--	35-91 in 90-230 cm	Enter patient height.
WEIGHT	Patient Weight	--	Neonatal: 0.22 - 44.09 lb 0.10 - 20.00 kg  Pediatric: 0.2 - 220.2 lb 0.1 - 99.9 kg  Adult: 55-551 lb 25-250 kg	Enter patient weight for respiratory mechanics calculations.


**ABG DATA ENTRY Screen**

PaCO <sub>2</sub>	Arterial Carbon Dioxide	40 mmHg (5.4 kPa or %)  ("--" displayed until an initial value is entered)	0-250 mmHg 0.0-20.0 kPa 0.0-20.0 %	Partial pressure of carbon dioxide in arterial blood. Enter this value for calculation of Vd alv (alveolar deadspace), Vd/Vt (deadspace to tidal volume ratio) and Vd phys (physiologic deadspace).
PaO <sub>2</sub>	Arterial Oxygen	FiO <sub>2</sub> (Pb-47 mmHg) ("--" displayed until an initial value is entered)	0-750 mmHg 0.0-99.5 kPa 0.0-99.5 %	Partial pressure of oxygen in arterial blood.  Enter a value for this parameter if desired; does not affect CO <sub>2</sub> /flow calculations.
Hb	Hemoglobin Concentration or Hematocrit	11.0 gm/dL 6.8 mmol/L 33 % ("--" displayed until an initial value is entered)	Hb: 5.0-20.0 gm/dL Hb: 3.1-12.4 mmol/L Hct: 0-60 %	Concentration of hemoglobin or hematocrit in the blood.  Enter a value for parameter if desired; does not affect CO <sub>2</sub> /flow calculations.
MARK ABG TIME	Time when ABG blood sample is drawn	Current Time	hh:mm (hours:minutes)	Enter time ABG is drawn. Only accepts time since ETCO <sub>2</sub> was first detected.
DEFAULT ABGs	Blood gas values	PaCO <sub>2</sub> : 5 mmHg (0.7kPa or %) above the measured ETCO <sub>2</sub> value <sup>a</sup>  PaO <sub>2</sub> : Based on barometric pressure and INSP O <sub>2</sub> value.  Hb: 11.0 gm/dL		Resets blood gas settings to default values.

a. 40mmHg (5.4 kPa or %) if ETCO<sub>2</sub> is not available.

**Entering Patient Data** To enter (or view) patient data:

- 1 Press the **DATA ENTRY** key to activate **DATA ENTRY**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 2 Highlight the desired data by turning the **KNOB**.
- 3 Select the highlighted data by pressing the **KNOB**.
- 4 Turn the **KNOB** to adjust the value as desired.
- 5 Press the **KNOB** to accept the value.
- 6 Repeat these steps for the other settings.



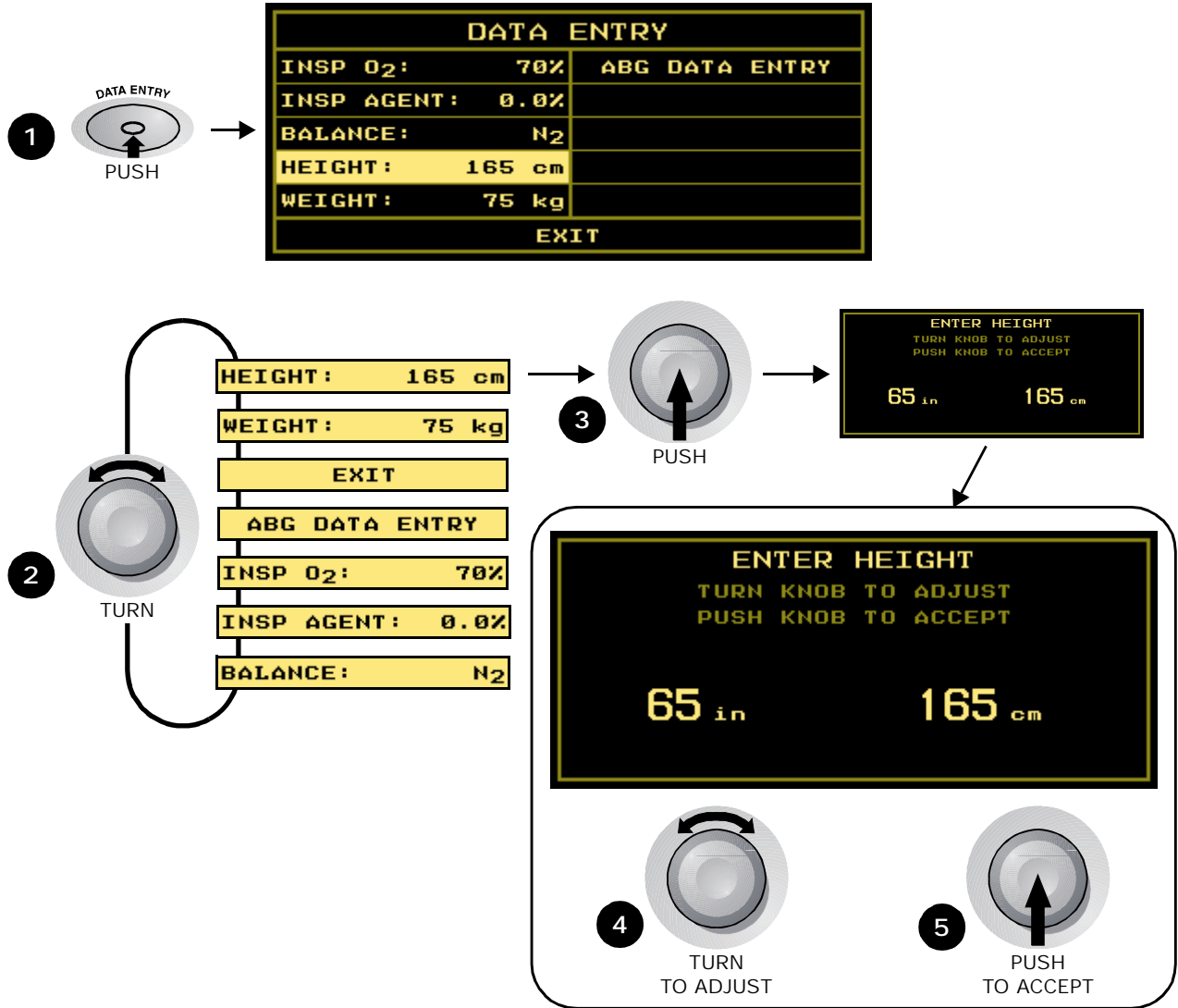


**Entering ABG Data**

To enter ABG data:

- 1 Press the **DATA ENTRY** key to activate **DATA ENTRY**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 2 Highlight **ABG DATA ENTRY** by turning the **KNOB**.

3 Select **ABG DATA ENTRY** by pressing the **KNOB**.

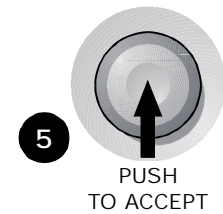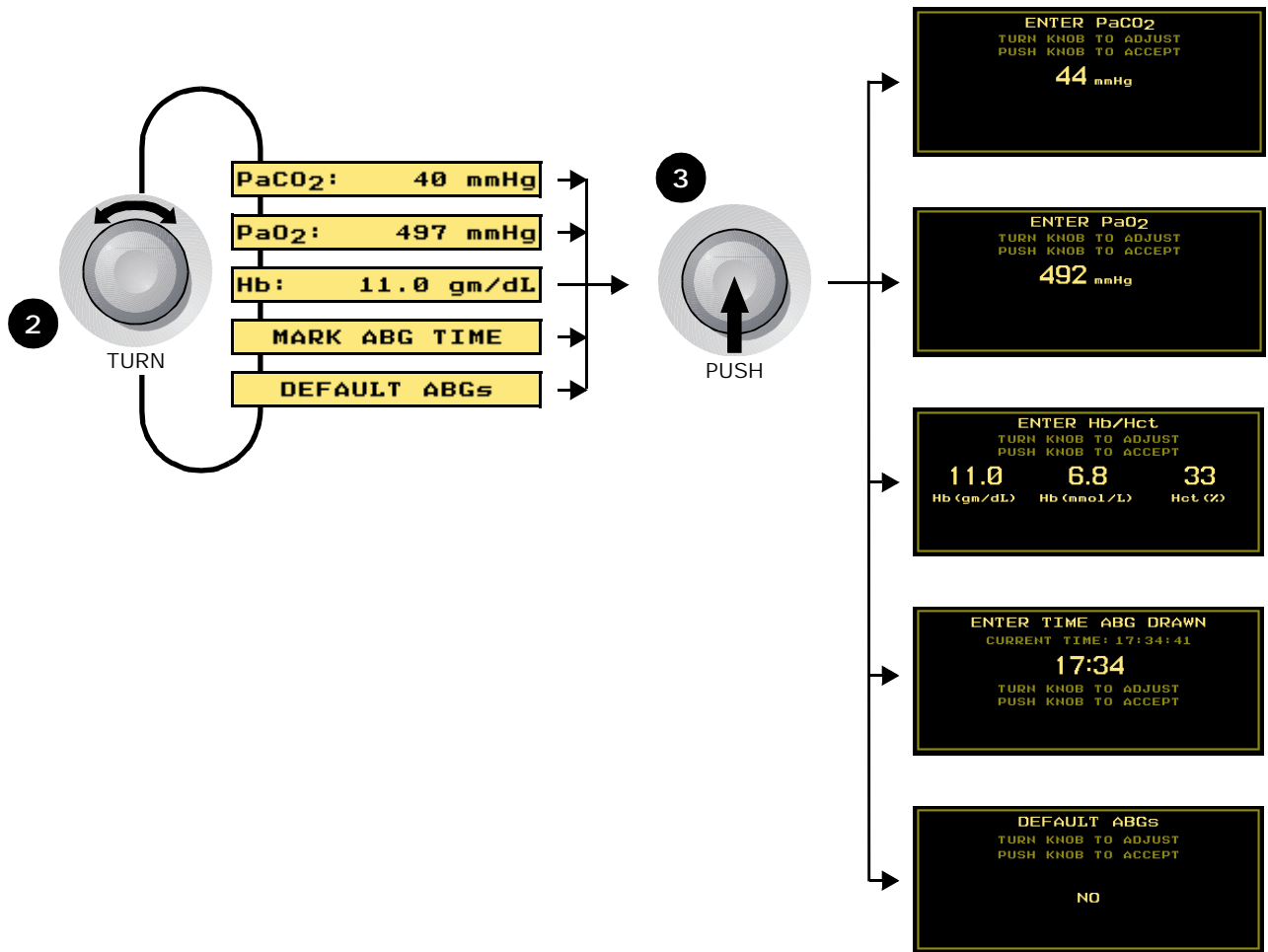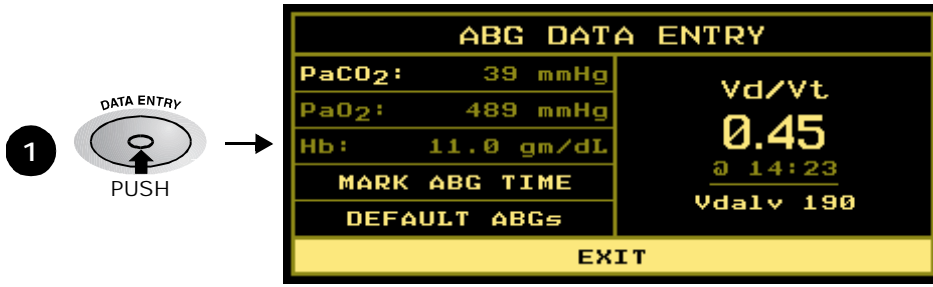

- 4 Turn and press the **KNOB** to select the PaCO<sub>2</sub>, PaO<sub>2</sub>, Hb/Hct, Mark ABG Time or Default ABG screen.

NOTE: PaO<sub>2</sub> and Hb selections are dimmed in Respiratory Mechanics mode to denote that, although they are available for data input, the data will not affect parameter calculations.

- 5 Turn **KNOB** to adjust value.
  - ABG Time is required for PaCO<sub>2</sub>.
  - Ensure that the NICO® clock is synchronized to the clock used to determine the draw time.
- 6 Press the **KNOB** to accept the value.
  - Updated Vd/Vt and Vdalv values are displayed.




## Monitoring and Setup Screens

### Cardiac Output mode

This section provides an overview of monitoring and setup screens in Cardiac Output mode.

#### C.O. Trend Screen

The Cardiac Output Trend Screen plots cardiac output over time.

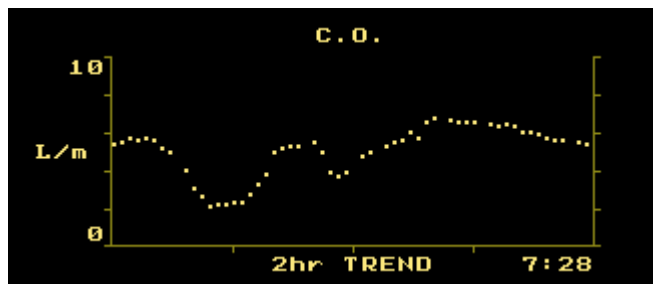



The **C.O. TREND** screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **C.O. TREND**.

- While viewing the **C.O. TREND** screen, each push of the **KNOB** will advance through the available 2, 4, 8 and 12 hour trend displays.
- Data is plotted from left (oldest data) to right (newest data). Once the display is filled, data shifts left so that the oldest data point on the left is pushed out to make room for the newest point entering from the right.
- The current time is shown in the lower right corner of the display and represents time at the right edge of the screen. See “Setup Screen” on page 45 to set the time.
- Each point on the trend (plotted or blank) represents the average C.O. value over a time period. The time periods are: 1 minute average for the 2 hour trend, 2 minutes for the 4 hour trend, 4 minutes for the 8 hour trend, and 6 minutes for the 12 hour trend.
- C.O. trends are automatically scaled to fit 0-5, 0-10, 0-15 and 0-20 L/min scales.
- A two-pixel wide dashed vertical line in the trend is used to denote a power-cycle where the NICO® monitor was turned off and back on.

#### Rebreathing Curves Screens

There are two Rebreathing Curves Screens—one displays the last completed NICO® cycle, and the other shows the current cycle. Both plot  $VCO_2$  and  $ETCO_2$  values over time and also provide numeric displays for  $VCO_2$ ,  $ETCO_2$  and pulse rate. Additionally, CO-f (C.O. fast-mode, unaveraged) is reported for the last completed cycle.

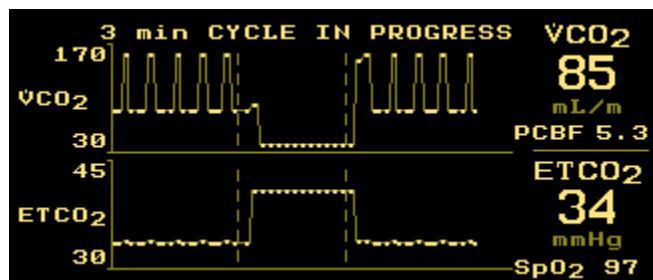

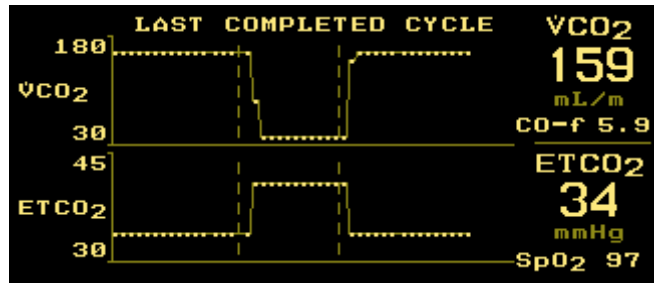


The Rebreathing Curves Screens can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **REBREATH CURVES**. To switch between the two curves, push the **KNOB**.

- Data is plotted from left (oldest data) to right (newest data).
- Points on the trend (plotted or blank) represent the VCO<sub>2</sub> or ETCO<sub>2</sub> value for each breath.
- The curves are automatically scaled to fit the display area.
- Two one-pixel wide dashed vertical lines are used to divide the curve into its baseline, rebreathing and stabilization phases.
- ETCO<sub>2</sub> and SpO<sub>2</sub> values will flash if an alert limit is exceeded.

### CO<sub>2</sub> and SpO<sub>2</sub> Waveform Screen

The CO<sub>2</sub> and SpO<sub>2</sub> Waveform Screen plots the capnogram and plethysmogram signals as well as providing numeric displays of ETCO<sub>2</sub>, Respiratory rate, SpO<sub>2</sub> and pulse rate (♥).

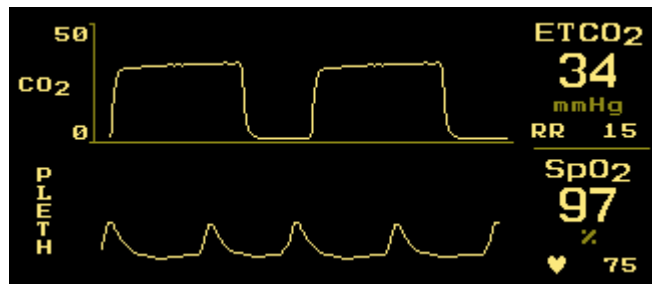


The CO<sub>2</sub>/SpO<sub>2</sub> screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **CO<sub>2</sub>/SpO<sub>2</sub>**.

- Information is updated in real time and reflects current patient status.
- The capnogram and plethysmogram are automatically scaled.
- The capnogram sweep speed is selectable in the **SETUP** menu.
- ETCO<sub>2</sub>, RR, SpO<sub>2</sub> and pulse rate values will flash if an alert limit is exceeded.

### Flow and Pressure Waveform Screen

The Flow and Pressure Waveform screen plots the airway flow and pressure signals over time as well as providing numeric displays of Vte, Vti, PIP and PEEP.

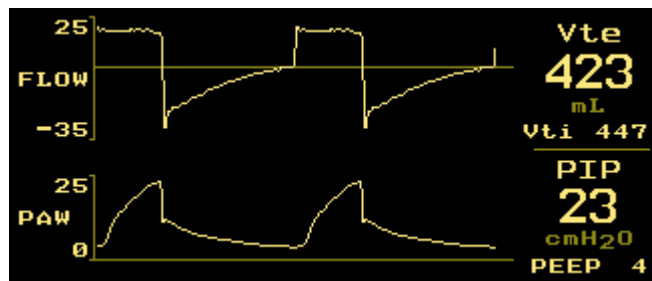

The **FLOW/PRESSURE** screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **FLOW/PRESSURE**.

- Information is updated in real time.
- The airway flow and pressure waveforms are automatically scaled.
- The waveform sweep speed is selectable in the **SETUP** menu.

### Numerics Screen

The **NUMERICS** screen display presents several monitored parameters together in one place.

ETCO <sub>2</sub> 34 mmHg	SpO <sub>2</sub> 97 %	MValv 5.2 L	PIP 23 cmH <sub>2</sub> O
RR 15	♥ 75	MV 6.2	MAP 9
VCO <sub>2</sub> 81 mL/m	Cdyn 27 mL/cmH <sub>2</sub> O	Raw 15 cm/(L/s)	PEEP 4 cmH <sub>2</sub> O

The **NUMERICS** screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **NUMERICS**.


- Information is updated in real time.
- The PEEP label is replaced with AUTO if Auto-PEEP (Intrinsic-PEEP) is detected.
- ETCO<sub>2</sub>, RR, SpO<sub>2</sub> and pulse rate values will flash if an alert limit is exceeded.

## Respiratory Screens (Optional)

NICO® offers four optional respiratory screens, displaying Respiratory Numerics, Flow Volume and Pressure Volume Loops, a Single Breath CO<sub>2</sub> Waveform, and the VCO<sub>2</sub>/MValv Trend. Press the **MENU** key and select **RESP SCREENS** from the **SELECT A SCREEN** menu by turning and then pressing the **KNOB**.

From the **ADD RESPIRATORY SCREENS** menu, highlight and then select which screens will appear in the base monitoring mode by turning and then pressing the **KNOB**; screen choices are dimmed until selected. When enabled, selected screens can be displayed by turning the **KNOB** while viewing any monitoring screen (See “KNOB Selectable Respiratory Screens” on page 12).

ADD RESPIRATORY SCREENS	
RESP NUMERICS	SELECTED RESPIRATORY SCREENS CAN BE ACCESSED BY KNOB
FLOW/PRES/VOL	
SINGLE BREATH CO <sub>2</sub>	
VCO <sub>2</sub> /MValv TREND	
EXIT	

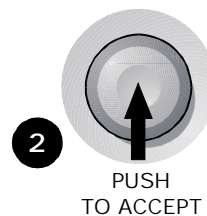

**Respiratory Numerics** The **RESPIRATORY NUMERICS** screen presents several monitored parameters together in one place.


Vte <b>423</b> mL	Vtalv <b>351</b> mL	MValv <b>5.2</b> L	Vd/Vt <b>0.52</b> @ 8:29
Vti 417	VdAW 71	MV 6.2	Vdalv 150
VCO <sub>2</sub> <b>160</b> mL/m	ETCO <sub>2</sub> <b>34</b> mmHg	PeCO <sub>2</sub> <b>23</b> mmHg	RSBI <b>0</b> br/m/L

The **RESPIRATORY NUMERICS** screen can be displayed, when enabled, by turning the **KNOB** while viewing any monitoring screen.

- Information is updated in real time.
- ETCO<sub>2</sub> value will flash if an alert limit is exceeded.

### Flow Volume and Pressure Volume Loops

The Flow Volume and Pressure Volume Loops screen displays a flow versus volume loop and a volume versus peak airway pressure loop.

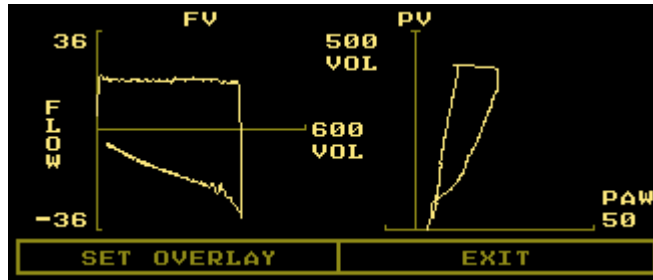


The Flow Volume and Pressure Volume Loops screen can be displayed, when enabled, by turning the **KNOB** while viewing any monitoring screen.

- The flow-volume loop is plotted in a clockwise direction, comparing flow versus volume for a single patient breath and providing information regarding the condition of the airways.
- The pressure-volume loop is plotted in a counter-clockwise direction; the slope from the beginning of inspiration to the end of inspiration depicts compliance while the width of the loop references resistance.
- The curve is automatically scaled to fit the display area.

### Set Overlay

The Set Overlay function “freezes” a single breath on the **FLOW/PRES/VOL** screen as a template. Subsequent loops will be drawn over that breath.

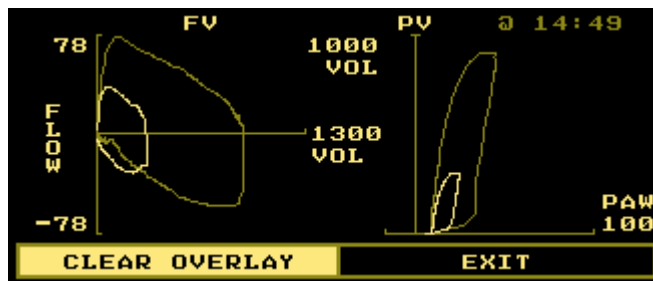

- Press the **KNOB** to highlight **SET OVERLAY** and again to select the breath.
- Press the **KNOB** to highlight **CLEAR OVERLAY** and again to erase the selected breath.
- The time stamp above the template wave corresponds to when the overlay was set.



### Single Breath CO<sub>2</sub>

The **SINGLE BREATH CO<sub>2</sub>** screen displays the CO<sub>2</sub> waveform for a single patient breath, as well as providing numeric displays of VCO<sub>2</sub>, Vte, ETCO<sub>2</sub>, Vtalv, Vti, and VdAw.



SBCO<sub>2</sub> presents the exhaled concentration of CO<sub>2</sub> versus tidal volume during a single expiration; useful for understanding the ventilation/perfusion relationship. It allows the clinician to detect relative changes in CO<sub>2</sub> production, deadspace, and effective ventilation by observing the shape of the graph.

- The curve is automatically scaled to fit the display area.

### Set Overlay


The Set Overlay function “freezes” a breath on the **SINGLE BREATH CO<sub>2</sub>** screen as a template. Subsequent waveforms will be drawn over that breath for comparison.



- Press the **NOB** to highlight **SET OVERLAY** and again to select the breath.
- Press the **NOB** to highlight **CLEAR OVERLAY** and again to erase the selected breath.
- The time stamp above the template wave corresponds to when the overlay was set.

### VCO<sub>2</sub>/MValv Trend

The **VCO<sub>2</sub>/MVALV TREND** screen displays trends for alveolar minute ventilation and CO<sub>2</sub> elimination. Monitoring spontaneous versus mechanical alveolar ventilation along with CO<sub>2</sub> elimination provides information on continued success or impending failure when weaning a patient off a ventilator.

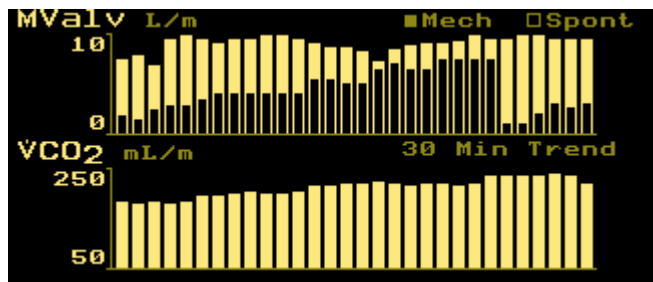

- While viewing the **VCO<sub>2</sub>/MV ALV TREND** screen, each push of the knob will advance through the available 30 minute, and 2, 4, 8, and 12 hour trend displays.
- Data is plotted from left (oldest data) to right (newest data). Once the display is filled, data shifts left so that the oldest data point on the left is pushed out to make room for the newest point entering from the right.


- Each bar on the trend represents the average MValv or VCO<sub>2</sub> value over a time period. The time periods are: 1 minute average for the 30 minute trend, 4 minute average for the 2 hour trend, 8 minutes for the 4 hour trend, 16 minutes for the 8 hour trend, and 24 minutes for the 12 hour trend.

### Systemic Vascular Resistance (SVR) Calculation

The **SVR CALCULATION** screen displays the Systemic Vascular Resistance formula and allows for entering the MBP, CVP and C.O. values for calculating the SVR value. The screen also displays the current time and last SVR value.

To calculate the SVR value:

- 1 Press the **MENU** key. The **SELECT A SCREEN** menu appears.
- 2 Highlight and then select **SVR CALCULATION** by turning and then pressing the knob.
- 3 The **SVR CALCULATION** screen is displayed.
- 4 Press then turn the **KNOB** to adjust the setting.
  - **MBP** - Mean Blood Pressure (25 to 300 mmHg).
  - **CVP** - Central Venous Pressure (-9 to 25 mmHg).
  - **C.O.** - Cardiac Output (0.5 to 19.9 L/m). This value will correspond to the displayed value (**AVERAGED** or **FAST**) and its value can be manually changed. If manually changed, the value will not be updated to reflect the displayed value.
  - **SVR** - Systemic Vascular Resistance (0 to 5000 dynes sec/cm<sup>5</sup>). Appears as "----" until the **MBP** value is entered.
  - **80** - A constant factor used to convert from Wood to VRU units.

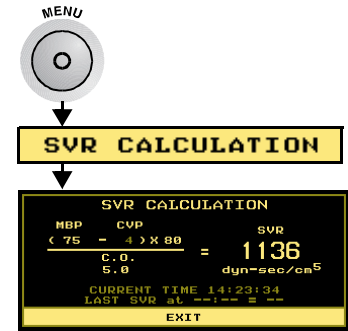

$$\frac{(\text{MBP} - \text{CVP}) \times 80}{\text{C.O.}} = \text{SVR} \text{ dyn-sec/cm}^5$$

- 5 Press the **KNOB** to accept the displayed value; turn to select the next setting.

### Tabular Data Screen

The **TABULAR DATA** screen displays the data collected for all NICO<sup>®</sup> parameters, in three-minute increments, in a table format.

The parameter displayed in each column will vary, depending on the setting chosen by the user.

- 1 Highlight and select the desired column by turning then pressing the **KNOB**.
- 2 Turn the **KNOB** to advance through all available NICO<sup>®</sup> parameters.
- 3 Press the **KNOB** to accept the displayed parameter; turn to select the next column.

TIME	CO-a	CI	SpO <sub>2</sub>	ETCO <sub>2</sub>	VCO <sub>2</sub>
11:40	6.1	2.4	97	35	156
11:43	5.7	2.3	97	35	156
11:46	5.7	2.3	97	35	156
11:47	5.3	2.1	97	39	74
11:48	6.1	2.1	97	34	157
11:52	5.7	2.3	97	35	156
11:55	5.7	2.3	97	35	156
EXIT					

### Time Column

With the **TIME** column selected (flashing), turn the **KNOB** to display NICO<sup>®</sup> parameters collected since the beginning of the monitoring session; the most recent records are at the bottom of the table. Three arrows ↓↓↓ will display at the bottom of the column when more records are available.

### Set Alerts Screen

The **SET ALERTS** screen displays the current patient values as well as the LOW and HIGH alert limit settings for various parameters. This screen also allows adjustment of alert limits, audible alert volume, and of the “No Respiration” alert. See “Alerts” on page 56 for details.

SET ALERTS	CURRENT	LOW	HIGH
CO (L/m)	6.5	🔔 OFF	🔔 OFF
ETCO <sub>2</sub> (mmHg)	34	🔔 OFF	🔔 OFF
SpO <sub>2</sub> (%)	97	🔔 OFF	🔔 OFF
RR (br/m)	15	🔔 OFF	🔔 OFF
♥ (bpm)	75	🔔 OFF	🔔 OFF
AUDIO: 10		🔔 NO RESP: 20s	
AUTO LIMITS		EXIT	

The **SET ALERTS** screen can only be displayed by pressing the **MENU** key and selecting **SET ALERTS**.

- **CURRENT** - the current patient value for a parameter, shown in real time  
NOTE: The **CURRENT** value will flash when a high or low alert limit is exceeded.
- **LOW** and **HIGH** - the values below and above which an alert will be generated
- 🔔 (Bell with slash) - the audible alert is disabled for this limit
- 🔔 (Bell) - the audible alert is enabled for this limit
- **AUDIO** - set the audible alert volume level
- **AUTO LIMITS** - have the NICO® monitor bracket alerts around current patient values
- **NO RESP:** - set the No Respiration alert delay timer

### Setup Screen

The **SETUP** screen allows the user to perform certain functions, such as a CO<sub>2</sub> Zero, erase trend data, change items like the waveform trace sweep speed and to set the time and date. Once set, the NICO® monitor will use these settings until again changed by the user.

To display the **SETUP** screen:

- 1 Press the **MENU** key to activate **SELECT A SCREEN**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 2 Highlight and then select **SETUP** by turning and then pressing the **KNOB**.
- 3 The **SETUP** screen is displayed.
- 4 Again, turn the **KNOB** to highlight an item and push the **KNOB** to select it.



SETUP	
CO <sub>2</sub> ZERO NOW	PULSE BEEP: OFF
PURGE NOW	SWEEP: MEDIUM
IABP MODE: OFF	C.O. MODE: AVERAGE
ERASE TRENDS	SET SPON THRESH
INPUT/OUTPUT	SET TIME & DATE
EXIT	

Items within the **SETUP** screen are described below:

Label	Settings/Range	Description
CO <sub>2</sub> ZERO NOW	Start or Cancel (default: Start)	Displays the <b>CO<sub>2</sub> ZERO NOW</b> screen. Place the CAPNOSTAT® CO <sub>2</sub> Sensor onto a clean and dry adapter. Place the adapter in room air and away from all sources of CO <sub>2</sub> . Select <b>START</b> to begin a CO <sub>2</sub> Zero or <b>CANCEL</b> to exit the selection and return to the <b>SETUP</b> menu.
PURGE NOW	n/a	The system immediately purges the NICO Sensor™ tubing. No messages are displayed. The Purge takes approximately 8 seconds. The flow and pressure waveform traces will return to zero during this period.


Label	Settings/Range	Description
IABP Intra-Aorta Balloon Pump	ON or OFF	Displays the <b>SET SpO<sub>2</sub> IABP MODE</b> screen. Turn the <b>KNOB</b> to select <b>OFF</b> (default setting) or <b>ON</b> to turn off the validator algorithm so that all pulsatile data, including the normally rejected artifact generated by the IABP, are allowed to influence the SpO <sub>2</sub> and Pulse Rate calculations. Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
ERASE TRENDS	Yes or No (default: No)	Displays the <b>ERASE STORED TRENDS?</b> screen. Turn the <b>KNOB</b> to select <b>NO</b> (the default setting) or <b>YES</b> . Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
INPUT/OUTPUT	ANALOG OUT 1-4 ANALOG CAL. RS232-2/RS232-3	Turn the <b>KNOB</b> to select <b>ANALOG OUT 1</b> through <b>ANALOG OUT 4</b> , <b>ANALOG CAL.</b> , <b>RS232-2</b> or <b>R232-3</b> . Push the <b>KNOB</b> to accept your selection and assign it the desired output parameter.
PULSE BEEP	OFF and 1-10 (default: OFF)	Displays the <b>SET PULSE BEEP</b> screen. Turn the <b>KNOB</b> to select a volume (1-10 and OFF) for the audible tone to accompany each detected pulse beat. Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
SWEEP	Slow, Medium, Fast (default: Medium)	Displays the <b>SET SWEEP SPEED</b> screen. Turn the <b>KNOB</b> to select how quickly the CO <sub>2</sub> , Flow and Pressure waveforms sweep across the display (Plethysmogram is unaffected). Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
C.O. MODE	AVERAGE or FAST	Displays the <b>SET C.O. MODE</b> screen. Turn the <b>KNOB</b> to select <b>AVERAGE</b> (CO-a, CObar displayed), <b>FAST</b> (CO-f, CObar not displayed). Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
SET SPON THRESHOLD	0-50 cmH <sub>2</sub> O	Displays the <b>SET SPONTANEOUS THRESHOLD</b> screen. The spontaneous threshold is the airway pressure chosen to differentiate between a spontaneous (patient-initiated) breath and a mechanical (ventilator) breath. Turn the <b>KNOB</b> to adjust the setting indicated by a dashed line. For optimal setting of the spontaneous threshold, the dashed line should be above any spontaneous breaths and below the mechanical breaths. Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
SET TIME & DATE	hh:mm dd mmm yyyy	Displays the <b>SET TIME / DATE</b> screen. Turn the <b>KNOB</b> to highlight the portion of the time/date to change. Push the <b>KNOB</b> to select that item—it begins to flash. Turn the <b>KNOB</b> to adjust the flashing item, and when correctly set, push the <b>KNOB</b> again to accept the value. Repeat for the other time/date entries. Finally, turn the <b>KNOB</b> to highlight <b>EXIT</b> and push the <b>KNOB</b> to return to the <b>SETUP</b> menu.

## Respiratory Mechanics mode

This section provides an overview of the various monitoring and setup screens in Respiratory Mechanics mode.

### CO<sub>2</sub> and SpO<sub>2</sub> Waveform Screen

The CO<sub>2</sub> and SpO<sub>2</sub> Waveform Screen plots the capnogram and plethysmogram signals.

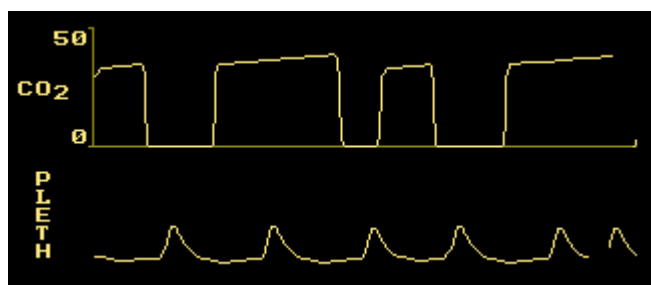


The  $\text{CO}_2/\text{SpO}_2$  screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting  $\text{CO}_2/\text{SpO}_2$ .

- Information is updated in real time.
- The capnogram and plethysmogram are automatically scaled.
- The capnogram sweep speed is selectable in the **SETUP** menu.

### Flow and Pressure Waveform Screen

The Flow and Pressure Waveform screen plots the airway flow and pressure signals over time.

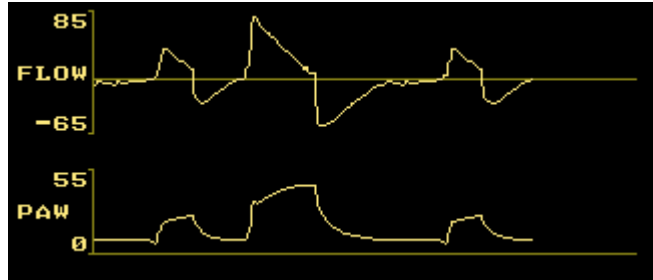


The **FLOW/PRESSURE** screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **FLOW/PRESSURE**.

- Information is updated in real time.
- The airway flow and pressure waveforms are automatically scaled.
- The waveform sweep speed is selectable in the **SETUP** menu.

### Flow Volume and Pressure Volume Loops screen

The Flow Volume and Pressure Volume Loops screen displays a flow versus volume loop and a volume versus peak airway pressure loop.

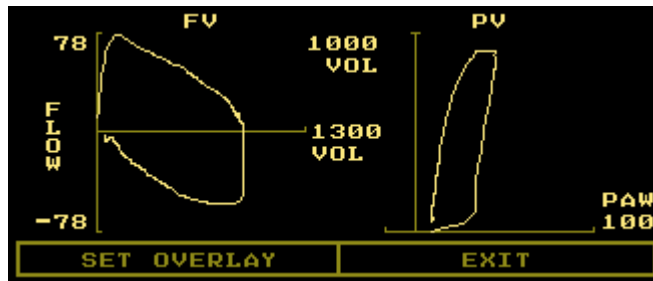


The **FLOW/PRES/VOL** screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **FLOW/PRES/VOL**.

- The flow-volume loop is plotted in a clockwise direction, comparing flow versus volume for a single patient breath and providing information regarding the condition of the airways.
- The pressure-volume loop is plotted in a counter-clockwise direction; the slope from the beginning of inspiration to the end of inspiration depicts compliance while the width of the loop references resistance.
- The curve is automatically scaled to fit the display area.

### Set Overlay


The Set Overlay function “freezes” a single breath on the **FLOW/PRES/VOL** screen as a template. Subsequent loops will be drawn over that breath.

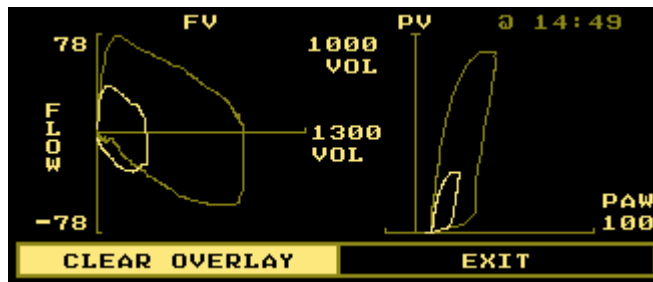

- Press the **KNOB** to highlight **SET OVERLAY** and again to select the breath.
- Press the **KNOB** to highlight **CLEAR OVERLAY** and again to erase the selected breath.
- The time stamp above the template waveform corresponds to when the overlay was set.

### VCO<sub>2</sub>/MValv Trend Screen


The **VCO<sub>2</sub>/MVALV** Trend screen displays trends for alveolar minute ventilation and CO<sub>2</sub> elimination. Monitoring spontaneous versus mechanical alveolar ventilation along with CO<sub>2</sub> elimination provides information on continued success or impending failure when weaning a patient off a ventilator.

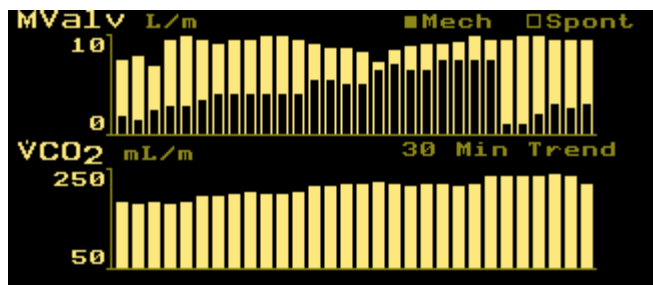

- While viewing the **VCO<sub>2</sub>/MVALV** Trend screen, each push of the knob will advance through the available 30 minute, and 2, 4, 8, and 12 hour trend displays.
- Data is plotted from left (oldest data) to right (newest data). Once the display is filled, data shifts left so that the oldest data point on the left is pushed out to make room for the newest point entering from the right.
- Each bar on the trend represents the average MValv or VCO<sub>2</sub> value over a time period. The time periods are: 1 minute average for the 30 minute trend, 4 minute average for the 2 hour trend, 8 minutes for the 4 hour trend, 16 minutes for the 8 hour trend, and 24 minutes for the 12 hour trend.

### Vt/Vd Trend Screen

The **Vt/Vd TREND** screen displays trends for mechanical and spontaneous airway deadspace and alveolar tidal volume. Monitoring spontaneous versus mechanical deadspace along with expired tidal volume provides information on continued success or impending failure when weaning a patient off a ventilator.

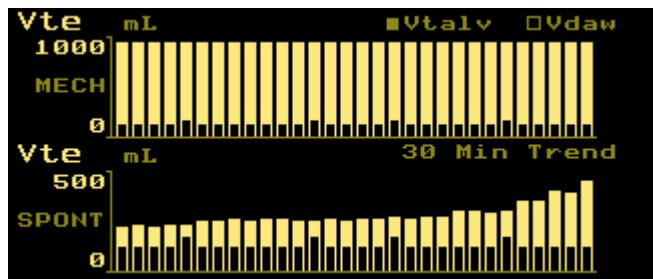

- While viewing the **Vt/Vd TREND** screen, each push of the knob will advance through the available 30 minute, and 2, 4, 8, and 12 hour trend displays.
- Data is plotted from left (oldest data) to right (newest data). Once the display is filled, data shifts left so that the oldest data point on the left is pushed out to make room for the newest point entering from the right.
- Each bar on the trend represents the average Vtalv or Vdaw value over a time period. The time periods are: 1 minute average for the 30 minute trend, 4 minute average for the 2 hour trend, 8 minutes for the 4 hour trend, 16 minutes for the 8 hour trend, and 24 minutes for the 12 hour trend.

### Single Breath CO<sub>2</sub> Screen

The **SINGLE BREATH CO<sub>2</sub>** screen displays the CO<sub>2</sub> waveform for a single patient breath, as well as providing numeric displays of Vd/Vt, Vtalv, Vdaw, PeCO<sub>2</sub>, and Vdalv.



SBCO<sub>2</sub> presents the exhaled concentration of CO<sub>2</sub> versus tidal volume during a single expiration; useful for understanding the ventilation/perfusion relationship. It allows the clinician to detect relative changes in CO<sub>2</sub> production, deadspace, and effective ventilation by observing the shape of the graph.

- The curve is automatically scaled to fit the display area.

### Set Overlay

The Set Overlay function “freezes” a breath on the **SINGLE BREATH CO<sub>2</sub>** screen as a template. Subsequent waveforms will be drawn over that breath.



- Press the **KNOB** to highlight **SET OVERLAY** and again to select the breath.
- Press the **KNOB** to highlight **CLEAR OVERLAY** and again to erase the selected breath.
- The time stamp above the template wave corresponds to when the overlay was set.

### Numerics Screen

The **NUMERICS** screen display presents several monitored parameters together in one place.

Vd/Vt <b>0.20</b> 12:38	Vtalv <b>546</b> mL	MValv <b>11.0</b> L	PIP <b>26</b> cmH <sub>2</sub> O
Vdalv 58	VdAW 64	MAP 16	Pplat --
PeCO <sub>2</sub> <b>32</b> mmHg	Cdyn <b>21</b> mL/cmH <sub>2</sub> O	Raw <b>10</b> cm/(L/s)	PEEP <b>10</b> cmH <sub>2</sub> O

The **NUMERICS** screen can be displayed by turning the **KNOB** while viewing any monitoring screen, or by pressing the **MENU** key and selecting **NUMERICS**.

- Information is updated in real time.
- The PEEP label is replaced with AUTO if Auto-PEEP (Intrinsic-PEEP) is detected.

### Tabular Data Screen

The **TABULAR DATA** screen displays the data collected for all NICO® parameters, in one-minute increments, in a table format. Cardiac output-related parameters will not be displayed and will remain dashed (--) and inactive in Respiratory Mechanics mode.

The parameter displayed in each column will vary, depending on the setting chosen by the user.

- 1 Highlight and select the desired column by turning then pressing the **KNOB**.
- 2 Turn the **KNOB** to advance through all available NICO® parameters.
- 3 Press the **KNOB** to accept the displayed parameter; turn to select the next column.

TIME	PIP	CO-a	SpO <sub>2</sub>	ETCO <sub>2</sub>	VCO <sub>2</sub>
11:12	32	--	97	39	430
11:13	32	--	97	39	430
11:14	32	--	97	39	430
11:15	32	--	97	39	430
11:16	32	--	97	39	430
11:17	32	--	97	39	430
11:18	30	--	97	39	430
EXIT					

### Time Column

With the **TIME** column selected (flashing), turn the **KNOB** to display NICO® parameters collected since the beginning of the monitoring session; the most recent records are at the bottom of the table. Three arrows ↓↓↓ will display at the bottom of the column when more records are available.



### Set Alerts Screen

The **SET ALERTS** screen displays the current patient values as well as the LOW and HIGH alert limit settings for various parameters (PIP is HIGH only). This screen also allows adjustment of alert limits, audible alert volume, and of the "No Respiration" alert. See "Alerts" on page 56 for details.

SET ALERTS	CURRENT	LOW	HIGH
PIP cmH <sub>2</sub> O	46		🔊 120
ETCO <sub>2</sub> (mmHg)	39	🔊 1	🔊 150
SpO <sub>2</sub> (%)	97	🔊 50	🔊 100
RR (br/m)	20	🔊 3	🔊 150
📈 (bpm)	75	🔊 31	🔊 249
AUDIO: 4	🔊 NO RESP: 20s		
AUTO LIMITS	EXIT		



The **SET ALERTS** screen can only be displayed by pressing the **MENU** key and selecting **SET ALERTS**.


- **CURRENT** - the current patient value for a parameter, shown in real time  
NOTE: The **CURRENT** value will flash when a high or low alert limit is exceeded.
- **LOW** and **HIGH** - the values below and above which an alert will be generated
-  (Bell with slash) - the audible alert is disabled for this limit
-  (Bell) - the audible alert is enabled for this limit
- **AUDIO** - set the audible alert volume level
- **AUTO LIMITS** - have the NICO® monitor bracket alerts around current patient values
- **NO RESP:** - set the No Respiration alert delay timer

## Setup Screen

The **SETUP** screen allows the user to perform certain functions, such as a CO<sub>2</sub> Zero, erase trend data, change items like the waveform trace sweep speed and to set the time and date. Once set, the NICO® monitor will use these settings until again changed by the user.

To display the **SETUP** screen:

- 1 Press the **MENU** key to activate **SELECT A SCREEN**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 2 Highlight and then select **SETUP** by turning and then pressing the **KNOB**.
- 3 The **SETUP** screen is displayed.
- 4 Again, turn the **KNOB** to highlight an item and push the **KNOB** to select it.

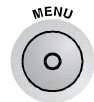

SETUP	
CO <sub>2</sub> ZERO NOW	PULSE BEEP: OFF
PURGE NOW	SWEEP: MEDIUM
IABP MODE: OFF	SET SPON THRESH
ERASE TRENDS	SET TIME & DATE
INPUT/OUTPUT	
EXIT	

Items within the **SETUP** screen are described below:

Label	Settings/Range	Description
CO <sub>2</sub> ZERO NOW	Start or Cancel (default: Start)	Displays the <b>CO<sub>2</sub> ZERO NOW</b> screen. Place the CAPNOSTAT® CO <sub>2</sub> Sensor onto a clean and dry adapter. Place the adapter in room air and away from all sources of CO <sub>2</sub> . Select <b>START</b> to begin a CO <sub>2</sub> Zero or <b>CANCEL</b> to exit the selection and return to the <b>SETUP</b> menu.
PURGE NOW	n/a	The system immediately purges the CO <sub>2</sub> /Flow Sensor tubing. No messages are displayed. The Purge takes approximately 8 seconds. The flow and pressure waveform traces will return to zero during this period.
IABP Intra-Aorta Balloon Pump	ON or OFF	Displays the <b>SET SpO<sub>2</sub> IABP MODE</b> screen. Turn the <b>KNOB</b> to select <b>OFF</b> (default setting) or <b>ON</b> to turn off the validator algorithm so that all pulsatile data, including the normally rejected artifact generated by the IABP, are allowed to influence the SpO <sub>2</sub> and Pulse Rate calculations. Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
ERASE TRENDS	Yes or No (default: No)	Displays the <b>ERASE STORED TRENDS?</b> screen. Turn the <b>KNOB</b> to select <b>NO</b> (the default setting) or <b>YES</b> . Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
INPUT/OUTPUT	ANALOG OUT 1-4 ANALOG CAL. RS232-2/RS232-3	Turn the <b>KNOB</b> to select <b>ANALOG OUT 1</b> through <b>ANALOG OUT 4</b> , <b>ANALOG CAL.</b> , <b>RS232-2</b> or <b>R232-3</b> . Push the <b>KNOB</b> to accept your selection and assign it the desired output parameter.

Label	Settings/Range	Description
PULSE BEEP	OFF and 1-10 (default: OFF)	Displays the <b>SET PULSE BEEP</b> screen. Turn the <b>KNOB</b> to select a volume (1-10 and OFF) for the audible tone to accompany each detected pulse beat. Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
SWEEP	Slow, Medium, Fast (default: Medium)	Displays the <b>SET SWEEP SPEED</b> screen. Turn the <b>KNOB</b> to select how quickly the CO <sub>2</sub> , Flow and Pressure waveforms sweep across the display (Plethysmogram is unaffected). Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
SET SPON THRESHOLD	0-50 cmH <sub>2</sub> O	Displays the <b>SET SPONTANEOUS THRESHOLD</b> screen. The spontaneous threshold is the airway pressure chosen to help NICO® differentiate between a spontaneous breath and a mechanical breath. Turn the <b>KNOB</b> to adjust the setting indicated by a dashed line. For optimal setting of the spontaneous threshold, the dashed line should be above any spontaneous breaths and below the mechanical breaths. Push the <b>KNOB</b> to accept your selection and return to the <b>SETUP</b> menu.
SET TIME & DATE	hh:mm dd mmm yyyy	Displays the <b>SET TIME / DATE</b> screen. Turn the <b>KNOB</b> to highlight the portion of the time/date to change. Push the <b>KNOB</b> to select that item—it begins to flash. Turn the <b>KNOB</b> to adjust the flashing item, and when correctly set, push the <b>KNOB</b> again to accept the value. Repeat for the other time/date entries. Finally, turn the <b>KNOB</b> to highlight <b>EXIT</b> and push the <b>KNOB</b> to return to the <b>SETUP</b> menu.



## Notes on Patient Monitoring

---

---

### Automatic Purging

A double lumen connecting line (tubing) connects the NICO® Flow and CO<sub>2</sub>/Flow sensors to the NICO® monitor. The NICO® monitor includes an automatic and manual purge feature which provides a flow rate of room air to keep the sensor tubing free from water condensation and patient secretions. This feature is available for the adult, pediatric, and neonatal modes.

#### Adult mode

The system automatically purges the sensor tubing every 10 minutes or less, depending on system conditions. In adult mode, the system will purge both sides of the line, one at a time, during each purge cycle. The higher the pressure, the more frequent the purging. This action anticipates increased moisture migration into the sensor tubing due to the increase in circuit pressure.

#### Neonatal and Pediatric modes

The automatic purge cycle used in the neonatal or pediatric mode is fixed at every 3 minutes regardless of circuit pressure. Only one side of the sensor tubing will be purged during each purge cycle. The purge will only occur during the exhalation portion of the ventilator cycle, regardless of exhalation time.

Unlike the adult purge mode, the neonatal or pediatric purge mode does not use the full force of the internal pump, but rather pressurizes an internal reservoir which is used for the purge. This minimizes the pressure delivered to the ventilator circuit, but does deliver a sufficient pressure to purge the sensor tubing.

### Manual Purging

Occasionally, purging may be required in between the automatic purge cycle. The manual purge may be used as often as needed, and may be used at all times **except** when a very low battery condition exists. During very low battery conditions, automatic and manual purging is not allowed.


To manually purge:

- 1 Press the **MENU** key.
- 1 Turn and press the **KNOB** to select **SETUP**.
- 1 Turn and press the **KNOB** to select **PURGE NOW**.
- 2 The purge cycle will begin:

In adult mode, the system will purge both sides of the line, one at a time, during each purge cycle.

In pediatric and neonatal modes, only one side of the sensor tubing will be purged during each purge cycle. The purge is synchronized to the exhalation phase of the ventilator cycle and will not exceed exhalation time.


A purge can also be initiated upon request if the flow waveform appears as though the lines are partially occluded (see example below) and the purge did not initiate automatically. To initiate a purge, see Manual Purging, above.

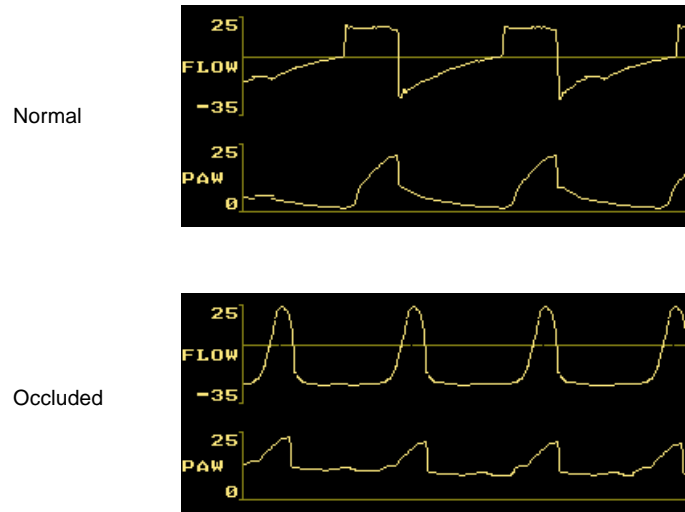

**NOTE:**

If the purge does not sufficiently clear the flow tubing lines, the flow sensor should be changed.

Below is an example of a waveform that exhibits a small leak in the breathing circuit. Replace the sensor; if problem persists, refer monitor to qualified service personnel.

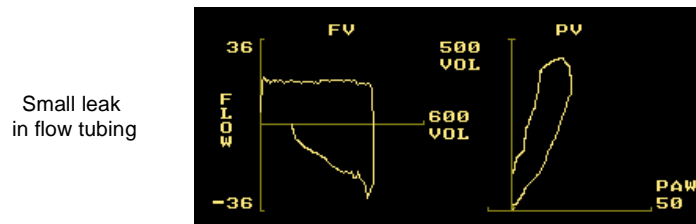

## Intra-Aortic Balloon Pump

NICO® uses advanced signal processing algorithms to distinguish valid pulsatile signals from signals generated by motion or other artifact. Motion artifact, very common in all but heavily sedated patients, can swamp the true pulsatile signal or distort it enough to produce significant errors in the SpO<sub>2</sub> and Pulse Rate calculations. The validator algorithms reject distorted plethysmographic signals or those that lack a regular rhythmic pattern; therefore, only valid (pulsatile) signals are allowed to affect the monitor's SpO<sub>2</sub> and Pulse Rate calculations. Rare conditions exist where the pulsatile waveform truly is distorted and lacks a fixed rhythm, specifically during use of an Intra-Aortic Pump (IABP).

During IABP procedures the pulsatile signal can be massively distorted without affecting the patient's SpO<sub>2</sub>. In order to accommodate these IABP procedures without compromising the monitor's superior artifact rejection algorithm, IABP MODE is available. IABP MODE allows the

user to turn off the validator algorithm so that all pulsatile data are allowed to influence the SpO<sub>2</sub> and Pulse Rate calculations.

<b>NOTE</b>	
<p>With IABP MODE turned ON, the clinician must exercise prudence in assessing the validity of the SpO<sub>2</sub> and Pulse Rate displays because any motion or other artifact—not just that associated with the IABP—can have a significant affect on the SpO<sub>2</sub> and Pulse Rate calculations.</p>	
<p>While in IABP MODE, the displayed Pulse Rate reflects true pulsatile signal-heart rate plus the IABP ratio (e.g. #1: heart rate = 120 bpm, IABP ratio = 1: 1, then displayed Pulse Rate should be 120 + (120/1)=240 beats/min. e.g. #2: heart rate = 120 bpm, IABP = 1: 3, then displayed Pulse Rate should be 120 + (120/3) = 180 beats/min). When in IABP MODE the Pulse Rate can be affected by motion or other artifact; the accuracy of the Pulse Rate can usually be used as an indicator of the quality of the SpO<sub>2</sub> display.</p>	

## Configuration Menu

Simultaneously press and hold the **MENU** and **DATA ENTRY** keys for 3 seconds to access the Configuration Menu. Turn and press the **KNOB** to adjust and accept settings.

Parameter	Range/Units	Description	Factory Default
CO <sub>2</sub> UNITS	mmHg, %, kPa	Select the desired units for the capnogram, PeCO <sub>2</sub> , ETCO <sub>2</sub> and PaCO <sub>2</sub> values.	mmHg
ETCO <sub>2</sub> AVG	10 sec, 20 sec, 1 Breath	Select the interval from which the displayed value of end tidal CO <sub>2</sub> (ETCO <sub>2</sub> ) is calculated.	10 sec
VCO <sub>2</sub> AVG	10 min, 7 min, 5 min, 3 min, 1 min, 8 Breath	Select the averaging time for the displayed value of CO <sub>2</sub> elimination (VCO <sub>2</sub> ).	1 min
ALLOW AUDIO OFF	N/A	When the ALLOW AUDIO OFF setting is set to YES, the audible alert tone may be silenced permanently by pressing and holding the SILENCE key.  If the ALLOW AUDIO OFF is set to NO, the audible alert tone may not be silenced for more than 2 minutes.	Yes
LANGUAGE	N/A	Select the desired language	English

## Reference Handbooks

For a discussion on waveform interpretations, refer to the Novamatrix Reference Handbooks on capnography, respiratory mechanics, and pulse oximetry. Contact Novamatrix Customer Service or your local sales representative for more information.



# Alerts

This section describes NICO® monitor alerts.

## NICO® Alert Priorities

The NICO® monitor prioritizes alert notifications. This prioritization allows an alert condition for which an immediate user response is required to take precedence over a lesser alert for which a less urgent response is acceptable. Alert notifications may include on-screen display messages and audible tones, and may result from violations of parameter limit settings, or from monitor or sensor related error conditions.

### High Priority Alert

- Action: Immediate user response
- Audible: 3 consecutive tones, repeated every 5 seconds (if enabled)
- Visual: the **SILENCE** key indicator flashes red, and a screen message is displayed
- Example: **LOW C.O.**



### Medium Priority Alert

- Action: Prompt user response
- Audible: 2 consecutive tones, repeated every 10 seconds (if enabled)
- Visual: Screen message
- Example: **HIGH RESP**

### Low priority Alert

- Action: User awareness
- Audible: a single tone, repeated every 15 seconds (if enabled)
- Visual: Screen message
- Example: **EXPAND LOOP**

### Status Messages

- Action: Informational, no urgency
- Audible: none
- Visual: Screen message
- Example: **ALERTS OFF**

## Responding to Alert Audio

The **SILENCE** key is used to mute/prevent audible alerts. It also visually indicates the presence of a “High Priority Alert”. The Silence feature operates in two modes; a temporary “2 Minute Silence” mode and an “Audio Disabled” mode.



- 2 Minute Silence — Press and release to activate or deactivate the two minute silence. The key's icon illuminates amber when active and audible alerts will be muted for two minutes, after which the icon turns off and any active audible alert will sound.
- Audio Disabled — Press and hold for two seconds to prevent or allow any audible alerts. The key's icon illuminates and flashes amber to indicate that all audible alerts are being suppressed.
- High Priority Alerts — The **SILENCE** key's icon illuminates and flashes red to indicate High Priority Alert is active. The icon alternately flashes red and amber if the audio is disabled and a High Priority Alert is active.

## Parameter Limit Alerts

The NICO® monitor allows for the establishment of high and low limit alerts for Cardiac Output (C.O., Cardiac Output mode only), Positive Inspiratory Pressure (PIP, Respiratory Mechanics mode only), End-Tidal Carbon Dioxide (ETCO<sub>2</sub>), Oxygen Saturation (SpO<sub>2</sub>), Respiratory Rate (RR), and Pulse Rate (♥). These alerts provide a visual (and audible, if desired) indication to the user that a parameter has violated the limit settings

Cardiac Output mode			Respiratory Mechanics mode		
Parameter	Range	Units	Parameter	Range	Units
C.O.	OFF, 0.1 - 19.9	L/m	PIP	OFF, 5 - 120	cmH <sub>2</sub> O
ETCO <sub>2</sub>	OFF, 1 - 150	mmHg	ETCO <sub>2</sub>	OFF, 1 - 150	mmHg
	OFF, 0.1 - 19.9	kPa		OFF, 0.1 - 19.9	kPa
	OFF, 0.1 - 19.9	%		OFF, 0.1 - 19.9	%
SpO <sub>2</sub>	OFF, 50 - 100	%	SpO <sub>2</sub>	OFF, 50 - 100	%
RR	OFF, 3 - 150	br/min	RR	OFF, 3 - 150	br/min
♥ (Pulse Rate)	OFF, 31 - 249	bpm	♥ (Pulse Rate)	OFF, 31 - 249	bpm

- All alert limit values are retained each time the monitor is turned off.
- The user can select individual limit values or have the monitor automatically assign limit values based on current patient values with the **AUTO LIMITS** feature.
- Audible alarms can be enabled or disabled for each individual high and low limit setting.
- Limits cannot be adjusted to provide less than a 5 digit separation between the high and low limit values.
- Violating the 5-digit minimum separation constraint may cause previously set alerts to reset or turn off.

### NOTE

The message **ALERTS OFF** will appear if the monitor is powered up with all alert limits and audible alerts set to **OFF**.

To cancel the **ALERTS OFF** message, adjust any individual limit value and activate the audible alert (🔔). Both the alert limit and audible alert must be activated.

### View Alert Settings

To view the existing alert settings:

- 1 Press the **MENU** key. The **SELECT A SCREEN** menu appears.
- 2 Highlight and then select **SET ALERTS** by turning and then pressing the knob.
- 3 The **SET ALERTS** screen is displayed.

### Adjust Alert Limits

To adjust **LOW** or **HIGH** alert limit values:

- 1 Display the **SET ALERTS** screen. (Press **MENU**, then select **SET ALERTS**.)
- 2 Highlight and then select the **LOW** or **HIGH** limit you wish to adjust (PIP is **HIGH** only).
- 3 Turn the knob to adjust the limit, then press the knob to accept the displayed value.
- 4 If desired, enable the audio for the alerts which have been enabled. (See below.)

### Enable/Disable Audible Alerts


To enable or disable the audible alert tone(s) for a particular alert:

- 1 Display the **SET ALERTS** screen. (Press **MENU**, then select **SET ALERTS**.)
- 2 Highlight the bell icon for the alert you wish to adjust.
- 3 Press the knob to switch between audio enabled (🔔) and disabled (🔕) settings.

### Adjust Alert Audio Volume

To adjust the Alert Audio Volume:

- 1 Display the **SET ALERTS** screen. (Press **MENU**, then select **SET ALERTS**.)
- 2 Highlight and then select **AUDIO** by turning and then pressing the knob.
- 3 Turn the knob to adjust the volume setting.
  - Volume ranges from 1 to 10 (loudest).
  - A representative tone sounds at each setting.
- 4 Press the knob to accept the displayed value.

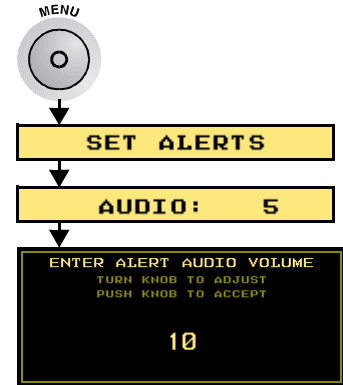

#### NOTE

Make sure that the audible alert volume is not set too low to be heard over ambient noise levels.

### Auto Alerts

Alert limits for **LOW** and **HIGH** C.O., PIP (**HIGH** only), ETCO<sub>2</sub>, SpO<sub>2</sub>, RR and Pulse rate can be set automatically for alerts which have been enabled. Alert limits are bracketed about the current patient values. Auto Alerts do not affect the status of the audible alert bell icon.

## No Respiration Alert

The NICO® monitor incorporates a No Respiration Alert (**NO RESP**). Once the monitor detects respiration, any loss of that signal starts the "No Resp" timer. A visual (and audible, if desired) alert occurs if the monitor does not detect a respiratory rate signal before the timer reaches its user set limit (20 seconds by default).



If a No Respiration alert occurs:

- The message **NO RESP: x:xx** is displayed (unless a higher priority alert is in progress).
- The counter shows the elapsed time (minutes:seconds) since the alert occurred.
- Press the **SILENCE** key to cancel the alert.
- The alert (message, audio and flashing **SILENCE** key) is automatically cancelled after ten minutes.

#### WARNING




The NICO® monitor is not intended to be used as an apnea monitor.

### Adjust the No Resp. Alert Limit

To adjust the No Respiration Alert Limit:

- 1 Display the **SET ALERTS** screen. (Press **MENU**, then select **SET ALERTS**.)
- 2 Highlight and then select **NO RESP** by turning and then pressing the knob.
- 3 Turn the knob to adjust the limit setting.
  - The limit timer is selectable from 10-60 seconds.
- 4 Press the knob to accept the displayed value.

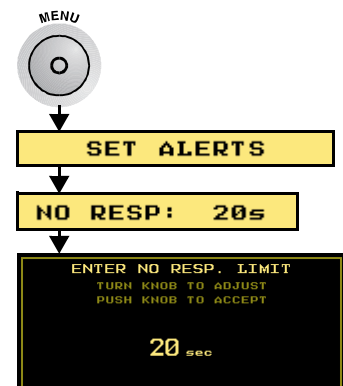






## NICO Sensor™

This section provides information regarding Disposable NICO Sensors™ and their use with the CAPNOSTAT® CO<sub>2</sub> Sensor and the NICO® monitor.

### Disposable NICO Sensors™

A NICO Sensor™ incorporates a rebreathing valve, NICO Loop™ (an adjustable rebreathing volume) and an adult CO<sub>2</sub>/Flow Sensor. It is disposable and intended for Single Patient Use only. The NICO Sensor™ is not for pediatric use.

A NICO Sensor™ may be connected to and removed from the NICO® monitor while the monitor is turned off or on.

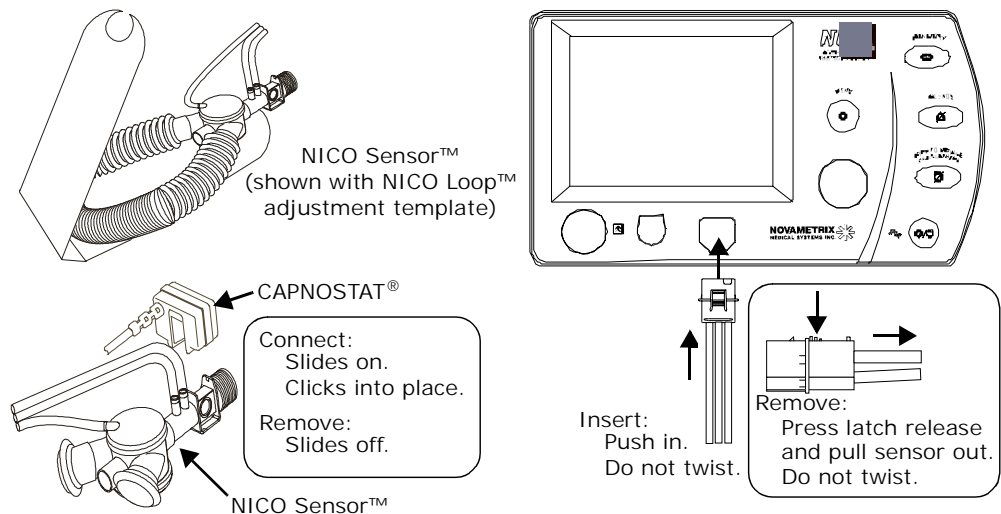

**WARNING**

The Disposable NICO Sensor™ is intended for SINGLE PATIENT USE ONLY (X).


Re-use, including disassembly, cleaning, disinfecting, sterilizing, and other efforts made in an attempt to re-use a NICO Sensor™, may compromise system performance and may cause a potential patient hazard. Performance is not guaranteed if reused.

Before use, verify the sensor is physically intact, with no broken or damaged parts. Do not use a broken or damaged sensor or one with wet, contaminated, or corroded connectors.

Do not allow the NICO Sensor™ to remain in the ventilator circuit when not connected to the NICO® monitor. A ventilator circuit leak will occur.

**CAUTION**

Connect only a Novamatrix NICO Sensor™, Catalog Number 8950-00, 8951-00 or 8952-00 to the NICO® monitor. Do not connect any other sensor in place of a NICO Sensor™.

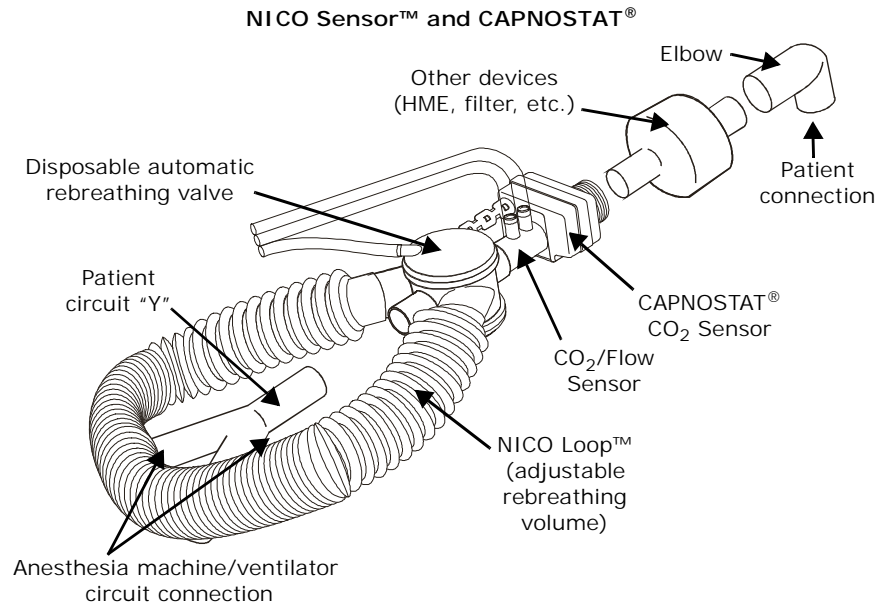


**Choosing a NICO Sensor™ size**

Choose a NICO Sensor™ size based on the tidal volume ranges listed below.

Sensor Size	Tidal Volume Range	Part Number
Small	For use with ventilator set tidal volumes of 200-500 ml	8950-00
Standard	For use with ventilator set tidal volumes of 400-1000 ml	8951-00
Large	For use with ventilator set tidal volumes of 750-1500 ml	8952-00

**WARNING**



 The NICO Sensor™ is not for pediatric use.

The NICO Sensor™ increases airway deadspace by 35 cc (minimum). At low tidal volumes, compensatory changes to ventilation protocol should be considered.



**Connecting a NICO Sensor™**

To connect a NICO Sensor™ into the ventilator circuit:

- 1 Remove a new NICO Sensor™ and loop adjustment template from its package. Open and inspect the sensor. Do not use if damaged.
- 2 Using the Initial Adjustment Template as a guide, adjust the NICO Loop™ to match the ventilator's tidal volume setting. Discard the template.

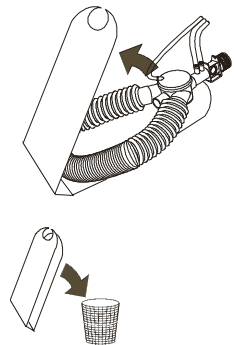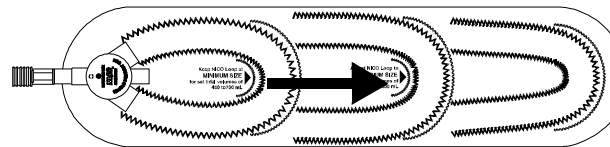


- 3 Insert the connector on the NICO Sensor™ into NICO® monitor's front panel. Verify a **NICO SENSOR IDENTIFIED** message is briefly displayed.
- 4 Attach the CAPNOSTAT® to the NICO Sensor™ and connect the NICO Sensor™ to ventilator circuit. For optimal results, place the NICO Sensor™ into the ventilator circuit between the endotracheal tube and the ventilator circuit wye.
  - Place other devices (HME, filters, etc.) between NICO Sensor™ and the patient connection.
  - Placement of a sidestream gas analyzer sampling port between the NICO Sensor™ and the patient connection may reduce NICO® accuracy at low tidal volumes.
  - Sidestream or mainstream gas analyzers placed between the NICO Sensor™ and the patient circuit "Y" may report diluted readings during the rebreathing phase of the NICO® cycle.
  - Place the sensor so that the triple lumen tubing lines exit from the top of the sensor (to help keep them clear and dry).
  - Keep the sensor clear of accumulations by proper circuit maintenance.
- 5 To begin monitoring, press the **STOP / CONTINUE REBREATHING** key.

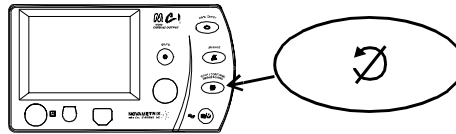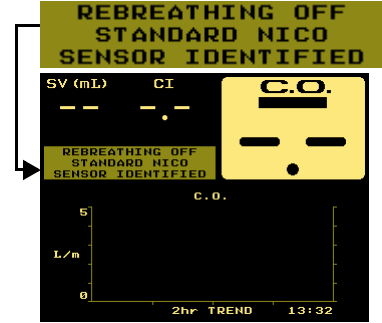

- To enhance accuracy, enter the respiratory gas composition and the arterial blood gas values whenever possible (by pressing the **DATA ENTRY** key).


### Status Messages

While monitoring, if the **EXPAND** or **RETRACT LOOP** message appears, adjust the loop 3-6 inches (8-15 cm) or until the message is removed.

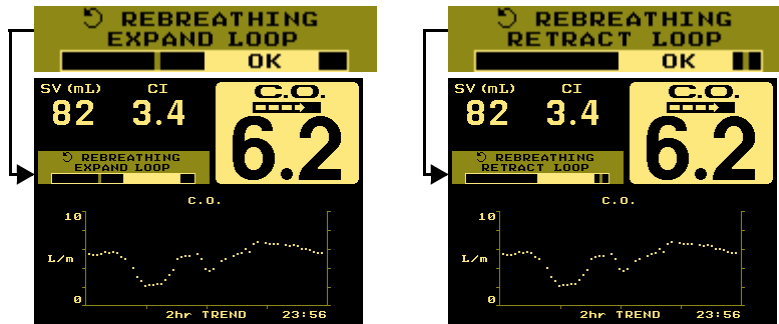

If the **EXPAND** or **RETRACT LOOP** message appears for more than three consecutive rebreathing cycles, and resizing the NICO Loop™ was not effective, the NICO® monitor will suggest a different sized sensor, with a larger or smaller loop, to correct the condition.


- Small NICO Sensor™: the monitor displays **CONSIDER USING A STANDARD NICO SENSOR** when ventilator set tidal volume is greater than 500 ml.
- Standard NICO Sensor™: the monitor displays **CONSIDER USING A SMALL NICO SENSOR** when ventilator set tidal volume is less than 300 ml.
- Standard NICO Sensor™: the monitor displays **CONSIDER USING A LARGE NICO SENSOR** when ventilator set tidal volume is greater than 1000 ml.
- Large NICO Sensor™: the monitor displays **CONSIDER USING A STANDARD NICO SENSOR** when ventilator set tidal volume is less than 1000 ml.

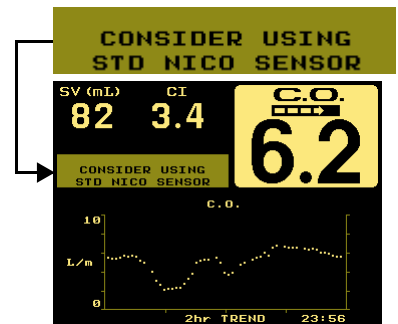

## CO<sub>2</sub>/Flow Sensors

This section details the use of combined CO<sub>2</sub>/Flow sensors with the NICO<sup>®</sup> monitor. This section also explains how to connect a sensor to the monitor, and how to apply the sensor to the patient. NICO<sup>®</sup> automatically operates in Respiratory Mechanics mode whenever a CO<sub>2</sub>/Flow sensor is connected to the monitor. Cardiac output and associated parameters are no longer available.

Only NICO<sup>®</sup> CO<sub>2</sub>/Flow sensors are compatible with the NICO<sup>®</sup> monitor. Sensors may be connected and removed with the monitor off or on.

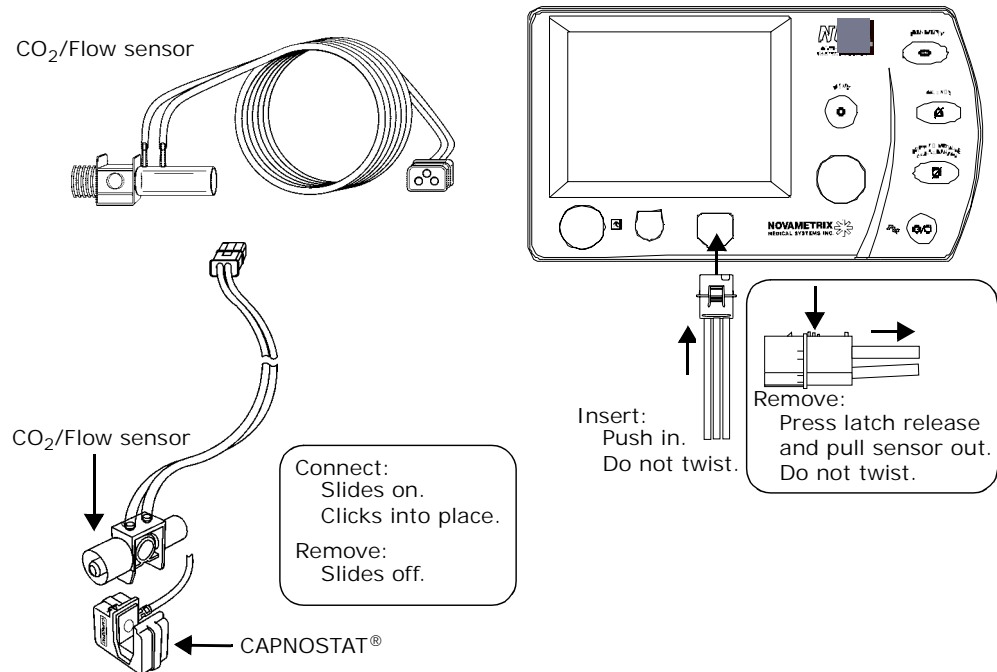


### WARNING

The Disposable CO<sub>2</sub>/Flow sensor is intended for SINGLE PATIENT USE ONLY (X).

Re-use, including disassembly, cleaning, disinfecting, sterilizing, and other efforts made in an attempt to re-use a CO<sub>2</sub>/Flow sensor, may compromise system performance and may cause a potential patient hazard. Performance is not guaranteed if reused.




Before use, verify the sensor is physically intact, with no broken or damaged parts. Do not use a broken or damaged sensor or one with wet, contaminated, or corroded connectors.

Do not allow the CO<sub>2</sub>/Flow sensor to remain in the ventilator circuit when not connected to the NICO<sup>®</sup> monitor. A ventilator circuit leak will occur.

	<b>CAUTION</b>
Connect only Novamatrix NICO® CO <sub>2</sub> /Flow sensors, Catalog Number 9765-00, 9766-00, or 9767-00, to the NICO® monitor. Do not connect any other sensor in place of a NICO CO <sub>2</sub> /Flow sensor.	

## Choosing a CO<sub>2</sub>/Flow Sensor

Select the appropriate combined CO<sub>2</sub>/Flow sensor based on endotracheal tube size (ETT), volume and flow rate. Ranges for each sensor are listed in the chart below:

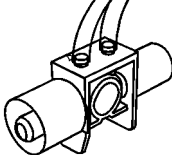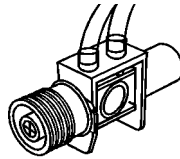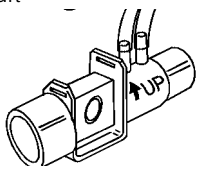
Sensor	Range			
	ETT (mm)	Volume (ml)	Flow	Deadspace
Neonatal 	2.5-4.0	1-100	0.25-25 LPM 4-417 ml/sec	Less than 1 ml
Pediatric 	3.5-6.0	30-400	0.5-100 LPM 8-1600 ml/sec	Less than 4 ml
Adult 	>5.5	200-3000	2.0-180 LPM 33-3000 ml/sec	Less than 8.5 ml

## Connecting a CO<sub>2</sub>/Flow Sensor

To connect a combined CO<sub>2</sub>/Flow sensor:

- 1 Remove a new CO<sub>2</sub>/Flow sensor from its package. Open and inspect the sensor. Do not use if damaged.
- 2 Insert the connector on the CO<sub>2</sub>/Flow sensor into NICO® monitor's front panel. Verify a **CO<sub>2</sub>/FLOW SENSOR IDENTIFIED** message is briefly displayed.
- 3 Attach the CAPNOSTAT® to the CO<sub>2</sub>/Flow sensor. The sensor will "click" when properly seated.

When a CO<sub>2</sub>/Flow sensor is replaced, the message **CHECK/CHANGE AIRWAY ADAPTER** may appear, this means a CO<sub>2</sub> zero is required to help the NICO® monitor "learn" the specific optical characteristics of the CO<sub>2</sub>/Flow sensor in use.

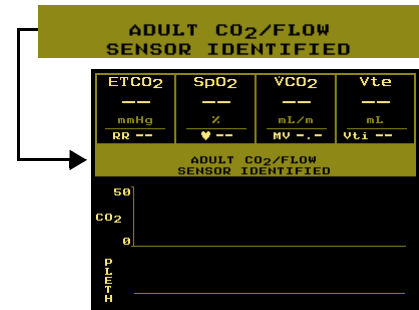

To perform an Adapter Zero:

- 1 Press the **MENU** key.
- 2 Highlight and then select **SETUP** by rotating and then pressing the knob.
- 3 Highlight and then select **CO<sub>2</sub> ZERO NOW**.
- 4 Place the CAPNOSTAT® onto a clean and dry CO<sub>2</sub>/Flow sensor (adapter) that is exposed to room air and away from all sources of CO<sub>2</sub>. Alternatively, the CO<sub>2</sub> adapter provided with the CAPNOSTAT® can be used for the zero procedure.
- 5 Highlight and then select **START**.
  - **ADAPTER ZERO IN PROGRESS PLEASE WAIT** is displayed during the 10 seconds needed to complete the process. Upon completion, **ADAPTER ZERO SUCCESSFUL** is briefly displayed before the monitor automatically returns to the **SETUP** screen.
- 6 Highlight and then select **EXIT** to return to the previous monitoring screen.



**CAUTION**

Position the airway adapter with its windows in a vertical and NOT a horizontal position: this helps keep patient secretions from “pooling” on the windows.

To prevent “rain-out” and moisture from draining into the airway adapter, do NOT place the airway adapter in a gravity dependent position.

Periodically check the CO<sub>2</sub>/Flow sensor and tubing for excessive moisture or secretion build up.



## CAPNOSTAT® CO<sub>2</sub> Sensor


This section provides information regarding the CAPNOSTAT® CO<sub>2</sub> Sensor and its use with NICO Sensors™ and the NICO® monitor.

### The CAPNOSTAT®



The CAPNOSTAT® CO<sub>2</sub> Sensor is a rugged, solid-state, mainstream sensor. It is factory calibrated and does *not* require further calibration during use.

### Connecting the CAPNOSTAT®


The CAPNOSTAT® may be connected to and removed from the NICO® monitor while the monitor is turned off or on.

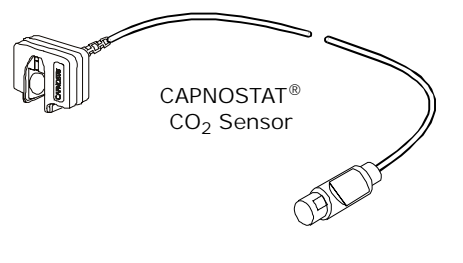


CAPNOSTAT®  
CO<sub>2</sub> Sensor

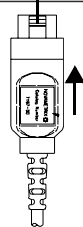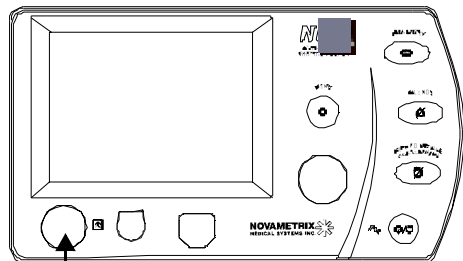

Insert:  
Push in.  
Do not twist.  
Remove:  
Pull out.  
Do not twist.

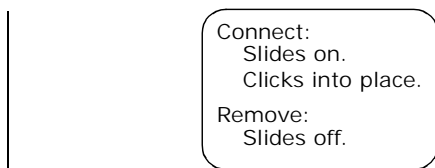



Connect:  
Slides on.  
Clicks into place.  
Remove:  
Slides off.

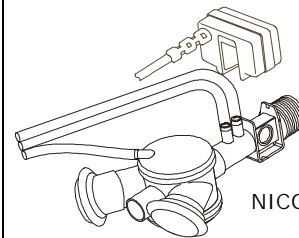

NICO Sensor™

Pediatric/adult  
CO<sub>2</sub>/Flow sensor

Pediatric or neonatal  
CO<sub>2</sub>/Flow sensor

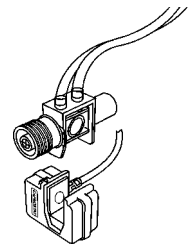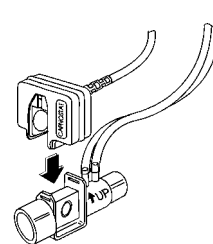

#### CAUTION

Connect only a Novametrix CAPNOSTAT® CO<sub>2</sub> Sensor, Catalog Number 9567-00, to the NICO® monitor. Do not use other CAPNOSTAT® types or other CO<sub>2</sub> sensors.



#### WARNING


Before use, verify the sensor is physically intact, with no broken/frayed wires or damaged parts. Do not use a broken or damaged sensor or one with wet, contaminated, or corroded connectors.

## CAPNOSTAT® Adapter Zero

An “Adapter Zero”, a quick 15 second process that lets the NICO® monitor “learn” the special characteristics of a particular CAPNOSTAT®, is necessary only when the NICO® monitor requests an Adapter Zero.

Such a request may occur the first time a particular CAPNOSTAT® is connected to a particular NICO® monitor—as is the case the first time you power up your NICO® monitor and CAPNOSTAT®—or if the monitor detects some change in the CAPNOSTAT®.

Once an Adapter Zero is performed, NICO® can be turned on and off and the CAPNOSTAT® can be unplugged and reconnected without having to calibrate or Adapter Zero sensor. However, if a second CAPNOSTAT® is connected in place of the first, the Adapter Zero process must be performed on the new sensor—and if at a later time, the first CAPNOSTAT® is reconnected, it too will then have to be put through the Adapter Zero process.

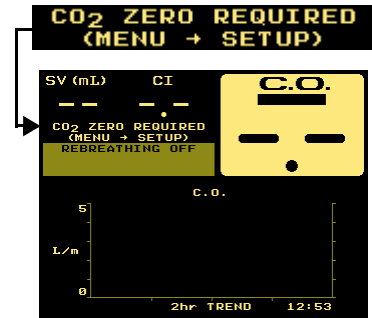

### Adapter Zero

To perform an Adapter Zero:

- 1 Press the **MENU** key.
- 2 Highlight and then select **SETUP** by rotating and then pressing the knob.
- 3 Highlight and then select **CO2 ZERO NOW**.
- 4 Place the CAPNOSTAT® onto a clean and dry NICO Sensor™ or CO<sub>2</sub>/Flow sensor (adapter) that is exposed to room air and away from all sources of CO<sub>2</sub>. Alternatively, a Single Patient Use Adult Airway Adapter (Cat. No. 6063-01) provided with the CAPNOSTAT® can be used for the zero procedure.
- 5 Highlight and then select **START**.
  - **ADAPTER ZERO IN PROGRESS PLEASE WAIT** is displayed during the 10 seconds needed to complete the process. Upon completion, **ADAPTER ZERO SUCCESSFUL** is briefly displayed before the monitor automatically returns to the **SETUP** screen.
- 6 Highlight and then select **EXIT** to return to the previous monitoring screen.






## Pulse Oximetry Sensors



This section details the various pulse oximetry sensors and accessories that can be used with the NICO<sup>®</sup> monitor. This section also explains how to connect a sensor to the monitor, and how to apply the sensor to the patient. NICO<sup>®</sup> uses the pulse oximeter portion of the monitor to enhance shunt corrections as well as to monitor patient oxygenation levels.

### Oximetry Sensors

NICO is compatible with a variety of Novamatrix pulse oximetry sensors including reusable Finger and Y-Sensors<sup>™</sup> and Single Patient Use sensors.

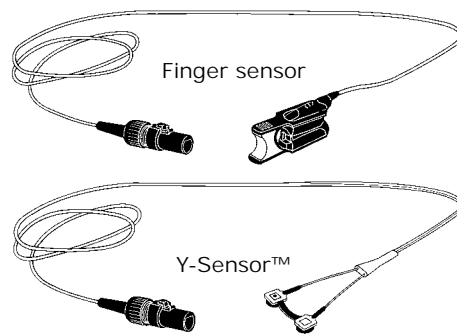


Sensors may be connected and removed with the monitor off or on.

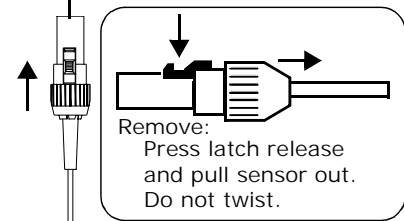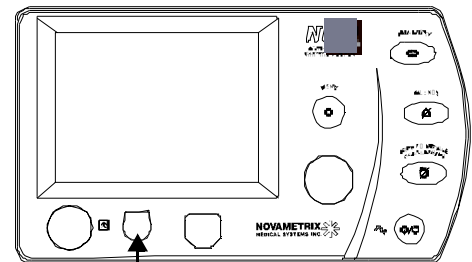


Insert:  
Push in. Do not twist.

**CAUTION**

 Connect only Novamatrix SuperBright<sup>™</sup> SpO<sub>2</sub> sensors, extension cables and accessories with the NICO<sup>®</sup> monitor. Do not use other SpO<sub>2</sub> sensors or accessories.

Overstretching the pulse oximeter finger sensor can damage the sensor and potentially affect pulse oximeter readings. Do not stretch the finger sensor open beyond the limit for which it was designed. Overstretching can be prevented: avoid opening the sensor by any means other than squeezing the grips; *DO NOT* force the sensor onto large objects such as a bedrail.

**WARNING**

 Before applying to the patient, verify the sensor is physically intact, with no broken/frayed wires or damaged parts. Do not use a broken or damaged sensor or one with wet, contaminated, or corroded connectors.

After applying to the patient, inspect the site often for adequate circulation—at least once every four hours. Do not wrap so tightly that circulation is restricted. Note the patient's physiological condition. For example, burn patients may be more sensitive to heat and pressure and require more frequent site checks.

### Sensor Quick Check

A quick functional check of basic sensor operation.

- 1 With the sensor connected to the monitor but not applied to the patient, position the sensor heads so they face each other (red light shines on the detector). Is **PULSE SEARCH** displayed?
- 2 Apply the sensor to your finger. Are reasonable SpO<sub>2</sub> and pulse rate values displayed?
- 3 A YES to BOTH #1 and #2 indicates the sensor is operational. Apply the sensor to the patient as instructed.
  - Note that the Quick Check tag on the sensor may refer to the message **Probe Off Patient**. This message does not apply to the NICO® monitor—**PULSE SEARCH** does.


## Finger Sensor

The reusable Finger Sensor is intended for adult and appropriately sized pediatric fingers and is not intended for neonatal applications.

To apply: Squeeze the grips. Position the fingertip as shown and release the grips.

To remove: Squeeze the grips. Slide the sensor from the finger and release the grips.

Caution: Overstretching can damage the sensor and affect oximetry readings. Do not force the sensor onto large objects such as bedrails.

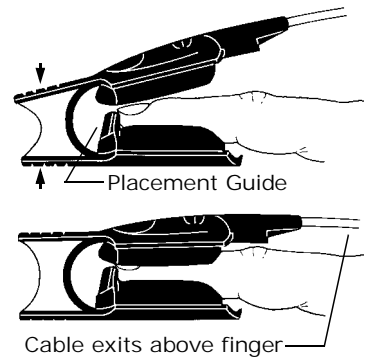


## Y-Sensor™

The reusable Y-Sensor™ is designed for use on all patients from adults to neonates.

### Y-Sensor™ configuration


The Y-Sensor's™ center strip may be carefully cut away if the distance between the sensor heads needs to be reduced to less than 25mm.

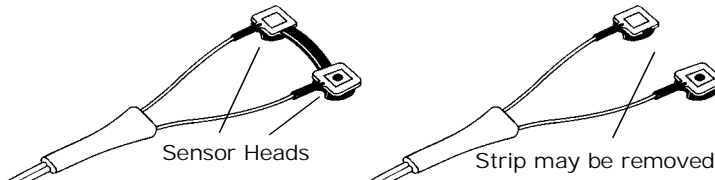


### Y-Sensor™ applicators

The flexible and versatile Y-Sensor™ is applied to the patient using a variety of adhesive and non-adhesive applicators.


Treat applicators (except ear-clip) in accordance with hospital protocol for single-patient use. Refer to instructions packaged with the various applicators.

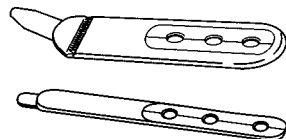


- 6929-00. Adhesive Foam Wraps, Large — Adult, pediatric or neonatal use.
- 6968-00. Adhesive Foam Wraps, Small — Neonatal or appropriately sized pediatric patient.

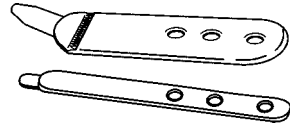

- 8836-00. Non-Adhesive Foam Wrap, Large  
— Adult, pediatric or neonatal use.
- 8943-00. Non-Adhesive Foam Wrap, Small  
— Neonatal or appropriately sized pediatric patient use.

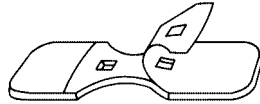


- 8831-00. 20mm Finger Style Y-Strip™ Tape (blue)  
— Appropriately sized pediatric or adult fingers.
- 8832-00. 25mm Finger Style Y-Strip™ Tape (green)  
— Adult fingers.

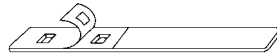

- 8828-00: Wrap Style Tapes 20mm (blue).  
— Neonatal foot, hand, pediatric toe, finger
- Catalog No. 8829: 25mm (green)  
— Neonatal foot, hand

Size refers to the distance between the holes in the tape.


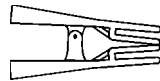


- 6131-00. Ear Clip  
— Adult or pediatric use.

## Using Foam Wraps






To use the foam wraps:

- 1 Press each sensor “button” through the blue side of the foam wrap.
  - Place the head with the red LED closest to the edge, and the other in either remaining hole (removing the Y-Sensor’s™ center strip as required).
- 2 If using an adhesive foam tape, remove both sides of the paper liner.
  - You can also remove the liners prior to buttoning the sensor into the foam wrap if desired.

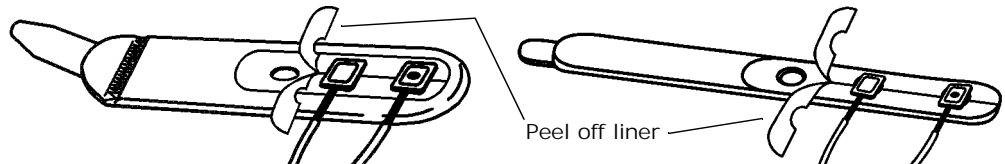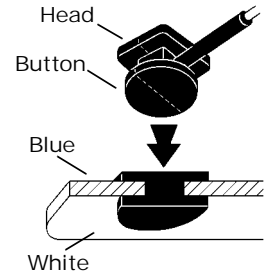


- 3 With the blue foam towards the patient, wrap around the site. Ensure the sensor heads are opposite each other through the tissue. Secure in place with the white plastic tab.

- The tab on the Small foam wrap is removable, allowing shortening for a better fit. Reapply the tab to secure the wrap in place.

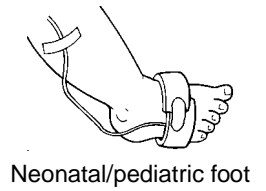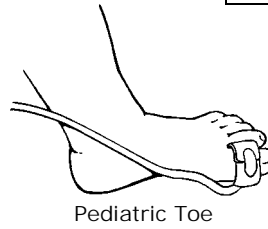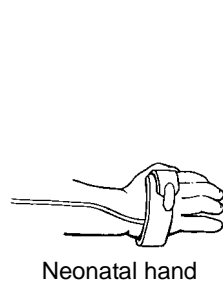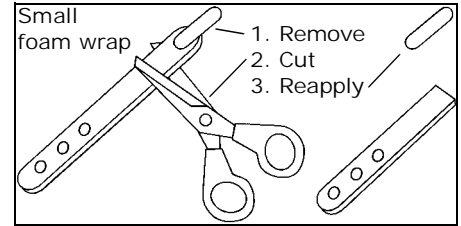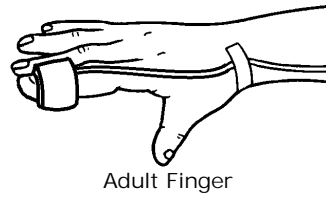

**Using Y-Strip™ Tapes** To use the Y-Strip™ Tapes;




- 1 Remove the release liner with the holes.

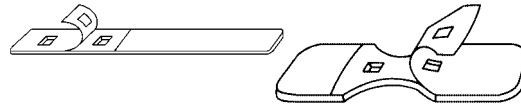

- 2 Press each sensor "button" through the adhesive side of the tape. (Remove the Y-Sensor's™ center strip if required.)

- 3 Remove the remaining release liner. Apply the sensor/ tape to the patient. Ensure the sensor heads are opposite each other.


Adult/Pediatric Finger

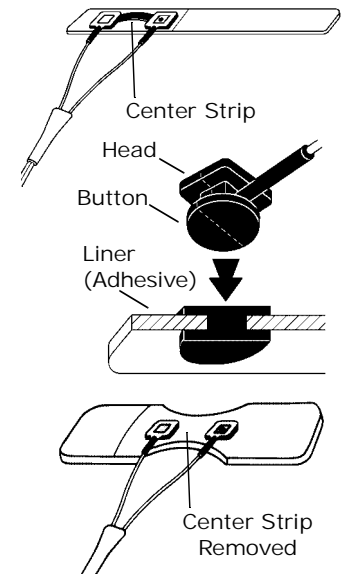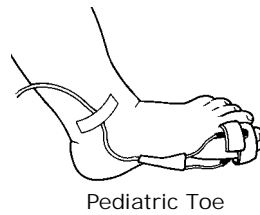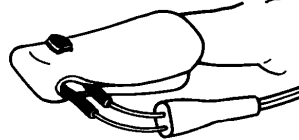

## Using the Ear Clip


To use the ear clip;

- 1 Remove the center strip from the Y-Sensor™.
- 2 Slide each Y-Sensor™ head into an Ear Clip receptacle with the blue button facing outwards.
- 3 Open the clip by squeezing its ends and apply it to the ear.
  - It may be necessary to rub the ear with your fingers in order to increase circulation prior to applying the sensor.
  - Adhesive Dots (8700-00) are included with the ear clip to help hold the clip to the ear.

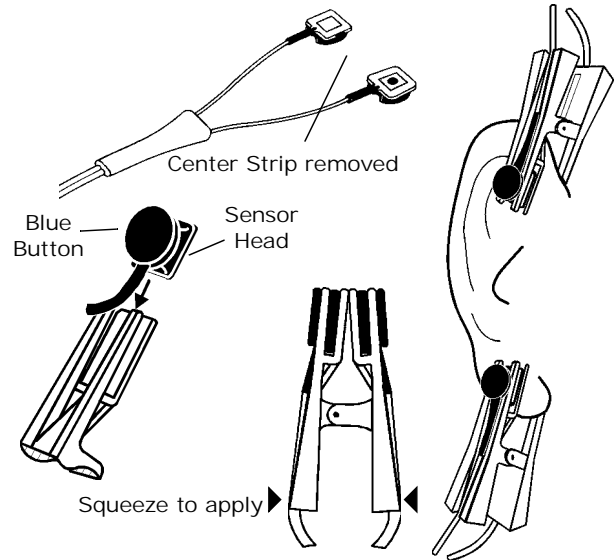

## Single Patient Use Sensors

These Single Patient Use Sensors are for use on appropriately sized patients.

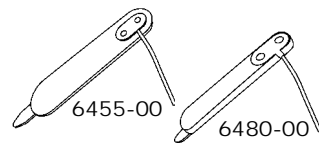

- 6455-00. Pediatric/Adult, Foam Wrap Style  
— Adult or appropriately sized pediatric patients.
- 6480-00. Neonatal/Pediatric, Foam Wrap Style  
— Neonatal or appropriately sized pediatric patients.


### CAUTION:





Single Patient Use SpO<sub>2</sub> sensors can be reapplied to a single patient as needed, but should not be used across multiple patients. Single Patient Use sensors should not be cleaned or disinfected. System performance may be compromised as a result. Replace sensor instead.

- 1 Select the appropriate size sensor based on the patient type.
- 2 With the blue foam towards the patient, wrap around the site. Ensure the sensor heads are opposite each other through the tissue. Secure in place with the white plastic tab.
  - Double-sided adhesive dots, included with the sensor, can be applied over the LED and detector components to help hold the sensor to the site.


- The tab on the Neonatal/Pediatric sensor is removable, allowing shortening for a better fit. Reapply the tab to secure the sensor in place.

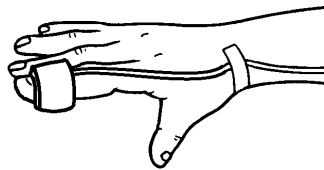


Adult/Pediatric Finger

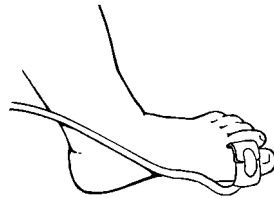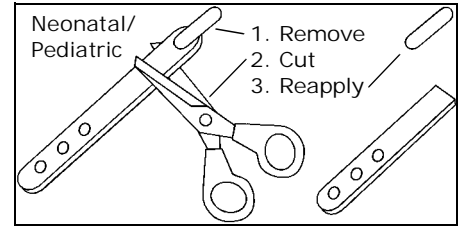


Pediatric Toe

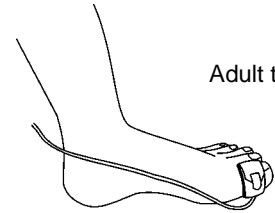

Adult toe

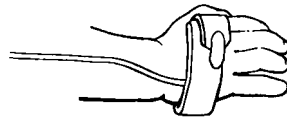

Neonatal hand

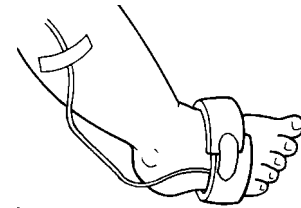

Neonatal foot




# Messages

## Message Areas - Cardiac Output mode

For Respiratory Mechanics messages, see page 78.

In Cardiac Output mode, the monitor uses two screen locations, the General Message Area and the C.O. Message Area, to convey information to the user.

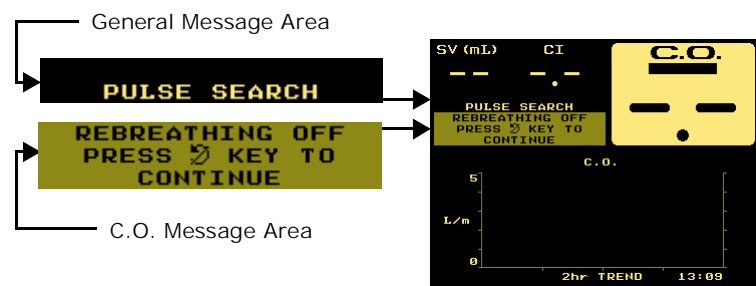

### General Message Area

The General Message Area in Cardiac Output mode displays system status, alert, and error conditions. It may be blank or:

- contain 1 single-line message
- contain 2 single-line messages
- contain 1 multi-line message



The General Message Area messages are listed below, in alphabetical order.

(Alert Class: H-High Priority, M-Medium Priority, L-Low Priority, S-Status Message. See "NICO® Alert Priorities" on page 56 for details.)

General Message Area	Message Description	Alert Class
ALERTS OFF	Displayed as a reminder that the default for all user selectable alerts is OFF. To cancel the message, adjust any individual limit value and activate the audible alert. Both the alert limit and audible alert must be activated.	S
AMBIENT LIGHT COVER SpO <sub>2</sub> PROBE	The monitor detects interference on the SpO <sub>2</sub> sensor from ambient light. This can be corrected by covering the SpO <sub>2</sub> sensor, or possibly by changing the sensor site.	S
CHECK / CHANGE AIRWAY ADAPTER	A change in the CO <sub>2</sub> adapter portion of the NICO Sensor™ is detected. Possible causes: <ul style="list-style-type: none"> <li>• CAPNOSTAT® CO<sub>2</sub> sensor off adapter</li> <li>• High level of moisture and/or secretions in the adapter. The moisture can be drained from the NICO Sensor™ and re-inserted in the circuit. Facing the pneumatic tubing upward as it exits the NICO Sensor™ can minimize this.</li> </ul>	S

General Message Area	Message Description	Alert Class
CO <sub>2</sub> SENSOR FAILURE REPLACE SENSOR	A problem with the CAPNOSTAT <sup>®</sup> CO <sub>2</sub> sensor has been identified. Replace the CAPNOSTAT <sup>®</sup> and return it to Novamatrix for exchange or repair.	M
CO <sub>2</sub> SENSOR?	The CAPNOSTAT <sup>®</sup> CO <sub>2</sub> sensor is not connected to the NICO <sup>®</sup> monitor.	S
CO <sub>2</sub> ZERO REQUIRED (MENU → SETUP)	The CAPNOSTAT <sup>®</sup> CO <sub>2</sub> sensor needs to be zeroed. Press the <b>MENU</b> key, then select <b>SETUP</b> , then <b>CO2 ZERO NOW</b> , and follow the instructions on the screen. NOTE: the CO <sub>2</sub> adapter provided with the CAPNOSTAT <sup>®</sup> can be used for the CO <sub>2</sub> zero procedure rather than a new NICO Sensor <sup>™</sup> .	S
DEMO MODE	The monitor is in demonstration mode and is not displaying patient data (all data is simulated). To exit demo mode, turn the monitor off, then back on.	S
ETCO <sub>2</sub> : XX mmHg	End tidal CO <sub>2</sub> is greater than 60 mmHg or exceeds the high alert limit. Appears in the general message area to supplement screens that do not display the ETCO <sub>2</sub> value.	L
HIGH C.O.	The displayed cardiac output value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
HIGH ETCO <sub>2</sub>	The displayed ETCO <sub>2</sub> value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
HIGH PULSE	The displayed pulse rate value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
HIGH RESP RATE	The displayed respiration rate value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
HIGH SpO <sub>2</sub>	The displayed SpO <sub>2</sub> value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
INCOMPATIBLE CO <sub>2</sub> SENSOR	The wrong part number CAPNOSTAT <sup>®</sup> is connected to the NICO <sup>®</sup> monitor. Use only a CAPNOSTAT <sup>®</sup> with part number 9567-00 (this can be distinguished from other CAPNOSTAT <sup>®</sup> part numbers by the yellow part number label on the CAPNOSTAT <sup>®</sup> 's electrical connector).	M
INCOMPATIBLE FLOW SENSOR	The wrong flow sensor is connected to the NICO <sup>®</sup> monitor. Use only a NICO Sensor <sup>™</sup> , part numbers 8950-00, 8951-00 or 8952-00 (the correct flow sensor is an integral part of the NICO Sensor <sup>™</sup> ).	M
INSP CO <sub>2</sub> : xx	(where xx is a numeric value with units of mmHg, kPa, or %). At least 3 mmHg, 0.1% or 0.1 kPa of CO <sub>2</sub> has been detected during inspiration (other than during rebreathing) for at least ten continuous seconds.	S
LOW C.O.	The displayed cardiac output value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
LOW ETCO <sub>2</sub>	The displayed ETCO <sub>2</sub> value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
LOW PULSE	The displayed pulse rate value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H




General Message Area	Message Description	Alert Class
LOW RESP RATE	The displayed respiration rate value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
LOW SpO <sub>2</sub>	The displayed SpO <sub>2</sub> value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
MONITOR INOPERABLE AIRWAY ZERO ERROR	The NICO® monitor detected a problem with its pneumatic flow and pressure sub-system. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE BARO PRESS ERROR	The NICO® monitor detected a problem with its internal barometric pressure sensor. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE CLOCK FAILURE	The NICO® monitor detected a problem with its internal clock. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE FLOW ZERO ERROR	The NICO® monitor detected a problem with its flow zeroing sub-system or related pneumatic components. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE NICO VALVE ZERO ERR	The NICO® monitor detected a problem with its internal rebreathing valve control circuitry or related pneumatic components. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE SpO <sub>2</sub> HDW ERROR	The NICO® monitor detected a problem with its pulse oximetry sub-system. Contact qualified personnel for monitor repair or exchange.	L
NO RESP: xx:xx	The time selected in the <b>SET ALERTS</b> screen for the NO RESP (no respiration) alert was exceeded since the end of the last detected breath (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
PULSE SEARCH	The pulse oximeter is not detecting a sufficient pulse. This could be due to: <ul style="list-style-type: none"> <li>• SpO<sub>2</sub> sensor is off of the patient</li> <li>• Insufficient perfusion at the site</li> <li>• Tissue at the site is too thick or too thin</li> </ul>	S
SpO <sub>2</sub> PROBE FAILURE	The pulse oximeter sensor is faulty. Replace the sensor and contact qualified service personnel.	M
SpO <sub>2</sub> PROBE?	<ul style="list-style-type: none"> <li>• The pulse oximeter sensor was not connected to the NICO® monitor when it was powered up.</li> <li>• (Flashing) SpO<sub>2</sub> sensor was disconnected from the monitor after it was powered up. Acknowledge by pressing the Silence key.</li> </ul>	S
WARMUP	The CAPNOSTAT® CO <sub>2</sub> sensor is not at proper operating temperature yet.	S



### C.O. Message Area

The C.O. Message Area is dedicated to cardiac output related information and may be blank or contain a message.

The C.O. Message Area messages are listed below, in alphabetical order.

(Alert Class: H-High Priority, M-Medium Priority, L-Low Priority, S-Status Message. See “NICO® Alert Priorities” on page 56 for details.)

C.O. Message Area	Message Description	Alert Class
BLOOD GASES HAVE BEEN RESET	A NICO® cycle has started after 20 minutes of a no respiration timeout. Data not available when no breaths are detected 20 minutes after the last detected breath.	S
CO <sub>2</sub> UNSTABLE CO <sub>2</sub> STABILIZING	The NICO® monitor was not able to calculate a cardiac output value. This can be due to: <ul style="list-style-type: none"> <li>• Spontaneous breaths or efforts</li> <li>• Surgeon moving the lungs</li> <li>• Ventilator adjustments</li> </ul>	L
CONSIDER USING SMALL NICO SENSOR	Resizing the NICO Loop™ was not effective because the ventilator set tidal volume is less than 300 ml. The NICO® monitor is suggesting a different sized sensor, with smaller loop, to correct the condition.	L
CONSIDER USING STANDARD NICO SENSOR	If a small NICO Sensor™ is currently in use: Resizing the NICO Loop™ was not effective because the ventilator set tidal volume is greater than 500 ml. The NICO® monitor is suggesting a different sized sensor, with a larger loop, to correct the condition.  If a large NICO Sensor™ is currently in use: Resizing the NICO Loop™ was not effective because the ventilator set tidal volume is less than 1000 ml. The NICO® monitor is suggesting a different sized sensor, with a smaller loop, to correct the condition.	L
CONSIDER USING LARGE NICO SENSOR	Resizing the NICO Loop™ was not effective because the ventilator set tidal volume is greater than 1000 ml. The NICO® monitor is suggesting a different sized sensor, with a larger loop, to correct the condition.	L
EXPAND LOOP	The NICO Loop™ (expandable rebreathing volume on the NICO Sensor™) needs to be expanded. Expand the loop by approximately 3-6 inches or until the message is removed. Note that the message is displayed only during the rebreathing phase of the NICO® cycle. If the loop is not appropriately sized by the end of the rebreathing phase, the message will be removed and displayed again during the next rebreathing phase. If this message persists with maximal expansion of the NICO Loop™, the tidal volume may be too large for the ventilatory conditions for NICO® to report accurate results.	L
NEXT  : xx:xx	There are x:xx minutes:seconds until the beginning of the next rebreathing period (provided as an indicator as to the current state of the NICO® cycle).	S
NICO SENSOR?	The NICO Sensor™ has not been connected to the monitor since it was last turned on.	S

C.O. Message Area	Message Description	Alert Class
 REBREATHING	The patient is currently rebreathing a portion of his/her tidal volume in order for NICO <sup>®</sup> to calculate cardiac output (provided as an indicator of the current state of the NICO <sup>®</sup> cycle). The rebreathing phase of the NICO <sup>®</sup> cycle lasts for 50 seconds.	S
REBREATHING OFF	Rebreathing and therefore C.O. measurements are currently disabled. The <b>STOP/CONTINUE REBREATHING</b> key is illuminated amber while rebreathing is off, and can be pressed to enable rebreathing and C.O. measurements. Rebreathing is off when: <ul style="list-style-type: none"> <li>The monitor is first turned on until the <b>STOP/CONTINUE REBREATHING</b> key is pressed</li> <li>The <b>STOP/CONTINUE REBREATHING</b> key is pressed while C.O. measurements are enabled</li> <li>The monitor detected a system fault or condition which warrants automatic disabling of C.O. measurements</li> </ul>	S
REBREATHING OFF ADAPTER DISCONNECT REMOVE FROM CIRCUIT	The NICO Sensor™ was disconnected from the monitor after C.O. measurements had been made since the last power-up. The NICO Sensor™ should be removed from the breathing circuit in order to avoid leaking breathing circuit gas through the NICO Sensor™'s connector.	S
REBREATHING OFF LARGE NICO SENSOR IDENTIFIED	A large size NICO Sensor™ was just connected to the NICO <sup>®</sup> monitor.	S
REBREATHING OFF SMALL NICO SENSOR IDENTIFIED	A small size NICO Sensor™ was just connected to the NICO <sup>®</sup> monitor.	S
REBREATHING OFF STANDARD NICO SENSOR IDENTIFIED	A standard size NICO Sensor™ was just connected to the NICO <sup>®</sup> monitor.	S
REBREATHING OFF NICO SENSOR FAILURE	A problem with the NICO Sensor™ was detected by the monitor. Discard the sensor and replace it. If the problem persists, contact qualified service personnel.	M
REBREATHING OFF PRESS  KEY TO CONTINUE	Rebreathing and therefore C.O. measurements are currently disabled and can be enabled by pressing the <b>STOP / CONTINUE REBREATHING</b> key.	S
REBREATHING PAUSED WAITING FOR ETCO <sub>2</sub> < XX	Rebreathing and therefore C.O. measurements have been temporarily paused automatically by the monitor, and will resume automatically once the indicated parameter is within the stated range (here, XX = 85 mmHg, 11.5 kPa, 11.5 %).	L
REBREATHING PAUSED WAITING FOR ETCO <sub>2</sub> > XX	Rebreathing and therefore C.O. measurements have been temporarily paused automatically by the monitor, and will resume automatically once the indicated parameter is within the stated range (here, XX = 15 mmHg, 2.0 kPa, 2.0 %).	L
REBREATHING PAUSED WAITING FOR RR < 60 br/m	Rebreathing and therefore C.O. measurements have been temporarily paused automatically by the monitor, and will resume automatically once the indicated parameter is within the stated range.	L

C.O. Message Area	Message Description	Alert Class
REBREATHING PAUSED WAITING FOR RR > 3 br/m	Rebreathing and therefore C.O. measurements have been temporarily paused automatically by the monitor, and will resume automatically once the indicated parameter is within the stated range.	L
REBREATHING PAUSED WAITING FOR VCO <sub>2</sub> > 20 mL/min	Rebreathing and therefore C.O. measurements have been temporarily paused automatically by the monitor, and will resume automatically once the indicated parameter is within the stated range.	L
REBREATHING PAUSED WAITING FOR Vt > 200 mL	Rebreathing and therefore C.O. measurements have been temporarily paused automatically by the monitor, and will resume automatically once the indicated parameter is within the stated range.	L
RETRACT LOOP	The NICO Loop™ (expandable rebreathing volume on the NICO Sensor™) needs to be retracted (made smaller). Retract the loop by approximately 3-6 inches or until the message is removed. Note that the message is displayed only during the rebreathing phase of the NICO® cycle. If the loop is not appropriately sized by the end of the rebreathing phase, the message will be removed and displayed again during the next rebreathing phase. If this message persists with the NICO Loop™ at its minimal size, the tidal volume may be too small for the ventilatory conditions for NICO® to report accurate results.	L

## Message Areas - Respiratory Mechanics mode

For Cardiac Output messages, see page 73.

In Respiratory Mechanics mode, the monitor uses one General Message Area to convey information to the user.

General Message Area


WARMUP  
PULSE SEARCH  
FLOW SENSOR?



### General Message Area

The General Message Area in Respiratory Mechanics mode displays all messages, including system status, alert, and error conditions. It may be blank or:

- contain 1, 2 or 3 single-line messages
- contain 1 multi-line message
- contain 1 multi-line message and 1 single-line message

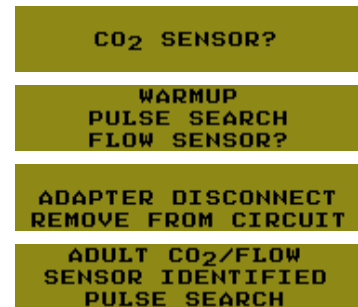

The General Message Area messages are listed below, in alphabetical order.

(Alert Class: H-High Priority, M-Medium Priority, L-Low Priority, S-Status Message. See "NICO® Alert Priorities" on page 56 for details.

General Message Area	Message Description	Alert Class
ADAPTER DISCONNECT REMOVE FROM CIRCUIT	The CO <sub>2</sub> /Flow sensor is not connected to the NICO® monitor. The CO <sub>2</sub> /Flow sensor should be removed from the breathing circuit in order to avoid leaking breathing circuit gas through the sensor's connector.  NOTE: If a two-line message will not fit in the General Message Area because other higher priority messages are present, the message "FLOW SENSOR?" will be used for the same condition.	S
ADULT CO <sub>2</sub> /FLOW SENSOR IDENTIFIED	An adult size CO <sub>2</sub> /Flow sensor was just connected to the NICO® monitor.	S
ALERTS OFF	Displayed as a reminder that the default for all user selectable alerts is OFF. To cancel the message, adjust any individual limit value and activate the audible alert. Both the alert limit and audible alert must be activated.	S
AMBIENT LIGHT COVER SpO <sub>2</sub> PROBE	The monitor detects interference on the SpO <sub>2</sub> sensor from ambient light. This can be corrected by covering the SpO <sub>2</sub> sensor, or possibly by changing the sensor site.	S
CHECK / CHANGE AIRWAY ADAPTER	A change in the CO <sub>2</sub> adapter portion of the CO <sub>2</sub> /Flow sensor is detected. Possible causes: <ul style="list-style-type: none"> <li>• CAPNOSTAT® CO<sub>2</sub> sensor off adapter</li> <li>• High level of moisture and/or secretions in the adapter. Replace if needed.</li> </ul>	S
CO <sub>2</sub> SENSOR FAILURE REPLACE SENSOR	A problem with the CAPNOSTAT® CO <sub>2</sub> sensor has been identified. Replace the CAPNOSTAT® and return it to Novamatrix for exchange or repair.	M
CO <sub>2</sub> SENSOR?	The CAPNOSTAT® CO <sub>2</sub> sensor is not connected to the NICO® monitor.	S
CO <sub>2</sub> ZERO REQUIRED (MENU → SETUP)	The CAPNOSTAT® CO <sub>2</sub> sensor needs to be zeroed. Press the <b>MENU</b> key, then select <b>SETUP</b> , then <b>CO<sub>2</sub> ZERO NOW</b> , and follow the instructions on the screen. NOTE: the CO <sub>2</sub> adapter provided with the CAPNOSTAT® can be used for the CO <sub>2</sub> zero procedure rather than a new sensor.	S
DEMO MODE	The monitor is in demonstration mode and is not displaying patient data (all data is simulated). To exit demo mode, turn the monitor off, then back on.	S
ETCO <sub>2</sub> > XX mmHg	End tidal CO <sub>2</sub> is greater than 60 mmHg or exceeds the high alert limit. Appears in the general message area to supplement screens that do not display the ETCO <sub>2</sub> value.	L
FLOW SENSOR?	The CO <sub>2</sub> /Flow sensor is disconnected.  NOTE: This message may also be used in place of ADAPTER DISCONNECT REMOVE FROM CIRCUIT, when a two-line message will not fit in the General Message Area.	S
HIGH ETCO <sub>2</sub>	The displayed ETCO <sub>2</sub> value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
HIGH PIP	The displayed PIP value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M

General Message Area	Message Description	Alert Class
HIGH PULSE	The displayed pulse rate value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
HIGH RESP RATE	The displayed respiration rate value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
HIGH SpO <sub>2</sub>	The displayed SpO <sub>2</sub> value is above the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
INCOMPATIBLE CO <sub>2</sub> SENSOR	The wrong part number CAPNOSTAT <sup>®</sup> is connected to the NICO <sup>®</sup> monitor. Use only a CAPNOSTAT <sup>®</sup> with part number 9567-00 (this can be distinguished from other CAPNOSTAT <sup>®</sup> part numbers by the yellow part number label on the CAPNOSTAT <sup>®</sup> 's electrical connector).	M
INCOMPATIBLE FLOW SENSOR	The wrong flow sensor is connected to the NICO <sup>®</sup> monitor. Use only a NICO <sup>®</sup> CO <sub>2</sub> /Flow sensor part number 9765-00, 9766-00, or 9767-00. If message persists, a hardware error may exist. Contact qualified personnel for monitor repair or exchange.	M
INSP CO <sub>2</sub> : xx	(where xx is a numeric value with units of mmHg, kPa, or %). At least 3 mmHg, 0.1% or 0.1 kPa of CO <sub>2</sub> has been detected during inspiration (other than during rebreathing) for at least ten continuous seconds.	S
LOW ETCO <sub>2</sub>	The displayed ETCO <sub>2</sub> value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
LOW PULSE	The displayed pulse rate value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
LOW RESP RATE	The displayed respiration rate value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	M
LOW SpO <sub>2</sub>	The displayed SpO <sub>2</sub> value is below the set alert limit in the <b>SET ALERTS</b> screen (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
MONITOR INOPERABLE AIRWAY ZERO ERROR	The NICO <sup>®</sup> monitor detected a problem with its pneumatic flow and pressure sub-system. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE BARO PRESS ERROR	The NICO <sup>®</sup> monitor detected a problem with its internal barometric pressure sensor. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE CLOCK FAILURE	The NICO <sup>®</sup> monitor detected a problem with its internal clock. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE FLOW ZERO ERROR	The NICO <sup>®</sup> monitor detected a problem with its flow zeroing sub-system or related pneumatic components. Contact qualified personnel for monitor repair or exchange.	L
MONITOR INOPERABLE SpO <sub>2</sub> HDW ERROR	The NICO <sup>®</sup> monitor detected a problem with its pulse oximetry sub-system. Contact qualified personnel for monitor repair or exchange.	L
NEONATAL CO <sub>2</sub> /FLOW SENSOR IDENTIFIED	A neonatal size CO <sub>2</sub> /Flow sensor was just connected to the NICO <sup>®</sup> monitor.	S

General Message Area	Message Description	Alert Class
NO RESP: xx:xx	The time selected in the <b>SET ALERTS</b> screen for the NO RESP (no respiration) alert was exceeded since the end of the last detected breath (press <b>MENU</b> key, then select <b>SET ALERTS</b> to view the alert limit settings).	H
PEDIATRIC CO <sub>2</sub> /FLOW SENSOR IDENTIFIED	A pediatric size CO <sub>2</sub> /Flow sensor was just connected to the NICO <sup>®</sup> monitor.	S
PULSE SEARCH	The pulse oximeter is not detecting a sufficient pulse. This could be due to: <ul style="list-style-type: none"> <li>• SpO<sub>2</sub> sensor is off of the patient</li> <li>• Insufficient perfusion at the site</li> <li>• Tissue at the site is too thick or too thin</li> </ul>	S
SpO <sub>2</sub> PROBE FAILURE	The pulse oximeter sensor is faulty. Replace the sensor and contact qualified service personnel.	M
SpO <sub>2</sub> PROBE?	<ul style="list-style-type: none"> <li>• The pulse oximeter sensor was not connected to the NICO<sup>®</sup> monitor when it was powered up.</li> <li>• (Flashing) SpO<sub>2</sub> sensor was disconnected from the monitor after it was powered up. Acknowledge by pressing the Silence key.</li> </ul>	S
WARMUP	The CAPNOSTAT <sup>®</sup> CO <sub>2</sub> sensor is not at proper operating temperature yet.	S



## External Devices

---

---

### GE Medical Systems Solar® Interface

The Solar® monitor from GE Medical Systems Information Technologies (GEMS-IT), connects to the NICO® via an Octanet and a DIDCA™ interface adapter.

Solar® interface software is compatible with NICO® release 2.1 or above. If your NICO® software is updated, please complete and fax the software notification found in your GEMS-IT service manual, so they can confirm continued compatibility.

#### Preparing for Use

To setup and connect the NICO® monitor to the GEMS-IT Solar® monitor, you will need the following:

- NICO® monitor with software revision 3.1 or greater
- GEMS-IT Solar® monitor, Model 7000, or 8000 with software version V6A or higher, or GEMS-IT Solar® monitor, Model 8000M with software version V1A or higher
- GEMS-IT Octanet module, (PN: OCTANET=A)
- DIDCA™ interface adapter, (PN: 420915-058)
- GEMS-IT Octanet cable (PN: 418335-00x). Cables come in various lengths and are available from GEMS-IT.
- GEMS-IT Solar® monitor cable (PN: 409752-00x or 700520-00x)
- Optional: GEMS-IT TRAM-Net Hub (PN: 410217-001) and TRAM-Net cable (PN: 409752-00x). Instructions for use of the TRAM-Net are detailed within Octanet and other service documents provided by GEMS-IT.

NOTE: All GEMS-IT parts are supplied by GEMS-IT.

#### Component Setup

Refer to the following steps for correct cable connections.

- 1** Connect the DIDCA™ interface adapter to the connector labeled RS232 1 on the back of the NICO®. The interface adapter is programmed specifically to work with a NICO® monitor.
- 2** Connect the Octanet cable between the interface adapter and one of the eight Octanet serial ports.
- 3** Connect the Solar monitor cable (PN: 409752-00x for Models 7000 and 8000, PN: 700752-00x for model 8000M) from the Octanet to the Solar® monitor

Complete the communication cabling by referring to the power up instructions found in the Octanet Connectivity Device Service Manual (PN: 418264-003).

- 4** Ensure Octanet serial port LED is steady green.

#### Transmitted Parameters

The following parameters are transmitted from the NICO® monitor to the Solar® monitor.

- The CO, CI, SV, CO confidence level and PCBF values will be displayed in the NICO parameter block.
- Four of the following sub parameters can be displayed in the RM parameter block: PEF, MV, MVs, MVm, TV, TVs, TVm, PIP, MAWP, PEEP, RR, RRs, RRm, I:E, CYDN, RAWe.



- The CO<sub>2</sub> inspired, CO<sub>2</sub> expired, and respiration rate (RR) values will be displayed in the CO<sub>2</sub> parameter block.
- The SPO<sub>2</sub> value and pulse rate (RATE) value will be displayed in the SPO<sub>2</sub> parameter block.

## Agilent Technologies VueLink Interface

The NICO® interface to the Agilent Technologies<sup>1</sup> VueLink module provides a pathway for NICO® data to be viewed on the Agilent Technologies Patient Monitoring System. The data from NICO® is available for display in conjunction with the other parameters configured for display on the Agilent Technologies Patient Monitor.

### Preparing for Use



To setup and connect the NICO® monitor to the Agilent VueLink module, you will need the following:

- NICO® monitor with software revision 3.1 or greater.
- Agilent Patient Monitor with software revision C (9.xx) or higher.
- Agilent VueLink module, Auxiliary Plus B with Open Interface (product number M1032 option A05).
- Agilent 4-meter VueLink Interface Cable (part number K6B)
- 25-pin female to 9-pin male null modem (crossover) cable

NOTE: The Agilent Patient Monitor, VueLink Module and Agilent interface cable are supplied by Agilent.


### Component Setup

- 1 Connect the CAPNOSTAT® CO<sub>2</sub> sensor, flow sensor, pulse oximeter probe and AC power cord to the NICO® monitor.
- 2 Connect the 25-pin male connector on the VueLink Interface Cable to the 25-pin female connector on the null modem cable, then connect the 9-pin male connector to the rear panel of the NICO®. Tighten screws.
- 3 Connect the other end of the VueLink cable to the VueLink module by grasping the end of the connector and pushing it into the receptacle on the VueLink module.
- 4 Plug the VueLink Module into an available slot in the Agilent monitoring system module rack.
  - The VueLink module may be stored in the rack when not in use.

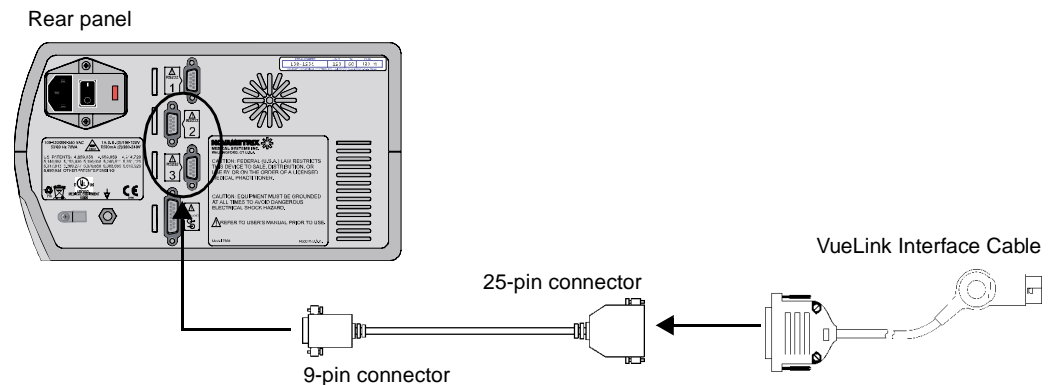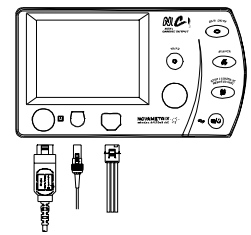

1. Agilent Technologies was formerly Hewlett-Packard. The information above also applies to existing HP products.

- 1 Press the **Operate/Standby** key to turn the monitor on and off.
  - NICO® can operate from its internal battery or from the AC Mains. (See “AC/Battery Operation” on page 4 for details.)
- 2 Press the **MENU** key to activate **SELECT A SCREEN**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 3 Highlight and then select **SETUP** by turning and then pressing the **KNOB**.
- 4 The **SETUP** screen is displayed.
- 5 Turn and press the **KNOB** to select **INPUT/OUTPUT**.
- 6 Turn and press the **KNOB** to select **RS232-2** or **RS232-3**.
- 7 Turn and press the **KNOB** to select **VUELINK**.

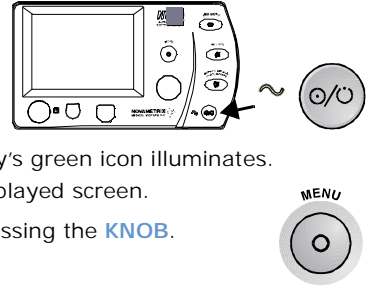

NOTE: If you need detailed installation and operating instructions for the Agilent VueLink Module and the Agilent Patient Monitor, refer to the Agilent Operator’s Manual.

### Transmitted Parameters

A total of twenty-four parameters are transmitted to the Agilent Patient Monitor. The parameter selection varies depending upon whether you are in Cardiac Output mode or Respiratory Mechanics mode. Seventeen of these parameters will display in the Agilent Main Screen, a maximum of six at a time. Seven parameters display in the Task Window only.

Label	Parameter	Units
AWRR	Respiratory Rate	rpm
Cdyn	Dynamic Compliance	ml/cmH <sub>2</sub> O
RAW e	Expired Airway Resistance (Dynamic)	cmH <sub>2</sub> O/L/S
PIP	Peak Inspiratory Pressure	cmH <sub>2</sub> O
PEEP	Positive End Expiratory Pressure	cmH <sub>2</sub> O
Pmean	Mean Airway Pressure	cmH <sub>2</sub> O
Pulse	Pulse Rate	bpm
SpO <sub>2</sub>	SpO <sub>2</sub>	%
ETCO <sub>2</sub>	ETCO <sub>2</sub>	*mmHg, %, kPa
PECO <sub>2</sub>	Mixed Expired CO <sub>2</sub>	*mmHg, %, kPa
PEF	Peak Expiratory Flow	L/min
TVin	Inspired Tidal Volume	ml
TVex	Expired Tidal Volume	ml
MV t	Total Minute Volume	L
MV s	Spontaneous Minute Volume	L
MV m	Mechanical Minute Volume	L
VCO <sub>2</sub>	CO <sub>2</sub> Elimination	ml/min
VD aw	Airway Deadspace	ml
Respiratory Mechanics Mode only:		
RAW i	Inspired Airway Resistance (Dynamic)	cmH <sub>2</sub> O/L/S
VT alv t	Total Alveolar Tidal Volume	ml

VT alv s	Spontaneous Alveolar Tidal Volume	ml
VT alv m	Mechanical Alveolar Tidal Volume	ml
MV alv	Alveolar Minute Volume (Total)	L
VdVtPh	Physiologic deadspace to tidal volume ratio	
Cardiac Output Mode only:		
CO-a	Average Cardiac Output	L/min
CO-f	Fast-mode Cardiac Output	L/min
CI	Cardiac Index	L/min/m <sup>2</sup>
PCBF	Pulmonary Capillary Blood Flow	L/min
SV	Stroke Volume	ml
SVR	Systemic Vascular Resistance	dynes sec/cm <sup>5</sup>

\*CO<sub>2</sub> parameters are reported in either mmHg, kPa, or %, depending upon the CO<sub>2</sub> units currently selected for the NICO<sup>®</sup> monitor.

**Transmitted Waveforms**


A total of five waveform are transmitted to the Agilent Patient Monitor; a maximum of two may viewed at one time.

Label	Waveform	Units
AWF	Airway Flow	L
AWP	Airway Pressure	cmH <sub>2</sub> O
CO <sub>2</sub>	CO <sub>2</sub> Capnogram	mmHg, %, kPa
PLETH	Plethysmogram	N/A
AWV	Volume	ml

**ASCII Output**

The NICO<sup>®</sup> monitor must be set to the appropriate interface.

- 1 Press the **Operate/Standby** key to turn the monitor on and off.
  - NICO<sup>®</sup> can operate from its internal battery or from the AC Mains. (See "AC/Battery Operation" on page 4 for details.)
- 2 Press the **MENU** key to activate **SELECT A SCREEN**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 3 Highlight and then select **SETUP** by turning and then pressing the **KNOB**.
- 4 The **SETUP** screen is displayed.
- 5 Turn and press the **KNOB** to select **INPUT/OUTPUT**.
- 6 Turn and press the **KNOB** to select **RS232-2** or **RS232-3**.
- 7 Turn and press the **KNOB** to select:
  - **ASCII OUTPUT** to a serial printer
  - **ASCII OUTPUT 2** for data collection using a PC
- 8 Connect the device to the appropriate rear panel connector.

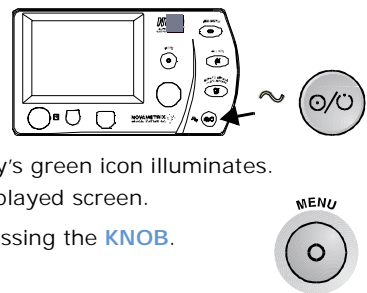

### Output to a Serial Printer

For output to a serial printer such a strip chart recorder or the Seiko DPU-414 printer, configure your serial printer to 9600 bps, no parity, 8 data bits, 1 stop bit.

An output sample is shown below:

JUL 6, 2000		9:47:14		
C.O.	C.I.	S.V.	C.O.-f	PCBF
5.2	2.8	69	5.2	4.9
VCO2	MV	Mvalv	Vti	Vte
0	0.0	0.0	0	0
PIP	MAP	PEEP	Cdyn	Raw
23	9	4	27	15
ETCO2	RR	SpO2	PR	
33.9	15	97	75	

### Output files for use on a PC

When configured to **ASCII OUTPUT 2**, NICO® will output a text file that can be imported into a spreadsheet program such as Microsoft Excel. Configure your serial device to 19200 bps, no parity, 8 data bits, 1 stop bit.

Data is output once every breath. An output sample is shown below:

10:40:04, PR=75, SPO2=97, RR=20, C.O.=0.0, C.O.-f=0.0, CI=0.0, SV=0, VCO2=429, ETCO2=39, PECO2=32, MV=12.2, MVALV=11.0, Vti=916, Vte=907, Vtalv=546, VdAW=63, PIP=45, MAP=22, PEEP=9, Cdyn=27, Raw=13

## Analog Output


NICO® directly supports interface to analog devices.

### Setup

The NICO® monitor must be set to the appropriate interface.

- 1 Press the **Operate/Standby** key to turn the monitor on and off.
  - NICO® can operate from its internal battery or from the AC Mains. (See "AC/Battery Operation" on page 4 for details.)
- 2 Press the **MENU** key to activate **SELECT A SCREEN**. The key's green icon illuminates.
  - Press the key again to return to the previously displayed screen.
- 3 Highlight and then select **SETUP** by turning and then pressing the **KNOB**.
- 4 The **SETUP** screen is displayed.
- 5 Turn and press the **KNOB** to select **INPUT/OUTPUT**.
- 6 Turn and press the **KNOB** to select **ANALOG OUT 1-4**.
- 7 Turn and press the **KNOB** to select the parameters you want from **ANALOG OUT 1-4** (maximum of 1 per), to be transmitted to an analog device like a strip chart recorder, etc.
 

CO, CI, SV, PCBF, ETCO<sub>2</sub>, SpO<sub>2</sub>, resp rate, and pulse rate value outputs are available, as well as CO<sub>2</sub>, plethysmogram, flow, and airway pressure waveform data.
- 8 Calibrate the recorder to the correct voltage levels using the **ANALOG CAL.** setting. Press **ZERO**, **HALF** and **FULL** to set the analog outputs.
  - **ZERO** — 0 volts
  - **HALF** — 0.50 volts
  - **FULL** — 1.00 volts

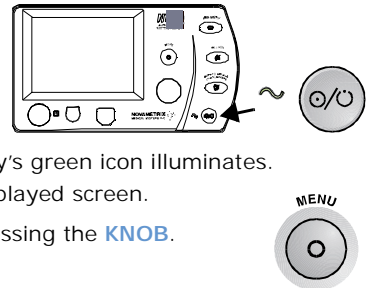

## Ranges and Units

- Analog Input/Output Port (selectable, 0 to 1 volt range):
  - C.O. - Cardiac Output, 0-20 L/m, 50mV/L/m
  - CI - Cardiac Index, 0-20 L/m, 50mV/L/m
  - SV - Stroke Volume, 0-20 L/m, 50mV/L/m
  - PCBF - 0-20 L/m, 50 mV/L/m
  - ETCO<sub>2</sub> - 0-150 mmHg, 0-20 kPa or %, 6.67mV/mmHg
  - SpO<sub>2</sub> - 0-100%, 10mV/%
  - Resp Rate - 0-150 br/min, 6.67mV/br/min
  - Pulse Rate - 0-250 bpm, 4mV/bpm
  - CO<sub>2</sub> Waveform - 0-150 mmHg, 0-20 kPa or %, 6.67mV/mmHg
  - Pleth Waveform - auto scaled
  - Flow Waveform - -125 to +125 L/m, 4mV/L/m
  - Airway Pressure Waveform - -20 to +105 cmH<sub>2</sub>O, 8mV/cmH<sub>2</sub>O

Pin#	Description	Pin #	Description
1	Ground	9	Ground
2	Channel 1 - input (not enabled)	10	Ground
3	Channel 2 - input (not enabled)	11	Channel 1 - output
4	Channel 3 - input (not enabled)	12	Channel 2 - output
5	Channel 4 - input (not enabled)	13	Channel 3 - output
6	Ground	14	Channel 4 - output
7	Ground	15	Input/Output Sense
8	Ground		

## Alerts and Messages

Alarms and messages are not transmitted from the NICO<sup>®</sup> monitor to the Agilent VueLink module; most alarms and messages are not transmitted from the NICO<sup>®</sup> monitor to the GEMS-IT Solar<sup>®</sup> monitor. All alarms and messages should be viewed directly from the NICO<sup>®</sup> monitor.



## Maintenance

---

---

This section details routine maintenance procedures for the NICO® monitor, its sensors and accessories.

### Cleaning and Sterilization

To clean and/or sterilize the monitor and its accessories:

#### *Single Patient Use NICO Sensor™*

- Treat the NICO Sensor™ in accordance with hospital protocol for single-patient use items.

#### *CO<sub>2</sub>/Flow Sensors*

- Treat CO<sub>2</sub>/Flow sensors in accordance with hospital protocol for single-patient use items.

#### *CAPNOSTAT® CO<sub>2</sub> Sensor*

- Do not immerse the sensor. Do not sterilize the sensor.
- The sensor can be cleaned and disinfected by wiping with solutions such as a 70% isopropyl alcohol, 2% glutaraldehyde, or 10% bleach solution. Then wipe down with a water dampened clean cloth to rinse. Dry before use.
- Make certain that the sensor windows are clean and dry before reuse.

#### *NICO® Monitor*

- Do not immerse the monitor. Do not sterilize the monitor.
- Turn the monitor off and unplug from the AC power source before cleaning.
- The monitor can be cleaned and disinfected by wiping with solutions such as a 70% isopropyl alcohol, 2% glutaraldehyde, or 10% bleach solution. Then wipe down with a water dampened clean cloth to rinse. Dry before use.

#### *SpO<sub>2</sub> Finger Sensor*

- Do not immerse the finger sensor. Do not sterilize the finger sensor.
- The sensor can be cleaned and disinfected by wiping with solutions such as a 70% isopropyl alcohol, 2% glutaraldehyde, or 10% bleach solution. Then wipe down with a water dampened clean cloth to rinse. Dry before use.
- Make certain that the finger sensor windows are clean and dry before reuse.
- After cleaning the finger sensor, perform a Quick Check to verify the sensor is functional (See "Sensor Quick Check" on page 68).

#### *SpO<sub>2</sub> Y-Sensor™*

- The Y-Sensor™ may be immersed up to—but not including—the connector, in a 2% glutaraldehyde solution, or 10% bleach solution. Refer to manufacturer's instructions and standard hospital protocols to determine recommended times for disinfection and sterilization.
- Rinse thoroughly with water and dry before use. (Do not rinse the connector).
- After cleaning or sterilizing the Y-Sensor™, perform a Quick Check to verify the sensor is functional (See "Sensor Quick Check" on page 68).

#### *SpO<sub>2</sub> Tapes and Foam Wraps*

- Treat tapes and foam wraps in accordance with hospital protocol for single-patient use items.

### Ear Clip

- Clean with a cloth dampened with 70% isopropyl alcohol. After cleaning, thoroughly wipe the ear clip with a clean, water dampened cloth.

## Monitor Maintenance Schedules

The NICO® monitor performs a diagnostic self-test at powerup that checks the internal electronics. If this self test fails, the normal monitoring display will not appear. Remove the NICO® monitor from use and contact qualified service personnel.

The NICO® monitor should undergo inspection and safety checks on a regular basis or according to institutional protocol. A Service Manual (Catalog No. 9226-90) containing information to assist qualified service personnel is available.

## Battery Maintenance

The monitor may not power up on battery power if the battery is not sufficiently charged. If the NICO® monitor has not been used or powered by the AC Mains for an extended time—3 months or more—allow the battery to charge for 12 hours before use. (The internal battery may slowly discharge over long periods of non-use.)

To charge the battery, connect the line cord to an AC source and set the rear panel power switch ON (I). Check that the AC Mains Power Indicator on the front panel is illuminated (green). Allow the battery to charge for 12 hours. (Refer battery replacement to qualified service personnel.)



## Specifications

---

### General

Specifications for the Novametrix NICO<sup>®</sup> Monitor, Model 7300, are listed for informational purposes only, and are subject to change without notice.

### Cardiac Output

- Measurement Frequency: Rebreathing cardiac output measurement made every three minutes, rebreathing period is 50 seconds.
- Cardiac Output Range: 0.5-19.9 liters/minute
- Cardiac Output Resolution: 0.1 liters/minute
- Pulmonary Capillary Blood Flow (PCBF) Range: 0.5-19.9 L/min, Resolution: 0.1 L/min
- Cardiac Index Range: 0-9.9 L/min/meter<sup>2</sup>, Resolution: 0.1 L/min/meter<sup>2</sup>
- Stroke Volume Range: 0-250 ml, Resolution: 1 ml
- Rebreathing Valve/sensor:
  - Valve type: dual diaphragm, pneumatically controlled
  - Return spring: automatically returns valve to normal position
  - Resistance: 3 cmH<sub>2</sub>O/L/min maximum
  - Rebreathed volume: normal position 35 ml; rebreathing position 150-450 ml (std.)
  - CO<sub>2</sub>/flow sensor: integrated into valve assembly
- Parameter limits for NICO<sup>®</sup> measurements:
  - VCO<sub>2</sub>: >20 ml/min
  - RR: >3, <60
  - Vt: >200 (small and standard), >400 (large)
  - ETCO<sub>2</sub>: >15, <85 mmHg (<100 mmHg during rebreathing)  
>2.0, <11.5 kPa or % (<13.5 kPa or % during rebreathing)

### CO<sub>2</sub>

- Principle of Operation: Non-Dispersive Infrared (NDIR) absorption, dual wavelength ratiometric-single beam optics, mainstream sensor.
- Response Time: Less than 60 ms
- Gas composition effects: Operator selectable
- CAPNOSTAT<sup>®</sup> CO<sub>2</sub> Sensor:
  - Weight: Less than 18 g without cable
  - Sensor Size: 1.3 x 1.67 x .85 inches (3.3 x 4.2 x 2.2 cm), 8 foot cable (2.44 m)
  - Construction: Durable high performance plastic, ultra-flexible cable. Shock Resistant: Sensor will withstand a 6 foot drop to a tile floor.
- End Tidal CO<sub>2</sub>:
  - Range: 0-150 mmHg, 0-20 kPa or % at Pb 760 mmHg
  - Accuracy: ± 2 mmHg for 0-40 mmHg, ± 5% of reading for 41-70 mmHg, ± 8% of reading for 71-150 mmHg
- Respiratory Rate:
  - Range: 2-150 breaths/min
  - Accuracy: ± 1 breath/min

### Flow

- Flow Range (L/min) at Pb 760 mmHg, room air, 35°C
  - Adult: 2 to 180
  - Pediatric: .5 to 100
  - Neonatal: .25 to 25
- Flow Accuracy: Greater of ± 3% reading or:
  - Adult: .5 L/min
  - Pediatric: .25 L/min
  - Neonatal: .125 L/min
- Tidal Volume Range (ml)
  - Adult: 200 to 3000
  - Pediatric: 30-400



- Neonatal: 1-100
- Airway Pressure Range (cmH<sub>2</sub>O): ± 120
- Accuracy: greater of 0.5 cmH<sub>2</sub>O or ± 2% reading

## SpO<sub>2</sub>

- Oxygen Saturation
  - Range: 0-100%
  - Accuracy: ± 2% for 80-100% (for 1 standard deviation or approximately 68% of readings), unspecified for 0-79%
  - Averaging Time: 2 seconds
- Pulse Rate:
  - Range: 30-250 beats per minute
  - Accuracy: ± 1% of full scale (for 1 standard deviation or approximately 68% of readings)
  - Averaging Time: 8 seconds

## Monitor Specifications

- Classification (IEC601-1): Class I/internal power source, type BF, continuous operating mode, enclosure protection rating IPX0.
- Operating Environment: 50-95° F (10-35° C), 0-90% relative humidity (non-condensing)
- Size: Height 6.5 in., Width 10.75 in., Depth 9.5 in.
- Weight: 9 lbs, 6 oz.
- Power: 100-240 VAC, 50-60 Hz, 70VA
- Fuse Rating: 100-240 VAC, 0.5 A 250 V Slo-Blo (x2); 200-240 VAC, T 250 mA/250 V (x2)
- Battery: Internal, Sealed lead-acid gel-cell, 45 minute life on full charge (on-screen life indicator), 12 hours recharge time.
- Display: 4.625 x 3.5 inch EL, 320x240 pixels
- Electromagnetic Emissions: Conforms to EMC Directive 89/336/EEC, CISPR Class A. Tested to EN55011 (1991) and CISPR11 (1990).
- Electromagnetic Immunity: Conforms to EMC Directive 89/336/EEC, EN50082-1 (1992). Tested to IEC801-3 (1984) Radiated Immunity. Conforms to Medical Device Directive 93/42/EEC EN60601-1 (1992). Tested to IEC801-2 (1991) ESD, IEC801-4 (1988) EFT, and IEC1000-4-5 (1995) Surge Immunity.

## RS232 Communications

- RS232 Communications Ports:

Pin #	RS232-1	RS232-2	RS232-3
2	Rx	Rx	Rx
3	Tx	Tx	Tx
5	Ground	Ground	Ground
7	n/a	RTSB	n/a
8	n/a	CTSB	n/a
9	n/a	Power	n/a



## NICO<sup>®</sup> Accessories

Catalog No.	Description
9226-00	<b>NICO<sup>®</sup> Non-Invasive Cardiac Output Monitor</b> , Model 7300 Includes: Monitor, CAPNOSTAT <sup>®</sup> CO <sub>2</sub> Sensor, SpO <sub>2</sub> Sensor, Power Cord and User's Manual. Warranty for NICO <sup>®</sup> Monitor and CAPNOSTAT <sup>®</sup> CO <sub>2</sub> Sensor is 2 years.
8950-00	<b>NICO Sensor<sup>™</sup></b> (10 per box) <b>Small size (for tidal volumes of 200 - 500 mL)</b>
8951-00	<b>NICO Sensor<sup>™</sup></b> (10 per box) <b>Standard size (for tidal volumes of 400 - 1000 mL)</b>
8952-00	<b>NICO Sensor<sup>™</sup></b> (10 per box) <b>Large size (for tidal volumes of 750 - 1500 mL)</b>
9567-00	<b>CAPNOSTAT<sup>®</sup> CO<sub>2</sub> Sensor - NICO</b>
6934-00	<b>Cable Management Straps</b> for use with the CAPNOSTAT <sup>®</sup> CO <sub>2</sub> Sensor. Organizes and holds multiple cables and tubings (package of 5)
8751-00	<b>CAPNOSTAT<sup>®</sup> CO<sub>2</sub> Sensor Cable Holding Clips</b> (50 per box)
8776-00	<b>SuperBright<sup>™</sup> Finger Sensor</b> (10 ft. sensor cable) 1 yr. warranty
8791-00	<b>SuperBright<sup>™</sup> Y-Sensor</b> (10 ft. sensor cable) 90 day warranty
9765-00	<b>NICO<sup>®</sup> CO<sub>2</sub>/Flow Sensors</b> (10 per box) <b>Neonatal size</b>
9766-00	<b>NICO<sup>®</sup> CO<sub>2</sub>/Flow Sensors</b> (10 per box) <b>Pediatric size</b>
9767-00	<b>NICO<sup>®</sup> CO<sub>2</sub>/Flow Sensors</b> (10 per box) <b>Adult size</b>
6063-01	<b>Single Patient Use Airway Adapter</b> - (1 piece) Adult size
6455-00	<b>Single Patient Use SpO<sub>2</sub> Sensor (10 per box) Pediatric/Adult</b>
6455-25	<b>Single Patient Use SpO<sub>2</sub> Sensor (25 per box) Pediatric/Adult</b>
6480-00	<b>Single Patient Use SpO<sub>2</sub> Sensor (10 per box) Neonatal/Pediatric</b>
6480-25	<b>Single Patient Use SpO<sub>2</sub> Sensor (25 per box) Neonatal/Pediatric</b>
4941-00	<b>Saturation Sensor Extension Cable</b> (4 feet)
4942-00	<b>Saturation Sensor Extension Cable</b> (6 feet)
4943-00	<b>Saturation Sensor Extension Cable</b> (10 feet)
5266-00	<b>Saturation Sensor Extension Cable</b> (25 feet)
6147-00	<b>Saturation Sensor Extension Cable</b> (50 feet)
8828-00	<b>20mm Wrap Style Y-Strip Taping System</b> (100 per box) Neonatal foot and hand, pediatric toe or finger, color coded blue
8829-00	<b>25mm Wrap Style Y-Strip Taping System</b> (100 per box) Neonatal foot and hand, color coded green
8831-00	<b>20mm Finger Style Taping System</b> (100 per box) Use on pediatric finger or on small adult finger, color coded blue
8832-00	<b>25mm Finger Style Taping System</b> (100 per box) Use on adult finger, color coded green
6929-00	<b>Adhesive Foam Wraps</b> , Large (25 per box)
6968-00	<b>Adhesive Foam Wraps</b> , Small (25 per box)
8836-00	<b>Non-Adhesive Foam Wraps</b> , Large (25 per box)
8943-00	<b>Non-Adhesive Foam Wraps</b> , Small (25 per box)
6131-50	<b>Ear Clips</b> (5 per box)
6131-25	<b>Ear Clips</b> (25 per box)
8700-00	<b>Adhesive Dots</b> (200 per box)

Catalog No.	Description
9861-00	<b>Flow Connector Cap</b> , 3 port (25 per bag)
601012	<b>NICO</b> <sup>®</sup> 9-pin to 9-pin null modem (crossover) cable
600026	<b>Power Cord</b> (included with monitor)
9226-23	<b>NICO</b> <sup>®</sup> User Manual
9226-90	<b>NICO</b> <sup>®</sup> Service Manual



## Warranty

---

---

Equipment manufactured or distributed by Novamatrix Medical Systems Inc., is fully guaranteed, covering materials and workmanship, for a period of one year from the date of shipment, except for certain disposable products and products with stated guarantees other than one year. Novamatrix reserves the right to perform guarantee service(s) at its factory, at an authorized repair station, or at the customer's installation.

Novamatrix' obligations under this guarantee are limited to repairs, or at Novamatrix' option, replacement of any defective parts of our equipment, except fuses, batteries, and calibration gasses, without charge, if said defects occur during normal service.

Claims for damages during shipment must be filed promptly with the transportation company. All correspondence concerning the equipment must specify both the model name and number, and the serial number as it appears on the equipment.

Improper use, mishandling, tampering with, or operation of the equipment without following specific operating instructions will void this guarantee and release Novamatrix from any further guarantee obligations.

### Service Department

For factory repair service:  
Call toll free: 1-800-243-3444  
To Call Direct: (203) 265-7701  
Facsimile (203) 284-0753  
<http://www.novamatrix.com>  
[techline@novamatrix.com](mailto:techline@novamatrix.com)

**Caution:** Federal (U.S.A.) law restricts this device to sale, distribution, or use by or on the order of a licensed medical practitioner.

Copyright © 2001. Novamatrix Medical Systems Inc.,  
5 Technology Drive, Wallingford, Connecticut, 06492.  
All rights reserved.